Imtron GmbH
 Media-Saturn-Strasse 1
 85053 Ingolstadt
 Germany
www.imtron.eu
www.koenic-online.com

IM_KGM 5424_231127_V04

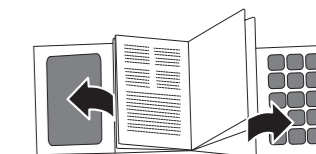
KOENIC

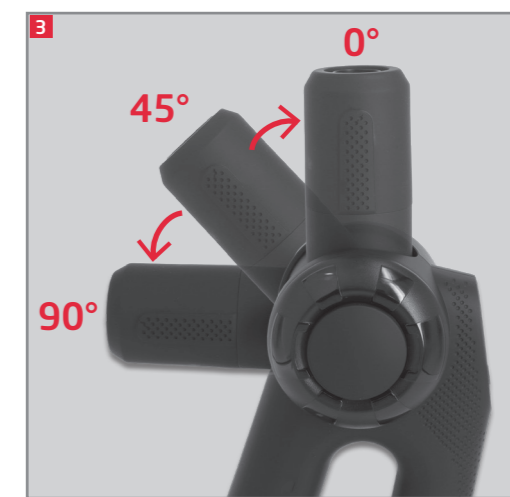
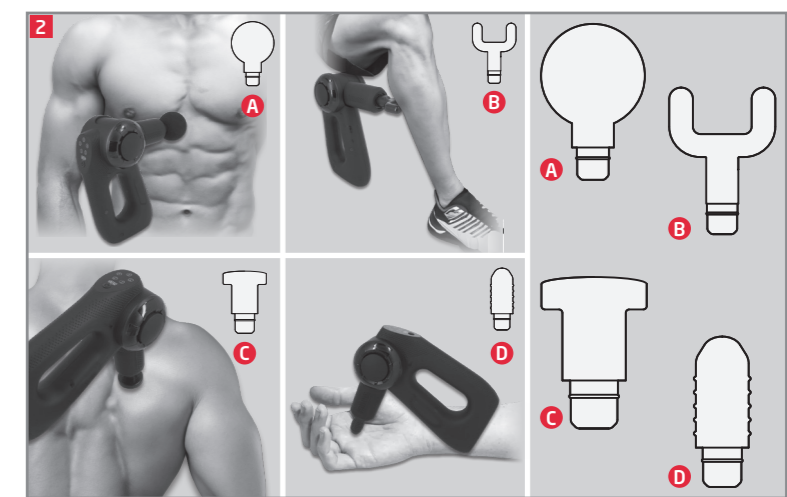
Massage Gun KGM 5424

- DE Gebrauchsanweisung
- EN User manual
- ES Manual de instrucciones
- FR Mode d'emploi
- HU Használati leírás
- IT Manuale dell'utente
- NL Gebruiksaanwijzing
- PL Instrukcja obsługi
- TR Kullanım Kılavuzu

Deutsch	3 - 8	—
English	9 - 14	—
Español	15 - 19	—
Français	20 - 25	—
Magyar	26 - 31	—
Italiano	32 - 37	—
Nederlands	38 - 43	—
Polski	44 - 49	—
Türkçe	50 - 54	—

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Safety instructions

-  **Warning!** Do not use this product in the vicinity of water or high levels of humidity, for example in a damp cellar or close to a swimming pool or bathtub.
- This product can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they are supervised or have received instruction concerning use of the product in a safe way and understand the hazards involved.
- Children shall not play with the product.
- Cleaning and user maintenance shall not be carried out by children unless they are supervised.
- **Warning!** Keep the product dry.
-  **Warning!** For the purposes of recharging the battery, only use the detachable supply unit IVP0500-2000U provided with this appliance. The polarity of the adaptor must match the polarity of the product .
-  The power supply unit shall not be used, if pins of the plug part are damaged.
- Never expose rechargeable batteries to excessive heat (e.g. bright sun, fire), and never throw them into fire. The rechargeable batteries could explode.
- The rechargeable battery is permanently built into the product and cannot be replaced.

Safety instructions

- When disposing of the product, please consult an electrician or your retailer to remove and dispose of the built-in batteries in an environmentally friendly manner.
- Never short-circuit the contacts of the rechargeable battery. Do not throw the rechargeable battery or the product into fire. There is a danger of fire and explosion!
- Never charge the rechargeable battery of the product unattended.

Safety instructions

- Read this user manual thoroughly before first use and pass it on with the product. Pay attention to the warnings on the product and in this user manual. It contains important information for your safety as well as for the use and maintenance of the equipment.
- Use the product only for its intended purpose and with recommended accessories and components. Improper use or wrong operation may lead to hazards.
- Never leave the product unsupervised while it is switched on.
- Never immerse electrical parts of the product in water during cleaning or operation. Never hold the product under running water.
- Do not under any circumstances attempt (e.g. damaged power cord, product has been dropped, etc.) to repair the product yourself. For servicing and repairs, please consult an authorised service agent.
- Use this product in dry indoor areas only and never outdoors.
- Never use, expose or put the product near to:
 - direct sunlight and dust;
 - fire (fireplace, grill, candles), water (water splashes, vases, bath tube) or heavy moisture.
- This product is not intended for commercial use. It is designed for household use only.
- Clean the product and all other accessories after each use.
- Follow the instructions in the chapter **Cleaning and care**.
- This product is intended for private use only and cannot replace physiotherapy or medical treatment. Do not use the product for therapeutic or medical purposes.
- The product must not be used:
 - on people with pacemakers, heart problems or pain of unknown origin.
 - for febrile colds, varicose veins, thromboses, phlebitis, jaundice, diabetes, nervous diseases such as sciatica, acute inflammation or during pregnancy.
 - on swollen, inflamed or injured parts of the body.
 - for massaging bones, especially the spine and head. Only use the product on muscles.
 - near the eye area or other sensitive parts of the body.
 - Around the abdomen or lower back for pregnant women.
 - on small children.
- If you are unsure, ask your doctor for advice before use.
- While using the product, stop and contact a physician or doctor if you experience unusual or unfamiliar pain.
- Do not use the product under a blanket or pillow. The product may overheat and increase the risk of fire.
- Only attach or replace a massage head when the product is turned off.
- **Important!** Check to ensure that the product is in a serviceable condition before use. The product shall not be used if it shows any signs of damage or if it's suspected to be defect.

English

Congratulations!

Thank you for your purchase of a **KOENIC** product. Please read this manual carefully and keep it for future reference.

Intended use

This device is only intended for massaging parts of the human body. This product is designed for household use only and not intended for commercial use. Any other use may result in damage to the product or injuries. The Imtron GmbH assumes no liability for damage to the product, for property damage or for personal injury due to careless or improper usage of the product, or usage of the product, which does not meet the manufacturer's specified purpose.

Before first time use

Carefully take the product and its accessories out of the original packaging. Retain the original packaging in case you would like to store the product away or transport it at a later date. If you choose to dispose of the original packaging, then do so according to the applicable legal regulations.

Check the packaging contents for completeness and damages. Should the contents in the packaging not be complete or damaged, do not use the product. Return it to your retail outlet immediately.

After unpacking, please refer to chapter

Cleaning and care.

Note:

- For first time use, charge the product for at least 3 hours. The general charging time is about 2 hours which also depends on the condition of the battery.
- When the battery is fully charged, the product can be used for 120-360 minutes. The time may vary (depending on the battery condition).

Content

- 1 x Massage device
- 1 x USB charger
- 1 x USB cable
- 4 x Massage attachments
- 1 x Storage case
- 1 x Storage bag
- 1 x User manual

Controls and components

- A** Storage case
- B** Storage bag
- C** Joint release button
- D** Massage level buttons with indicators
- E** Battery status indicator
- F** On/Off button
- G** Handle
- H** USB C charging socket
- I** USB charging cable
- J** Charger
- K** Small round head for all parts of the body
- L** Round head for joints
- M** Ball head for large muscle groups
- N** U-head for vertebrae
- O** Massage attachment socket

Preparation

Fig. 1 Connect the USB C plug with the USB C charging socket, the USB A plug with the charger and the charger with a suitable socket outlet.

When the unit is charging, the battery indicator will flash. When the battery is fully charged, the battery indicator will light up steadily. Charge the battery for at least 4 hours before first use.

Low battery warning: When the battery is less than 10% the indicator light flashes white.

Operation

Fig. 2 Replacing or attaching the massage attachment.

- 1 Select the desired massage attachment.
- 2 Remove any massage attachment that may be attached by pulling it off.
- 3 Insert the new massage attachment into the holder.

Please note the following information to select the right massage attachment.

A	Round head: Arm, waist, hip, thigh, calf and other large muscle groups
B	Double/U-head: Relaxes and shapes different muscle groups
C	Flat head: Neck, spine and Achilles tendon
D	Precision head: Deep tissue such as joints, palms and soles

Adjusting the massage angle

- Fig. 3
- Turn off the massager before adjusting the angle.
 - Press the joint release button and adjust the massage arm to the desired position and release the buttons again.
 - Three angles are available: 0°, 45° and 90°.

Power on and mode selection

- Fig. 4
- Press the on/off button to switch on the device. The device now runs at the lowest level 1.
 - Select a massage level from 1-5 by pressing the corresponding button.
 - Press and hold the on/off button for about 3 seconds to turn off the device.

Massage

- Hold the massager by the handle and place the massage head on the skin.
- Move the massage head slowly over the skin in circular or straight movements.
- Massage without applying pressure.
- Keep the skin as taut as possible to achieve optimal effect.

Fig. 5 To guide the massager, please note the series of images.

Note:

A complete massage cycle takes 15 minutes. The product switches off automatically after 15 minutes.

Do not use the product for more than 2 massage cycles (2 x 15 minutes = 30 minutes) continuously, but switch it off for at least 30 minutes before switching it on again.

English

Cleaning and care

Warning

When cleaning, never use solvents or abrasive materials, hard brushes, metallic or sharp objects. Solvents are harmful to human and can attack plastic parts, while abrasive cleaning mechanisms and tools may scratch the surface(s).

- Clean the product's exterior surfaces with a slightly damp cloth and then dry them thoroughly.
- Make sure that no water enters the product interior (for example, through the switch). Keep the product in a dry, cool location which is protected against moisture and direct sunlight. Keep this product out of the reach of children.

Technical data

Massage device

Rated input voltage / current : 5.0 V --- / 2.0 A
 IP code : IPX0
 Re-chargeable battery : Li-ion, 11.1 V --- , 2000 mAh, 22.2 Wh

Charger

Manufacturer : Xiamen Invcy Electronics Co., Ltd.
 Imported by : Imtron GmbH
 Media-Saturn-Strasse 1, 85053 Ingolstadt, Germany
 Commercial registration number : HRB 4580
 Model : IVPO500-2000U
 Rated input voltage : 100 - 240 V~
 Input AC frequency : 50/60 Hz
 Output voltage : 5.0 V ---
 Output current : 2.0 A
 Output power : 10.0 W
 Average active efficiency : 80.90 %
 Efficiency at low load (10 %) : 75.37 %
 No-load power consumption : 0.077 W
 Protection class : IP20

Troubleshooting

If there are any problems during operation, please refer to the following table. If the problems cannot be solved as indicated under **Solution**, please contact your sales agent.

Problem	Possible causes	Solution
Device does not work.	<ul style="list-style-type: none"> • Interrupted power supply • Batteries empty 	<ul style="list-style-type: none"> • Check the power supply • Re-charge batteries
LED indicator does not light on while charging batteries	<ul style="list-style-type: none"> • No power supply • Power cord is not connected to product or charging stand 	<ul style="list-style-type: none"> • Check the power supply • Check the power cord connection

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Üretici : Xiamen Invcy Electronics Co., Ltd.
İthal eden : Imtron GmbH