

**PHILIPS**

**sonicare**

Advanced Clean

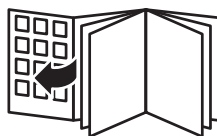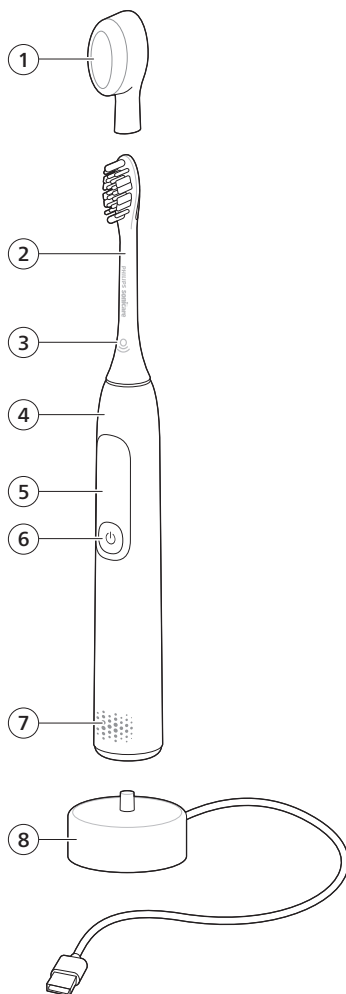English 6  
Azərbaycanca 25  
Čeština 44  
Dansk 63  
Deutsch 82  
Eesti 105  
Español 124  
Français 144  
Hrvatski 165  
Italiano 184  
Latviešu 205  
Lietuviškai 224  
Magyar 243  
Nederlands 263  
Norsk 284  
O'zbek 303  
Polski 323  
Português 344  
Română 364  
Shqip 383  
Slovenščina 403  
Slovensky 422  
Srpski 441  
Suomi 460  
Svenska 479  
Türkmençe 498  
Ελληνικά 517  
Български 539  
Кыргызча 560  
Македонски 581  
Русский 601  
Тоҷикӣ 623  
Українська 645

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## Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at **[www.philips.com/welcome](http://www.philips.com/welcome)**.

## Important safety information

Only use the product for its intended purpose. Read this information carefully before you use the product and its batteries and accessories, and save it for future reference. Misuse can lead to hazards or serious injuries. Accessories supplied may vary for different products.

## Warnings

- Only use original Philips accessories or consumables. If supplied, only use the original USB cable and USB power supply unit.
- This appliance contains no user-serviceable parts. If the appliance is damaged, stop using it and contact the Consumer Care Center in your country (see 'Warranty and support').
- If the USB cable or USB power supply unit is damaged, stop using it.
- Keep the charger away from water.
- This appliance can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
- Children should be supervised to ensure that they do not play with the product.
- Do not use adapters and chargers outdoors or near heated surfaces.
- Do not clean any part of the product in the dishwasher.
- After cleaning, make sure the charger is completely dry before you connect it.
- Rinse the handle, especially the brush head connection. Gently clean around the rubber seal.

- Replace the brush head every 3 months or sooner if signs of wear appear. Stop using a brush head with crushed or bent bristles.
- This appliance has only been designed for cleaning teeth, gums and tongue.
- Avoid direct contact with products that contain essential oils or coconut oil. Contact may result in bristles dislodging.
- The maximum altitude of use is 4500 meters.
- Information on REACH substances, contained in Philips products, can be found on [www.philips.com/REACH](http://www.philips.com/REACH)

## Medical warnings

- This appliance is a personal care device and is not intended to be shared by multiple users.
- Stop using the appliance and consult your dentist/doctor if excessive or prolonged bleeding occurs after use, or if you experience discomfort or pain during use.
- If you have a pacemaker or other implanted device, contact your doctor or the manufacturer of the implanted device prior to use.

## Battery safety instructions

- This appliance contains batteries that are only removable by skilled persons.
- Keep product and batteries away from fire and do not expose them to direct sunlight or high temperatures.
- If the product becomes abnormally hot or smelly, changes color or if charging takes longer than usual, stop using and charging the product and contact Philips.
- Do not place products and their batteries in microwave ovens or on induction cookers.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent batteries from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge batteries.

- This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of, see section on Recycling.
- Charge only with Philips WAA1001, WAA2001, safety approved IPX4 adapters with input ratings (100-240 V~; 50/60 Hz; 3.5 W) and output ratings (5 V dc; 2.5 W). For the list of adapters, please see **[www.philips.com/support](http://www.philips.com/support)**.
- Charge and use the product at a temperature between 0 °C and 40 °C (32°F and 104 °F).
- Power consumption:
  - Off mode: Not applicable
  - Standby mode: < 0.1 W (automatically reached within 1 minute)

## Storage and transportation conditions

Temperature -10°C to +60°C (14 °F to 140 °F).

## Electromagnetic fields (EMF)

This Philips appliance complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

## Radio Equipment Directive

Hereby, Philips declares that electric toothbrushes with radio equipment type (Bluetooth, NFC) are in compliance with Directive 2014/53/EU.

- The Bluetooth radio frequency interface in applicable products operates at 2.4 GHz.
- The maximum output power by the Bluetooth appliance is 3 dBm.
- The NFC radio frequency interface in applicable products operates at 13.56 MHz.
- The maximum RF power transmitted by the appliance is 30.16 dBm.

The full text of the EU declaration of conformity is available at the following internet address:

**[www.philips.com/support](http://www.philips.com/support)**

Note: The Features of different models may vary. Some models may not be equipped with Bluetooth or NFC.

## LIST OF SYMBOLS USED

The following symbols may appear on the product:  
This symbol means read the instructions for use, before using the appliance.

This symbol means to use the indicated power supply.

## Intended use

This Sonicare power toothbrush is intended to remove adherent plaque and food debris from the teeth to reduce tooth decay and improve and maintain oral health. This Sonicare power toothbrush is intended for consumer home use. Use by children should be with adult supervision.

## Your Philips Sonicare (Fig 1.)

- 1 Brush head cap
- 2 Smart brush head
- 3 BrushSync symbol
- 4 Handle
- 5 Smart display
- 6 Power/mode button
- 7 Pressure sensor indicator
- 8 USB charger (wall adapter not included)

Note: Accessories included may vary based on the model purchased.

## Getting started

### Attaching the brush head

- 1 Align the brush head so the bristles point in the same direction as the front of the handle.

- 2 Firmly press the brush head down onto the metal shaft until it stops.

Note: It is normal to see a slight gap between the brush head and the handle. This allows the brush head to vibrate properly.

### Using your Philips Sonicare toothbrush

If you are using your Sonicare toothbrush for the first time, it is normal to feel more vibration than when using a non-electric toothbrush. It is common for first time users to apply too much pressure. Only apply gentle pressure and let the toothbrush do the brushing for you. Please follow brushing instructions below for the best experience.

To help you transition to your new Sonicare power toothbrush, it comes with the EasyStart feature activated. This feature gradually increases the power over the first 14 brushing sessions to help you adapt to brushing with your Sonicare toothbrush.

### Brushing instructions


- 1 Wet the bristles and apply a small amount of toothpaste.


- 2 Press the power button to turn on your Philips Sonicare. Then select a brushing mode.

Note: The first 3 seconds of all modes start with a gentle vibration. This prevents toothpaste from splattering and lets you cycle through the modes to select your preferred brushing mode.

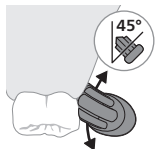
- 3 Place the toothbrush bristles against the teeth at a slight angle (45 degrees), pressing gently to make the bristles reach the gumline or slightly beneath the gumline.

Note: Keep the center of the brush in contact with the teeth at all times.

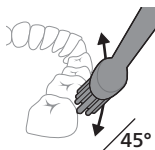
- 4 Gently keep the bristles placed on the teeth and in the gumline. Brush your teeth with small back and forth motion so the bristles reach between the teeth. Continue this motion throughout your brushing cycle.

Note: The bristles should slightly flare. It is not recommended to scrub as you would with a manual toothbrush.

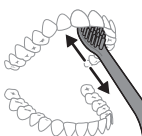
- 5 To clean the inside surfaces of the front teeth, tilt the brush handle semi-upright and make several vertical overlapping brushing strokes on each tooth.

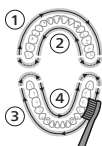- 6 Begin brushing in section 1 (outside upper teeth) and brush for 30 seconds before you move to section 2 (inside upper teeth). Continue brushing in section 3 (outside lower teeth) and brush for 30 seconds before you move to section 4 (inside lower teeth).

Note: Additional time is added when you use the Deep Clean, Gum Pro, or White modes.

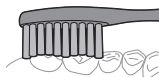- 7** After you have completed the brushing cycle, you can spend additional time brushing the chewing surfaces of your teeth and areas where staining occurs. You can also brush your tongue, with the toothbrush turned on or off, as you prefer.

Your Philips Sonicare toothbrush should be safe to use on braces (brush heads wear out sooner when used on braces) and dental restorations (fillings, crowns, veneers) if they are properly adhered and not compromised. If a problem occurs, please follow up with a dental professional.

**Note:** When the Philips Sonicare toothbrush is used in clinical studies, it must be used in Clean mode with the handle fully charged and the EasyStart feature deactivated.

## Modes

Your power toothbrush is equipped with five different modes: Clean, Deep Clean, Gum Pro, Sensitive, and White.

Based on the smart brush head you use; the optimal mode will be selected automatically by the BrushSync technology.

To manually change the mode, press the power button when the toothbrush is on. All Philips Sonicare brush heads work with each mode.

| Mode       | Clean                                                                               | Deep Clean                                                                          | Gum Pro                                                                             | Sensitive                                                                           | White                                                                               |
|------------|-------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|
| Mode image |  |  |  |  |  |
|            | Clean                                                                               | Deep Clean                                                                          | Gum Pro                                                                             | Sensitive                                                                           | White                                                                               |

| Mode                       | Clean                                                                                                                  | Deep Clean                                                                                                             | Gum Pro                                                                                                                                                                                       | Sensitive                                                                                                              | White                                                                                                                                                                                      |
|----------------------------|------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>Benefit</b>             | Plaque removal                                                                                                         | Superior plaque removal                                                                                                | Plaque removal and gum massage                                                                                                                                                                | Extra gentle mode for sensitive teeth and gums                                                                         | Plaque and surface stains removal                                                                                                                                                          |
| <b>Total brushing time</b> | 2 minutes                                                                                                              | 3 minutes                                                                                                              | 3 minutes                                                                                                                                                                                     | 2 minutes                                                                                                              | 2 minutes and 30 seconds                                                                                                                                                                   |
| <b>How to brush</b>        | Brush each segment for 30 seconds<br> | Brush each segment for 45 seconds<br> | Brush each segment for 30 seconds,<br><br>then brush along the gum line of each segment for 15 seconds each. | Brush each segment for 30 seconds<br> | Brush each segment for 30 seconds,<br><br>then brush the upper and lower front teeth for 15 seconds each. |

## BrushSync technology

BrushSync technology enables your brush head to communicate with your handle using a microchip. The  symbol at the bottom of the brush head indicates that the brush head is equipped with this technology.

BrushSync technology enables:

- Brush head replacement reminder
- Brush head mode pairing

Philips Sonicare offers a wide range of smart brush heads, equipped with BrushSync technology. To explore our full range of toothbrush brush heads, and find the brush head that is best for you, go to

[www.philips.com/toothbrush-heads](http://www.philips.com/toothbrush-heads) for more information.

## Features

- Pressure sensor
- SmarTimer
- Brushing reports
- QuadPacer
- EasyStart
- Brush head replacement reminder
- Brush head mode pairing

### Pressure Sensor

Your Philips Sonicare is equipped with an advanced sensor that measures the pressure you apply while brushing. If you apply too much pressure, the toothbrush provides immediate feedback to indicate that you need to adjust the pressure. When you apply too much pressure during brushing, you will feel a change in vibration and the bottom of the handle will light up in orange.

### Brushing reports

Your Philips Sonicare regularly provides a report on your brushing behavior.

After brushing the display shows a report with the brushing time.

- Brushing time of 2 minutes or more

- Brushing time of less than 2 minutes

- During brushing the display shows the current brushing mode and brushing time

## SmarTimer

The SmarTimer indicates that your brushing cycle is complete by automatically switching off the toothbrush at the end of the brushing cycle.

Dental professionals recommend brushing no less than 2 minutes twice a day.

Note: If you press the power button after you started the brushing cycle, the toothbrush will pause. After a pause of 15 seconds, the SmarTimer resets.

## QuadPacer

The QuadPacer is an interval timer that has a short beep and pause to remind you to brush the 4 sections of your mouth evenly and thoroughly. This Philips Sonicare toothbrush comes with the QuadPacer feature activated.

Note: Additional time is added when you use the Deep Clean, Gum Pro, or White modes.

## EasyStart

This Philips Sonicare model comes with the EasyStart feature activated. The EasyStart feature gently increases the power over the first 14 brushings to help you get used to the brushing with the Philips Sonicare.

When the Philips Sonicare toothbrush is used in clinical studies, it must be used at high intensity with the handle fully charged and the EasyStart feature deactivated.

To deactivate EasyStart (see 'Activating or deactivating features').

## Brush head replacement reminder

Your Philips Sonicare is equipped with BrushSync technology that tracks the wear of your brush head. This

feature gives you the guarantee for the best cleaning and care of your teeth.

**Note:** This feature only works with Philips Sonicare smart brush heads with BrushSync technology.

- 1 When a new smart brush head is attached for the first time, the handle recognizes that you have a Philips brush head with BrushSync technology and the display shows the brush head symbol.
- 2 Over time, based on the pressure you apply and how long the brush head has been used, the handle will track brush head wear in order to determine the optimal time to change your brush head. When the life of the brush head ends, the display shows the brush head replacement reminder symbol and indicates that you should replace the brush head as soon as possible.

**Note:** Replace Philips Sonicare brush heads every 3 months to achieve optimal results. Use only Philips Sonicare replacement brush heads.

## Brush head mode pairing

Brush heads with BrushSync technology are equipped with a microchip which communicates with the handle and automatically pairs it with the recommended mode and intensity. Most brush heads with BrushSync technology will auto pair to Clean Mode. Brush heads with specific benefits will pair with the corresponding mode (when available), examples include:

G3 Premium Gum Care: Gum Pro Mode

S2 Sensitive: Sensitive mode


W Optimal White: White Mode

**Note:** If you choose to change mode from the recommended setting, the brush head will remember the chosen setting for future brushing sessions.

## Charging and battery status

### Charging your Philips Sonicare

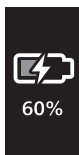- 1 Attach the USB charger to a wall adapter and plug the wall adapter into an electrical outlet.
- 2 Place the handle on the plugged-in charger.
  - The battery icon with lightning symbol is shown on the display and you hear two beeps. The pressure sensor indicator light will also light up white for 5 seconds.

Note: It can take up to 16 hours to fully charge the battery, but you can use the Philips Sonicare toothbrush before it is fully charged.

Note: This Philips Sonicare toothbrush is designed to provide at least 28 brushing sessions, each session lasting 2 minutes (14 days if used twice per day).

### Battery status

The LCD screen will show the battery level after each brushing session, when placing the handle on the charger, and when removing it from the charger. While charging, press the power button to check the battery level.

The battery icon changes color with different battery levels: it turns yellow if the battery level is below 40%; while low battery is indicated by a red icon flashing with two short beeps.

| Battery status | Battery icon color | Battery percentage |
|----------------|--------------------|--------------------|
| Low            | Red                | 0% - 20%           |
| Medium         | Yellow             | 21% - 40%          |
| Full           | Blue               | 41% - 100%         |

### Activating or deactivating features

You can activate or deactivate the following features of your toothbrush:

- EasyStart
- Brush head replacement reminder

To activate or deactivate these features, follow the instructions below:

## EasyStart

- 1 Place the handle on the plugged-in charger.
  - 2 Keep the power button pressed until you hear a single short beep (after 3 seconds).
  - 3 Release the power button.
- Triple tone of low-medium-high means the EasyStart feature has been activated. The icon on the display shows that the EasyStart has been activated.

- Triple tone of high-medium-low means the EasyStart feature has been deactivated. The icon on the display shows that the EasyStart has been deactivated.

**Note:** To achieve clinical efficacy, EasyStart needs to be deactivated.

## Brush head replacement reminder

**Note:** You can only activate or deactivate the brush head replacement reminder when a smart brush head is on the handle.

- 1 Place the handle on the plugged-in charger.
  - 2 Press and hold the power button while the handle remains on the charger.
  - 3 Keep the power button pressed until you hear two short beeps (after 6 seconds).
  - 4 Release the power button.
- Triple tone of low-medium-high means the brush head replacement reminder feature has been activated. The icon on the display shows that the brush head replacement reminder has been activated.


- Triple tone of high-medium-low means the brush head replacement reminder feature has been deactivated. The icon on the display shows that the brush head replacement reminder has been deactivated.

## Cleaning

You should clean your Philips Sonicare toothbrush regularly to remove toothpaste and other residue. Failure to clean your Philips Sonicare toothbrush may result in an unhygienic product and damage to it may occur.

**Warning: Do not clean the product or accessories with sharp objects or in the dishwasher, microwave, with chemicals or in boiling hot water**

### Toothbrush handle

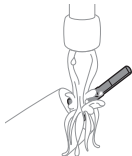- 1 Remove the brush head and rinse the metal shaft area with warm water. Make sure you remove any residual toothpaste on the handle and around the buttons at least once a week, but preferably after each use.


**Caution: Do not push on the rubber seal on the metal shaft or around the buttons with sharp objects, as this may cause damage.**

- 2 Wipe the entire surface of the handle with a damp cloth.

**Caution: Do not tap the handle on the sink or counter to remove excess water.**

- 3 Dry the handle with a soft cloth or towel.

### Brush head

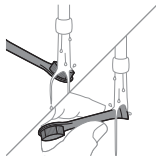- 1 Rinse the brush head and bristles after each use and allow to air dry.
- 2 Remove the brush head from the handle and rinse the brush head connection with warm water at least once a week.

## Charger

- 1 Unplug the charger before you clean it.
- 2 Wipe the surface of the charger with a damp cloth.

## Storage

If you are not going to use the product for a long time, unplug it, and clean it (see chapter 'Cleaning'). Store it in a cool and dry place away from direct sunlight.

Charge the toothbrush regularly to prevent the battery from discharging completely. A completely discharged battery may fail to charge again.

## Locating the model number

Look on the bottom of the Philips Sonicare toothbrush handle for the model number (HX379x).

## Troubleshooting

This chapter summarizes the most common problems you could encounter with the appliance. If you are unable to solve the problem with the information below, visit [www.philips.com/support](http://www.philips.com/support) and insert your model number for a list of frequently asked questions or contact the Consumer Care Center in your country.

| Problem                                                                                 | Possible Cause | Solution                                                                                                                                                    |
|-----------------------------------------------------------------------------------------|----------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------|
| I cannot fully attach the brush head. There is a gap between the brush head and handle. |                | This gap is normal and is needed for the brush head to vibrate appropriately. The brush head must be able to move to create the right amount of vibrations. |

| Problem                                                                | Possible Cause                                                                                   | Solution                                                                                                                                                                                                                                                                                            |
|------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| My Philips Sonicare toothbrush vibration is less powerful than before. | The brush head is too close to the body of the toothbrush or the toothbrush needs to be charged. | Remove the brush head from the handle and re-set it on the metal shaft leaving a small gap. (1-2mm).<br>When charging, make sure the charger is plugged into a power supply and placed on a flat surface before placing the toothbrush on the charger. Allow the toothbrush to charge for 24 hours. |
| I charged my toothbrush, but it only works for a short time.           | The base of the toothbrush is not placed flat on the charger when attempting to charge.          | Make sure the charger is plugged into a power supply and placed on a flat surface before placing the toothbrush on the charger. Allow the toothbrush to charge for 24 hours.                                                                                                                        |

## Warranty and support

If you need information or support, please visit **[www.philips.com/support](http://www.philips.com/support)** or read the international warranty leaflet.

## Warranty restrictions

The terms of the international warranty do not cover the following:

- Brush heads.
- Damage caused by use of unauthorized replacement parts.
- Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
- Normal wear and tear, including chips, scratches, abrasions, discoloration or fading.

## Recycling

- This symbol means that electrical products and batteries shall not be disposed of with normal household waste.
- Follow your country's rules for the separate collection of electrical products and batteries.

## Removal of built-in rechargeable battery

**CONTAINS LITHIUM-ION RECHARGEABLE CELL.  
MUST BE RECYCLED OR DISPOSED OF PROPERLY.**

The built-in rechargeable battery must be removed by a qualified professional when the product is discarded. Before removing the battery, make sure that the battery is completely empty.

**Take any necessary safety precautions when you handle tools to open the appliance and when you dispose of the rechargeable battery.**

**When you handle batteries, make sure that your hands, the product and the batteries are dry.**


**If batteries are damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.**

**To avoid accidental short-circuiting of batteries after removal, do not let battery terminals come into contact with metal objects (e.g. coins, hairpins, rings). Do not wrap batteries in aluminum foil. Tape battery terminals or put batteries in a plastic bag before you discard them.**


Before removing the battery, make sure that the battery is completely empty.

To remove the rechargeable battery, you need a flat-head (standard) screwdriver.


- 1 To deplete the rechargeable battery of any charge, remove the handle from the charger, turn on the Philips Sonicare and let it run until it stops. Repeat this step until you can no longer turn on the Philips Sonicare.
- 2 Remove the film at the bottom of handle, then release the screw.

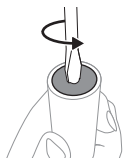
- 3 Remove the end cap from the toothbrush handle.

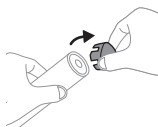
- 4 Holding the handle upside down, press the shaft down on a hard surface.

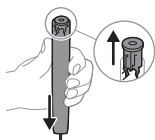
- 5 Wedge the screwdriver between the battery and the white frame at the bottom of the internal components. Then pry the screwdriver away from the battery to break the bottom of the white frame.

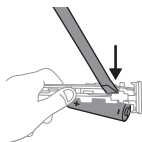- 6 Insert the screwdriver between the bottom of the battery and the white frame to break the metal tab connecting the battery to the green printed circuit board. This will release the bottom end of the battery from the frame.

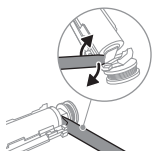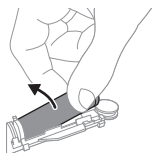- 7 Grab the battery and pull it away from the internal components to break the second metal battery tab.

**Caution: Be aware of the sharp edges of the battery tabs so as to avoid injury to your fingers.**

- 8 Cover the battery contacts with tape to prevent any electrical short from residual battery charge. The rechargeable battery can now be recycled and the rest of the product discarded appropriately.

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©2025 Koninklijke Philips N.V. (KPNV). All rights reserved. Philips and the Philips shield are trademarks of KPNV. Sonicare and the Sonicare logo are trademarks of Philips Oral Healthcare, LLC and/or KPNV.