
ok.

OPS 15422 B

Personenwaage // Personal Scale //  
Báscula para personas // Pèse-personne

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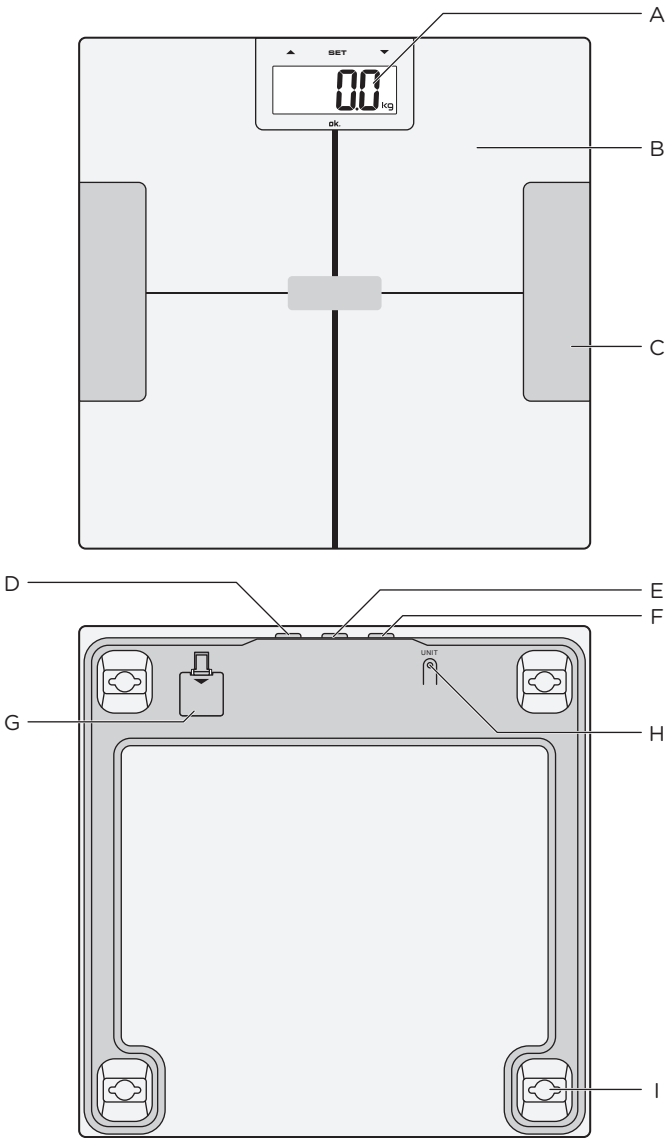EN

CONGRATULATIONS

Thanks for your purchase of an **ok.** product. Please read this manual carefully and keep it for future reference.

**Important safety instructions. READ CAREFULLY AND KEEP FOR FUTURE REFERENCE.**

- This product can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they are supervised or have received instruction concerning use of the product in a safe way and understand the hazards involved.
- Children shall not play with the product.
- Place this product on a horizontal and even surface only.
- Do not place the product on moist or wet surfaces, as there is a danger that the scale will slide away or slip when you stand on it.
- Caution!** Be careful not to damage the glass surface. Do not use this product if the glass surface is damaged or if the glass is cracked.
- Do not stand on the scales with wet or damp feet. Do not stand on the scales if you are wet or not dried off (e.g. after a shower). There is a danger that you may slip and injure yourself.
- Danger of injury!** Do not stand too close to the edges of the scale as this could cause the product to tip over.
- Danger of injury!** Do not jump on or off the product as the glass may shatter.
- Not suitable for pregnant women.
- Not suitable for persons with pace makers.
- Not suitable for children up to 8 years.
- Remove the batteries from the product before disposal.
- Follow the instructions in the chapter **Cleaning and care.**


- This product is not intended for commercial use. It is designed for household use only.
- Read this user guide carefully before first use. It contains important safety information as well as instructions for the use and care of the device. Save the guide for later reference and include it with the device when transferring it to a third party.
- To avoid dangerous situations, do not use the product for any purpose other than those described in the instructions. Improper use is dangerous and voids the warranty.
- Risk of injury!** Do not use this product if it is damaged or has fallen on the floor.

Battery warnings

- Caution!** There is danger of explosion from improper replacement of the battery. Replace the battery only with the same or equivalent type of battery.
- Make sure the battery is correctly inserted. Observe the polarity (+) plus / (-) negative on the battery and the housing. An improperly inserted battery can lead to leakage or, in extreme cases, a fire or an explosion.
- Never expose batteries to excessive heat (e.g. bright sun, fire), and never throw them into fire. The batteries could explode.
- The connections must not be short-circuted.
- Remove used batteries. Remove the batteries if you know that the product will not be used for a longer period of time. Otherwise, the batteries could leak and cause damage.
- Keep batteries out of reach of small children. If batteries have been swallowed, seek immediate medical attention.
- Non-rechargeable batteries must not be recharged.
- If the batteries leak, remove them with a cloth and dispose of them accordingly. Prevent battery liquid from coming into contact with skin and eyes. If you do get battery liquid in your eyes, flush them thoroughly with lots of water and consult a physician immediately. If battery liquid comes into contact with your skin, wash the affected area with lots of water and soap.

Warning!

- Do not swallow the battery due to risk of chemical burn.
- This product contains a coin cell battery. If the coin cell battery is swallowed, it can cause severe internal burns within 2 hours and lead to death.

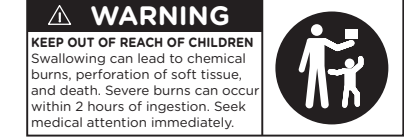- Keep new and used batteries out of reach of children.
- Do not use the product if the battery compartment cannot be closed securely and keep it out of reach of children.
- If you think that batteries have been swallowed or are in a part of the body, seek medical help immediately.
- Do not allow children to change batteries without supervision.
- Do not short-circuit batteries. Do not over-discharge the battery.
- Remove used batteries from the device immediately and dispose of them properly.
- Do not heat batteries. Do not weld or solder directly to batteries. Do not disassemble batteries. Do not deform batteries.
- A lithium battery with a damaged case should not come into contact with water.
- Do not encapsulate and / or modify batteries.

INTENDED USE

This product is intended as a personal scale. The scale has a weighting range of 2 kg up to 150 kg. Any other use may result in damage to product or injuries. Imtron GmbH assumes no liability for damage to the product, for property damage or for personal injury due to careless or improper usage of the product, or usage of the product which does not meet the manufacturer's specified purpose.  
**Note:** The measured values displayed on the scale are for reference only. For accurate measured values, please consult your physician or a qualified consultant or medical personnel.

CONTROLS AND COMPONENTS

- A. Display
- B. Platform
- C. Electrode
- D.  button
- E. **SET** button
- F.  button
- G. Battery compartment
- H. **UNIT** button (kg/lb)
- I. Anti-slip feet

PICTOGRAM IN DISPLAY

|                                                                                       |                          |                                                                                       |                                                                                               |
|---------------------------------------------------------------------------------------|--------------------------|---------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------|
|  | Age                      |  | Body fat %                                                                                    |
|  | Male                     |  | Hydration %                                                                                   |
|  | Female                   |  | Bone %                                                                                        |
|  | User profile P-01 – P-12 |  | Muscle %                                                                                      |
|  | Body fat indicator       |  | Weight unit                                                                                   |
|  | Calorie                  |  | Height unit                                                                                   |
|  | BMI parameter            |  | Percentage indicator when displaying body fat and water content, body bone and muscle content |

ERROR CODE

|                                                                                       |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                           |
|---------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|  | Overload                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                  |
|  | Low power                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                 |
|  | Wrong measuring<br> Step off and back on the scale, to measure again.<br> Another person should step on the scale to check the measuring.<br> If  &  still not working, it means sensor damaged. |
|  | Body fat less than 4%                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     |
|  | Body fat more than 60%                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                    |

BEFORE FIRST TIME USE

Remove the product and accessories carefully from the original packaging. Inspect the delivery contents for completeness and damages. Should the delivery contents be incomplete or damaged, contact your sales outlet immediately. Remove the battery protection before first use.

INSERT OR REPLACE BATTERIES

Fig.  Open the battery compartment. Insert 1 x 3 V battery (type CR2032) in the battery compartment. Make sure the polarity is correct as indicated. Close the battery compartment.  
**Note:**  
The scale has a battery low indicator. If the battery is low, than the display will show **Lo**. Replace the exhausted battery with a new battery of the same type.

#### ADJUST UNIT

Fig.  When the scale is on, press the **UNIT** button to switch the unit kg/lb.

#### OPERATION

##### Measuring body weight

Fig.  Always place the scale down on a solid, flat horizontal surface. Touch the platform and the scale will turn on. The display will show 0.0 kg. Step onto the platform and the scale starts to measure. The weight will be displayed. Do not move, moving will lead to discrepancies.
**Note:** The scale switches off automatically and turns on when steps on.

##### Set User profile

Fig.  Press **SET** and  /  buttons to choose an user profile P-01 to P-12. Confirm with **SET** and adjust the setting with  /  buttons Choose and adjust user profile P-01 to P-12 > Set male  or female  > set body height (100-250 cm) > Set age (10-100).

##### Body weight, body fat, muscles, bones and water content

Fig.  Press **SET** and  /  buttons to choose an user profile. Wait for 6 seconds or tap on the scale platform, 0.0 will be displayed. Step barefoot onto the platform and do not move. Ensure standing on the 2 electrodes. The measurement starts and the following values will be shown on the display:

**Body weight > Percentage of body fat content & fat status indication > Percentage of Water content > Percentage of muscles > Percentage of bones > calorie (KCAL) > BMI**

##### Auto User Recognition

Fig.  When the scale is in off- or standby mode, please step on the scale to get a weighing value. The display will lock the reading for one second. If the scale can find a similar pre-stored weighing value (difference less than ±2 kg), the display will depict the value, flashing once, and show the user profile. After automatically recognizing the current user, the profile will flash once and will then move to measuring other data, like body fat, hydration, etc.

##### Note:

- After searching users with similar weighing value and flashing the weighing value once, the scale cannot enter into body fat measuring if the user does not properly touch the metal electrodes. The scale will show the gender and current weighing value of the user before automatically turning itself off.
- If the scale finds two or more similar pre-stored weighing values, it will select the user profile with the closest weighing value. If the scale finds the same weighing values twice, it will select the user profile with the lower number.

##### Recall Function

- Press **SET** to switch on.
- Parameters of the last user profile or default parameters are shown.
- Press  /  to select user profile.
- Wait for 6 seconds or tap on the scale platform. 0.0 is shown in display.
- Press  to enter recall mode.
- Display the last test record of the selected user profile. Press  for quick page turning.
- The display loop of showing values repeats 3 times and then turns off automatically.

—The scale switches off automatically after these values have been displayed:—
In the following table you can find reference value of the scale.

| Age (years) | Female, values in % |               | Male, values in % |               | Status             | Indicator  |
|-------------|---------------------|---------------|-------------------|---------------|--------------------|------------|
|             | Body fat            | Water content | Body fat          | Water content |                    |            |
| ≤30         | 4.0-16.0            | 66.0-57.8     | 4.0-11.0          | 66.0-61.2     | severe underweight | <b>i--</b> |
|             | 16.1-20.5           | 57.7-54.7     | 11.1-15.5         | 61.1-58.1     | underweight        | <b>i-</b>  |
|             | 20.6-25.0           | 54.6-51.6     | 15.6-20.0         | 58.0-55.0     | normal             | <b>i</b>   |
|             | 25.1-30.5           | 51.5-47.8     | 20.1-24.5         | 54.9-51.9     | fat                | <b>i+</b>  |
|             | 30.6-60.0           | 47.7-27.5     | 24.6-60.0         | 51.8-27.5     | obese              | <b>i++</b> |
| >30         | 4.0-20.0            | 66.0-55.0     | 4.0-15.0          | 66.0-58.4     | severe underweight | <b>i--</b> |
|             | 20.1-25.0           | 54.9-51.6     | 15.1-19.5         | 58.3-55.3     | underweight        | <b>i-</b>  |
|             | 25.1-30.0           | 51.5-48.1     | 19.6-24.0         | 55.2-52.3     | normal             | <b>i</b>   |
|             | 30.1-35.0           | 48.0-44.7     | 24.1-28.5         | 52.2-49.2     | fat                | <b>i+</b>  |
|             | 35.1-60.0           | 44.6-27.5     | 28.6-60.0         | 49.1-27.5     | obese              | <b>i++</b> |

| Age     | Female   |        | Male     |        |
|---------|----------|--------|----------|--------|
|         | Muscle % | Bone % | Muscle % | Bone % |
| 16 – 50 | >31      | >8     | >38      | >10    |

| BMI parameter |           |            |     |
|---------------|-----------|------------|-----|
| Low           | Normal    | Overweight | Fat |
| <18,5         | 18.5 - 25 | 25 - 30    | >30 |

| Age          | Body height (cm) / Body weight (kg) |                   | Reference for energy intake in kcal/day |        |               |        |               |        |
|--------------|-------------------------------------|-------------------|-----------------------------------------|--------|---------------|--------|---------------|--------|
|              |                                     |                   | PAL-Value 1.4                           |        | PAL-Value 1.6 |        | PAL-Value 1.8 |        |
|              | male                                | female            | male                                    | female | male          | female | male          | female |
| <b>15-18</b> | 178.2 cm/ 69.2 kg                   | 165.5 cm/ 59.5 kg | 2600                                    | 2000   | 3000          | 2300   | 3400          | 2600   |
| <b>19-24</b> | 179.4 cm/ 70.8 kg                   | 165.8 cm/ 60.5 kg | 2400                                    | 1900   | 2800          | 2200   | 3100          | 2500   |
| <b>25-50</b> | 179.2 cm/ 70.7 kg                   | 165.1 cm/ 60.0 kg | 2300                                    | 1800   | 2700          | 2100   | 3000          | 2400   |
| <b>51-64</b> | 176.7 cm/ 68.7 kg                   | 162.6 cm/ 58.2 kg | 2200                                    | 1700   | 2500          | 2000   | 2800          | 2200   |
| <b>≥ 65</b>  | 174.2 cm/ 66.8 kg                   | 161.1 cm/ 57.1 kg | 2100                                    | 1700   | 2500          | 1900   | 2800          | 2100   |

PAL = physical activity level
Source: https://www.dge.de/wissenschaft/referenzwerte/energie/

PAL-Value 1.4: no/little physical activity
PAL-Value 1.6: medium physical activity

PAL-Value 1.8: heavy physical activity

#### CLEANING AND CARE

- Warning!** Never immerse electrical parts of the product in water or any other fluid during cleaning or operation. Never hold the product under running water.
- Warning!** When cleaning, never use solvents or abrasive materials, hard brushes, metallic or sharp objects. Solvents are harmful to human health and can attack plastic parts, while abrasive cleaning mechanisms and tools may scratch the surface(s).
- Clean the product's exterior surfaces with a slightly damp cloth and then dry them thoroughly.
- Make sure that no water enters the product interior (for example, through the display).
- Keep the product in a dry, cool location which is protected against moisture and direct sunlight. Keep this product out of the reach of children.

#### TECHNICAL DATA

Battery operated : 1 x 3 V  (Type CR2032) included

Loading capacity : Min. 2 kg - max. 150 kg

Devision (weight) : 100 g

#### TROUBLESHOOTING

If there are any problems during operation, please refer to the following table. If the problems cannot be solved as indicated under **Solution**, please contact your manufacturer.

| Problem                    | Possible causes                                    | Solution                                                     |
|----------------------------|----------------------------------------------------|--------------------------------------------------------------|
| Product does not work      | Battery fully discharged                           | Replace battery                                              |
| Incorrect measuring result | Scale isn't placed on even surface or not levelled | Reposition scale to even surface and ensure proper levelling |

#### DISPOSAL

 Do not dispose of this product as municipal waste. Return it to a designated collection point for the recycling of electrical and electronic equipment. Contact your retailer or local authorities for more information.

 Please dispose of batteries in an environmentally friendly manner Do not put batteries into the household waste. Please use the return and collection systems in your community or contact the dealer where you purchased the product.

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