
medisana®

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
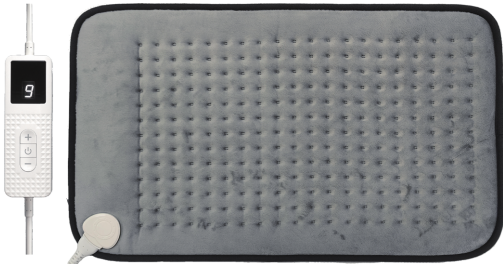
medisana®


GB Instruction manual
Heating pad HP 700

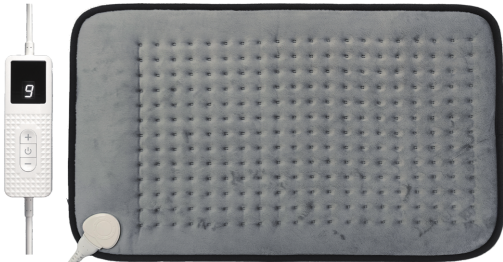
Explanation of symbols

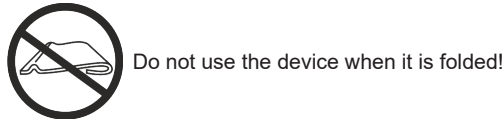
Do not use the device when it is folded!

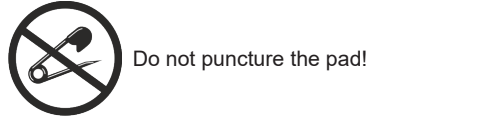
Do not puncture the pad!

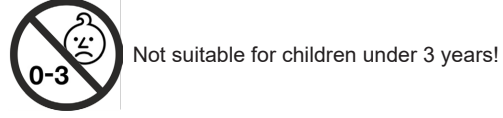
Not suitable for children under 3 years!

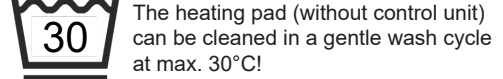
The heating pad (without control unit) can be cleaned in a gentle wash cycle at max. 30°C!

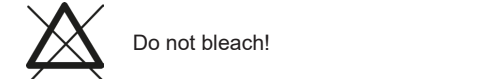
Do not bleach!

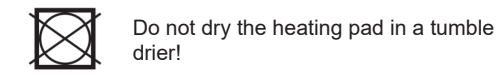
Do not dry the heating pad in a tumble drier!

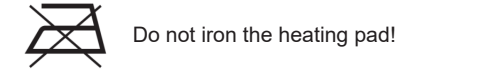
Do not iron the heating pad!

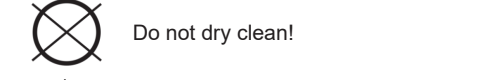
Do not dry clean!

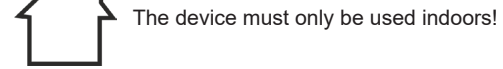
The device must only be used indoors!

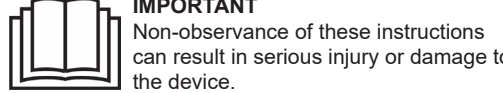
IMPORTANT
Non-observance of these instructions can result in serious injury or damage to the device.

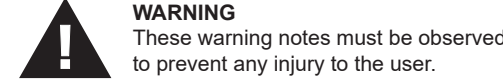
WARNING
These warning notes must be observed to prevent any injury to the user.

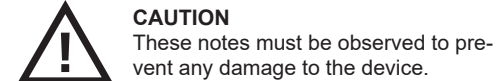
CAUTION
These notes must be observed to prevent any damage to the device.

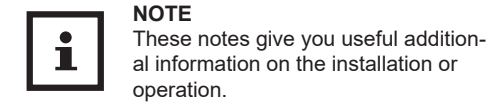
NOTE
These notes give you useful additional information on the installation or operation.

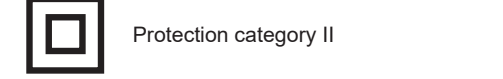Protection category II

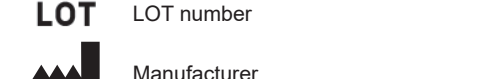LOT number

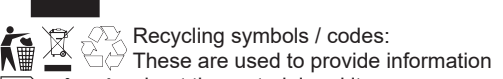Manufacturer

Recycling symbols / codes:
These are used to provide information about the material and its proper use and recycling.

GB

Important Instructions Retain for future use

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.

- **When this device is likely being used for an extended period of time the control unit must be set to level 1.**
- **This device is not intended for medical to be used for medical purposes in hospitals.**
- **This appliance must not be used by young children over the age of 3 years unless they have been given sufficient instruction on how to operate the controls safely.**
- **This device can be used by young children over the age of 3 years and younger than 8 years under supervision if the control unit is set to the lowest temperature level.**

- The device must only be used for its intended purpose as described in the instruction manual. Using the appliance for any other purpose invalidates the warranty.
- The device is not to be used for commercial use or medical purposes.

Safety Instructions – About the power supply

- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- The control switch should not be covered or placed on or below the pad while the unit is in operation.
- Keep the mains lead away from hot surfaces.
- Never carry, pull or turn the heating pad by the mains lead and do not let the lead get tangled.
- Never touch a device that has fallen into water. Unplug the unit from the mains outlet immediately.
- The controls and leads must not be exposed to any kind of moisture.
- The heating pad may only be operated with the control unit (SS43A) supplied.

Safety Instructions – Specific groups of people

- This appliance must not be used by persons insensitive to heat and other very vulnerable persons who are unable to react to overheating.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Children below the age of 3 years are not allowed to use this device, because they are not able to react on overheating.
- Should you have any concerns about health issues, consult your doctor before using the heating pad.
- Do not use the heating pad on parts of the body that are swollen, injured or inflamed.
- Should you experience persistent pains in the muscles and joints, please inform your doctor. Persistent pains can be a symptom of serious illness.
- If you wear a pacemaker, please be aware that the electromagnetic fields emanating from this electric product may, under certain circumstances, interfere with the functions of your pacemaker. Therefore, ask your doctor and the manufacturer of the pacemaker for advice before using this product.
- If you experience any pain or discomfort while using the unit, stop using it immediately.
- The cable and control unit of the appliance can give rise to risks of entanglement, strangulation, tripping or treading if not correctly arranged. The user shall make sure that excess ties and electric cords shall be arranged in a safe way.

Safety Instructions – Operating the device

- Carefully check the heating pad before each use, if you notice wear or damage.
- Do not use if you notice wear, damage or signs of improper use on the heating pad, the switch or the cables or if the device is not working.
- Do not fold or crease the heating pad during use.
- Neither pierce the pad with, nor attach safety pins or other spiky or sharp objects to the pad.
- Never use the pad unsupervised.
- Do not fall asleep while the heating pad is switched on. Prolonged use of the heated cushion at a high setting could lead to skin burns.
- Do not cover the pad with any other cushions.
- Do not sit on the heating pad. Instead, rest the pad on or against the part of the body to be treated. The intended use of the is to treat only back and neck.
- Do not use the heating pad when wet and only use it in dry environments (not in the bathroom or similar environments).

Safety Instructions – Maintenance and cleaning

- You are only permitted to clean the unit. To avoid hazards, never repair the device yourself. Contact the customer service.
- Cleaning and user maintenance shall not be made by children without supervision.
- If the power supply cable is damaged, it may only be replaced by **medisana**, an authorised dealer or by qualified personnel in order to avoid hazards.
- A damaged power supply cable must be replaced with a power cable of the same type.
- If a fault occurs, do not attempt to repair the unit yourself. Repairs must only be carried out by an authorized specialist dealer or other suitably qualified personnel.
- When storing the appliance, allow it to cool down before folding.
- Do not crease the appliance by placing items on top of it during storage.

Scope of delivery

Please check first whether the unit is complete and does not show any damage.

If in doubt, do not operate the unit and contact your dealer or service centre.

The scope of delivery includes:

- 1 **medisana** Heating pad **HP 700** with removable control unit SS43A
- 1 Quick Start Guide (short instructions)

If you notice any transport damage when unpacking, please contact your dealer immediately.

Packaging is recyclable or can be returned to the raw material cycle. Please dispose of packaging material that is no longer required properly.

WARNING
Please ensure that the packaging material is kept away from the reach of children! Risk of suffocation!

Control Unit SS43A

- ① LED display of the heating level (available: levels 1...9)
- ② Button to increase heating level
- ③ On/Off button
- ④ Button to reduce heating level

Using the device

Make sure, the connection cable of the control unit is properly connected to the plug of the heating pad. Plug the unit into the mains outlet and press the On/Off button ③. The device starts at lowest level (level 1), which is displayed at ①.

After a few minutes you will be able to feel the pad getting noticeably warmer. If you require a higher temperature setting, press the “+” button ② several times until you reach the required level. (The maximum level is level 9.) If the heating pad gets too hot for you, press the “-” button ④ several times until you reach the required level. After a short time, the heating pad will cool down.

The heating pad may cool down after a short period of time. This is normal. The fast re-heating according to the setting will follow immediately.

To switch the heating pad off, press the On/Off button ③. The LED display ① will go out, indicating that the unit is switched off. Unplug the unit from the mains outlet.

The heating pad will switch off automatically after approx. 90 minutes of continuous use. Unplug the unit from the mains outlet once you have finished using it.

Troubleshooting

If there is no function at all after switching on the heating pad (no heating level is displayed at ①), please check if the mains plug is correctly plugged in the mains outlet. If this is the case, the heating pad is defect. Then please do not use the heating pad and contact the service center.

Cleaning and maintenance

- Before cleaning the pad, unplug it from the mains outlet and let it cool down for at least ten minutes.
- The heating pad is equipped with a removable control unit (**SS43A**). Unplug the control unit from the heating pad by pulling out the connection cable from the plug on the heating pad.
- You may wash the heating pad gently by washing machine.
- However, the best way is to put the heating pad in a bath with lukewarm water and some mild detergent and then gently squeeze it out. Rinse the heating pad several times to remove all the detergent.
- Never use strong detergents or hard brushes.
- To dry, lay the pad flat on a moisture-absorbent towel or mat.
- Only use the blanket again once it is completely dry.
- Reconnect the control unit to the heating pad in the correct manner.
- You should untangle the cable from time to time.
- Store the unit in a clean and dry place, laid out flat with nothing covering or resting on top of it.

Technical Data

Name and model: **medisana** Heating Pad **HP 700** (with Control Unit SS43A)
Power supply: via Control Unit SS43A – input: 220–240V~ 50/60 Hz
Heat output: approx. 100 W (at max. heating level)
Heating levels: 9 levels.

Max. operating time: approx. 90 minutes, with automatic switch-off
Dimensions (L×W×H): approx. 60 × 35 cm
Weight: approx. 0.39 kg
Operating conditions: Use only in dry indoor areas
Storage conditions: Laid out flat and dry
Item number: 61154
EAN number: 40 15588 61154 4

Disposal

 This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner.
 Consult your municipal authority or your dealer for information about disposal.

Warranty and repair terms

Your statutory warranty rights are not restricted by our guarantee below. Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

1. The warranty period for **medisana** products is 3 years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

medisana GmbH, Carl-Schurz-Str. 2, 41460 NEUSS, GERMANY

In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.

The current version of this instruction manual can be found under <https://docs.medisana.com/61154>.
Information about service can be found here: www.medisana.com/servicepartners

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