# Shark CryoGlow

BRUGERVEJLEDNING  
BENUTZERANLEITUNG  
USER GUIDE  
GUÍA DEL USUARIO  
MANUEL D'UTILISATION  
GUIDA UTENTE

HANDLEIDING  
BRUKERVEILEDNING  
INSTRUKCJA OBSŁUGI  
KÄYTTÖOPAS  
ANVÄNDARGUIDE

**Intended use:**

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
## TABLE OF CONTENTS

GET TO KNOW YOUR SHARK CRYOGLOW MASK .....	38
WHAT'S IN THE BOX .....	38
CONTRAINDICATIONS .....	39
WARNINGS & PRECAUTIONS .....	41
LED PARAMETERS .....	43
GUIDANCE & MANUFACTURER'S DECLARATION - EMC .....	43
GETTING STARTED .....	45
FITTING THE MASK .....	46
GET TO KNOW THE REMOTE .....	47
USING THE REMOTE .....	48
LED SKIN ROUTINES .....	49
UNDER-EYE COOLING .....	50
ERROR MESSAGES .....	51
SYMBOLS .....	52
WARRANTY .....	53

SCAN QR CODE or visit [sharkclean.eu](http://sharkclean.eu)  
to access the full User Guide.
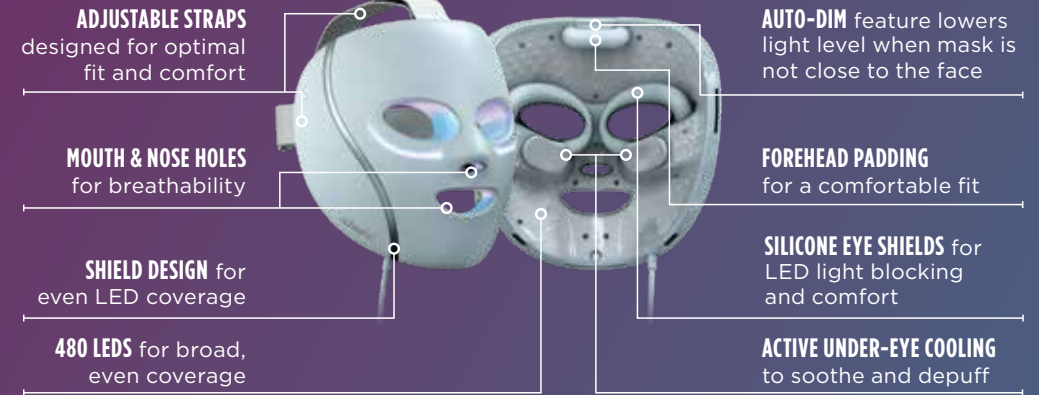
# GET TO KNOW YOUR SHARK CRYOGLOW MASK

## WHAT'S IN THE BOX

**Shark CryoGlow Mask + Straps**  
Pre-Assembled

-  Remote Control with LCD Screen + Control Buttons
-  Quick Start Guide + User Guide
-  Universal USB-C Charging Cable + Charging Plug
-  Optional clip-on chill pads
-  Storage Bag\*

## CONTRAINDICATIONS

### Do not use the Shark CryoGlow Mask if you are taking:

- Medication or supplements that cause sensitivity to light (photosensitivity), for example, medications for severe inflammatory acne (i.e. isotretinoin or doxycycline).

### Immediately stop using the device if you experience any discomfort and consult a medical professional.

### Do not use the Shark CryoGlow Mask if you have any of the following medical conditions:

- Retinal and/or macular disease, ocular albinism, or any eye condition that causes sensitivity to light
- Pregnant, nursing, or maybe pregnant
- Pacemaker or other electronic implanted device
- Implants in the face
- Lack of sensation (e.g., numbness)
- Light sensitive skin disorder (e.g., lupus, porphyria)
- Extreme light sensitivity of the eyes (i.e., photophobia)
- Light-induced headaches
- Albinism
- Melasma (i.e., condition that causes patches and spots)
- Sensitivity or allergy to plastic, aluminum, or self-fastening material

### Do not use the Shark CryoGlow Mask if you have any of the following conditions on the face, hairline, ears, or neck:

- Recent facial surgery (until you have been cleared by your doctor)
- Suspicious lesions or skin cancer
- Eczema
- Psoriasis
- Irregularly shaped, new, or changing brown spots
- Birthmarks
- Broken skin, skin rash, or blisters
- Tissue inflammation
- Skin infection, cold sore, or warts
- Active current outbreak of herpes simplex virus (cold sores)
- History of herpes outbreaks without receiving preventative treatment
- Facial treatment where the skin has not fully healed, including neurotoxin, dermal filler, micro needling, laser, or chemical peel

## CONTRAINDICATIONS, CONTINUED

### Do not apply the Shark CryoGlow Mask's cooling treatments on the eyeball or upper eyelid.

### Do not use the Shark CryoGlow Mask's cooling treatments if you have any of the following medical conditions:

- Cold hypersensitivity
- Cold urticaria (i.e., hives brought on by cold temperatures)
- Circulatory insufficiency

### PRECAUTIONS: Consult a medical professional before using the Shark CryoGlow Mask if:

- You have bruising on the face.
- You are experiencing skin sensitivity or irritation with retinol or benzoyl peroxide.

## WARNINGS & PRECAUTIONS

To reduce the risk of fire, thermal and chemical burns, electric shock, serious injury, or death, follow the precautions in this section and the procedures in this manual. The remote control contains a rechargeable lithium-ion battery. If damaged, the battery can overheat, catch fire, or leak a chemical that can burn skin and eyes.

### SUPERVISION AND INTENDED USERS:

- This equipment is not intended for use by persons with reduced physical, sensory, or mental capabilities, unless they have been given supervision or instruction concerning use of the equipment. Please refer to the contraindications before use.
- Do not use if any of the conditions in “Contraindications” (pages 39-40) apply to you.
- Do not leave the equipment unattended when it is powered on.
- This equipment is not intended for use by children. Children should be supervised to ensure they do not play with the equipment.
  - To reduce the risk of entanglement and potential strangulation, keep cords out of reach of children.
  - Small parts can be inhaled or swallowed by children—if any small parts break off or detach from the product, dispose of safely.
- This equipment is only intended to be used by one person. This person is the intended operator of the equipment.

### KEEP DRY.

- Keep the remote control dry. If it gets wet, do not charge or use—the battery may be damaged. Follow disposal instructions.
- Water and electricity are a dangerous combination!
  - Do not use this equipment in wet surroundings (e.g., in a bathtub, shower, or swimming pool).
  - Do not let water run into the device.
  - It is very important to keep the device dry. Any contact with water could damage the device. Do not use or store the device in the bathroom or near water.
  - Do not charge or use the device if it has been left outdoors.

### CHARGE SAFELY.

- Follow the charger manufacturer’s instructions.
- Only charge the equipment with the provided USB-C cable and charging brick: Dongguan Yinli Electronics, model YLJXA-T050200. Do not use any other USB-C cable or charging brick.
- Check if the voltage indicated on the equipment corresponds to the local mains voltage before you connect the equipment to avoid risk of electrocution or permanent damage to the equipment.
- If the wall socket used to power the equipment has poor connections, the plug will become hot. Make sure you plug the equipment into a properly installed wall socket.
- Do not charge the device in the bathroom or near water.
- Do not charge the remote after contact with water.
- Always position the equipment so the charger is easy to reach while charging.
- Store USB-C cable properly after charging.

### PROTECT THE REMOTE’S BATTERY FROM DAMAGE.

- Do not try to open the remote or replace the battery.
- Do not crush, drop, or damage the remote.
- Stop use or disconnect charging cable if the remote feels unusually hot or emits smoke or a bad smell.
- If liquid leaks from the remote, avoid contact. If liquid touches skin or eyes, wash with large amounts of water and seek medical advice.
- Follow disposal instructions.

## WARNINGS & PRECAUTIONS, CONTINUED

### USE SAFELY.

- The mask limits your vision. Do not wear while doing anything (such as driving) where limited vision poses a risk.
- CAUTION: The light emitted from this device may be harmful to the eyes. Infrared emitted from this device may cause eye irritation. Do not stare at the light source.
- Do not use the device when it has been left outdoors.
- If you are travelling by air, the mask should be carried in carry-on baggage and not checked in. The device must be completely switched off.
- Do not use the mask outdoors. Only use mask indoors under the following environmental conditions:
  - Temperature: between (10°C) and (30°C)
  - If the equipment is showing any signs of degradation, discontinue use immediately.
  - A relative humidity range of 15% to 90%, non-condensing, but not requiring a water vapor partial pressure greater than 50 hPa.
  - An atmospheric pressure range of 700 hPa to 1060 hPa.
- Do not store or charge the mask face down.
- Skin dryness, irritation, or redness may occur when using an LED routine. Apply a facial moisturizer if dryness occurs.
- To avoid the risk of temporary dry skin or hyperpigmentation, do not run back-to-back LED routines.

### CHECK THE PRODUCT FOR DAMAGE BEFORE EACH USE.

- Inspect the inside of the mask and strap for cleanliness, cracks, or sharp edges.
- Check the chill pads for signs of corrosion.
- These conditions can irritate or damage your skin.

### CLEAN PROPERLY.

- Clean the mask with a damp cloth.
- Avoid using chemical or alcohol based cleaners.
- Do if the equipment is showing any signs of degradation, discontinue use immediately.

### DO NOT REPAIR OR MODIFY.

- Do not attempt to repair or modify the mask or remote.
- Do not crush or dispose of in fire.

### DISPOSAL

- Dispose of the device according to local laws and requirements ensuring to consider lithium-ion battery disposal.

### POTENTIAL ELECTROMAGNETIC INTERFERENCE.

- WARNING: Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.
- WARNING: Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.
- WARNING: Portable radio frequency communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30cm to any part of the mask system, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.
- WARNING: Do not leave mask exposed in direct sunlight.

## TRANSPORT AND STORAGE CONDITIONS

Only transport and store the device in the following environmental conditions:

- Keep between -20°C and 45°C. For storage longer than 3 months, keep between 0°C and 25°C.
- Keep between 5°C and 32°C at a relative humidity up to 85%, non-condensing.
- Keep between 35°C and 60°C at a water vapour pressure up to 50 hPa.
- Atmospheric pressure range: 700 hPa to 1060 hPa.

## LED PARAMETERS

Mode	Red + Infrared Light Therapy	Blended Blue Light Therapy	Skin Sustain
Power (mW/cm <sup>2</sup> )	Red 64 ± 5mW/cm <sup>2</sup> , IR 64 ± 5 mW/cm <sup>2</sup> <b>Total mode</b> = 128 mW/cm <sup>2</sup>	<b>Step 1:</b> Blue 64 ± 8mW/cm <sup>2</sup> , IR 64 ± 5 mW/cm <sup>2</sup> <b>Step 2:</b> Blue 55 ± 6mW/cm <sup>2</sup> , Red 73 ± 5 mW/cm <sup>2</sup> <b>Step 3:</b> Red 73 ± 5mW/cm <sup>2</sup> , IR 55 ± 5 mW/cm <sup>2</sup>	<b>Total mode</b> = 128 mW/cm <sup>2</sup> <b>Total mode</b> = 128 mW/cm <sup>2</sup> <b>Total mode</b> = 128 mW/cm <sup>2</sup>

## GUIDANCE & MANUFACTURER'S DECLARATION - ELECTROMAGNETIC EMISSION

The device is Risk Group 1 as defined in IEC 62471.

The <b>Shark CryoGlow Mask</b> is intended for use in the electromagnetic environment specified below. The customer or the user of The <b>Shark CryoGlow Mask</b> should assure that it is used in such an environment.		
EMISSIONS TEST	COMPLIANCE	ELECTROMAGNETIC ENVIRONMENT - GUIDANCE
RF emissions	CISPR 11	Home healthcare environment
RF emissions	CISPR 11	Home healthcare environment
Harmonic emissions	IEC 61000-3-2	
Voltage function / flicker emissions	IEC 61000-3-3	

The <b>Shark CryoGlow Mask</b> is intended for use in the electromagnetic environment specified below. The customer or the user of The <b>Shark CryoGlow Mask</b> should assure that it is used in such an environment.			
IMMUNITY TEST	IEC 60601 TEST LEVEL	COMPLIANCE LEVEL	ELECTROMAGNETIC ENVIRONMENT - GUIDANCE
Electrostatic discharge (ESD) IEC 61000-4-2	Up to 8 kV contact Up to 15 kV air	Up to + 8 kV contact Up to ± 15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrostatic transient / burst IEC 61000-4-4	±2kV for main power lines	±2kV for main power lines	Mains power quality should be that of a typical commercial or home environment.
Surge IEC 61000-4-5	±0.5kV 1.5kV differential mode	±0.5 kV to 1.5kV differential mode	Mains power quality should be that of a typical commercial or home environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0% UT; 0.5 cycle <sup>(a)</sup> At 0°, 45°, 90°, 135°, 180°, 225°, 270°, & 315° <sup>(a)</sup>  0% UT; 1 cycle and 70% <sup>(a)</sup> ; UT; 25/30 cycles <sup>(b)</sup> Single phase: at 0°  0% UT; 250/300 cycle <sup>(b)</sup>	0% UT; 0.5 cycle <sup>(a)</sup> At 0°, 45°, 90°, 135°, 180°, 225°, 270°, & 315° <sup>(a)</sup>  0% UT; 1 cycle and 70% <sup>(a)</sup> ; UT; 25/30 cycles <sup>(b)</sup> Single phase: at 0°  0% UT; 250/300 cycle <sup>(b)</sup>	Mains power quality should be that of a typical commercial or home environment. If the user of the <b>Shark CryoGlow Mask</b> requires continued operation during power mains interruptions, it is recommended that the <b>Shark CryoGlow Mask</b> be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz magnetic field) IEC 61000-4-8	30 A/m at 50-60Hz	30 A/m at 50-60Hz	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or home environment.

NOTE U<sub>i</sub> is the a.c. mains voltage prior to application of the test level.

IMMUNITY TEST	IEC 60601 TEST LEVEL	COMPLIANCE LEVEL	ELECTROMAGNETIC ENVIRONMENT - GUIDANCE
Conducted RF IEC 61000-4-6	150KHz - 80MHz - Test level: 3V & 6V - Modulation: AM 1KHz 80%	- 150KHz - 80MHz - Test level: 3V & 6V - Modulation: AM 1KHz 80%	Portable and mobile RF communications equipment should be used no closer to any part of the <b>Shark CryoGlow Mask</b> , including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.  Recommended separation distance: $d = \left[ \frac{3.5}{V_1} \right] \sqrt{P}$ $d = \left[ \frac{12}{V_2} \right] \sqrt{P}$ $d = \left[ \frac{3.5}{E_1} \right] \sqrt{P}$ $d = \left[ \frac{7}{E_1} \right] \sqrt{P}$ 80 MHz to 800 MHz      800 MHz to 2.7 GHz
Radiated RF IEC 61000-4-3	80MHz - 2.7GHz - Test level: 10V/m - Modulation: AM 1KHz 80%	- 80MHz - 2.7GHz - Test level: 10V/m - Modulation: AM 1KHz 80%	where p is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m). <sup>b</sup>  Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey <sup>a</sup> . Should be less than the compliance level in each frequency range. <sup>b</sup> Interference may occur in the vicinity of equipment marked with the following symbol:  NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.  NOTE 2 These guidelines may not apply in all situations. Electromagnetic is affected by absorption and reflection from structures, objects and people.

**a)** The ISM (industrial, scientific and medical) bands between 150 kHz and 80 MHz are 6,765 MHz to 6,795 MHz; 13,553 MHz to 13,567 MHz; 26,957 MHz to 27,283 MHz; and 40,66 MHz to 40,70 MHz. The amateur radio bands between 0.15 MHz and 80 MHz are 1.8 MHz to 2.0 MHz, 3.5 MHz to 4.0 MHz, 5.3 MHz to 5.4 MHz, 7 MHz to 7.3 MHz, 10.1 MHz to 10.15 MHz, 14 MHz to 14.2 MHz, 18,07 MHz to 18,17 MHz, 21,0 MHz to 21.4 MHz, 24,89 MHz to 24,99 MHz, 28,0 MHz to 29.7 MHz and 50,0 MHz to 54,0 MHz.

**b)** Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the **Shark CryoGlow Mask** is used exceeds the applicable RF compliance level above, the TheraFace Mask should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the **Shark CryoGlow Mask**.

**c)** Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3V/m.

Recommended separation distances between portable and mobile RF communications equipment and the EQUIPMENT or SYSTEM - for EQUIPMENT and SYSTEMS

Recommended separation distances between portable and mobile RF communications equipment and the <b>Shark CryoGlow Mask</b>				
The <b>Shark CryoGlow Mask</b> is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the <b>Shark CryoGlow Mask</b> can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the <b>Shark CryoGlow Mask</b> as recommended below, according to the maximum output power of the communications equipment				
Rated maximum output of transmitter (W)	Separation distance according to frequency of transmitter (m)			
		150 kHz to 80 MHz outside ISM and amateur radio bands $d = \left[ \frac{3.5}{V_1} \right] \sqrt{P}$	150 kHz to 80 MHz outside ISM and amateur radio bands $d = \left[ \frac{12}{V_2} \right] \sqrt{P}$	80 MHz to 800 MHz $d = \left[ \frac{3.5}{E_1} \right] \sqrt{P}$
0.01	0.12	.020	0.035	0.07
0.1	0.38	0.63	0.11	0.22
1	1.2	2.00	0.35	0.70
10	3.8	6.32	1.10	2.21

For transmitters rated at a maximum output power not listed above the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

Please read the Warnings, Instructions, and Contraindications before first use.

### Charge fully before first use

Plug small end of USB-C cable into charging port on top of remote.

Plug large end of USB-C into power adapter. Plug power adapter into wall outlet.

Allow remote to charge fully (**this should take about 3 hours**). The remote screen will show your charging progress.

Once charged, unplug power adapter and disconnect USB-C cable from remote.

**NOTE:** Unplug before using face mask.  
Mask will not work when plugged into wall outlet.

### Powering ON/OFF

**ON:** To turn on your mask and remote, **long-press** the remote's dial/center button for 2 seconds.

**OFF:** To turn off, **long-press** the dial/center button.

**AUTOMATIC SHUT-OFF:** If the device is not manually turned off, the mask and remote will automatically shut off after completing a session.

The mask fans will continue to run for 60 seconds while the mask is powering down after a session.

**NOTE:** Demo video will run upon first use of the remote—please watch it for further fit and comfort instructions. You can also rewatch this video by heading to the remote's settings.

### Step 1 / Adjust the back strap

Start by adjusting the back strap (with the buckles). To tighten or loosen, **slide the buckles simultaneously** together or apart. This strap may cover the top of your ears—adjust until this feels comfortable.

Tighten

Loosen

### Step 2 / Adjust the top strap

To adjust the mask up or down on your face, loosen or tighten the self-fastening strap on the top of your head. Ensure the strap sits in the center of your head (following a center part).

### Step 3 / Adjust chill pad fit (optional)

The chill pads should sit under your eyes, on top of your cheekbones. If the chill pads have poor or inconsistent contact with your under-eye area, **add the optional clip-on chill pads to improve fit and comfort**

Remove the mask and take the optional clip-on chill pads out of the box. They are marked "L" and "R"—make sure the letters are pointing toward the chin of the mask. Align the optional clip-on chill pads with the original chill pads (that are permanently attached to the mask) and snap them in place. You'll feel a click when they are properly attached.

**NOTE:** LEDs work best when they are not in contact with your skin.  
This is why the mask sits at a slight distance from your face.

### Pro Tips:

- 1 The back strap should sit parallel to the floor, and your eyes should be centered in the eye holes.
- 2 Use one hand to hold the back strap in position and the other hand to adjust the top self-fastening strap until the mask is positioned comfortably and securely on your head.

Get comfortable! We recommend using the mask in a seated or reclined position. Take this time to relax, and enjoy your routine!

Please visit [sharkclean.eu](https://sharkclean.eu) for more helpful images and videos.

### Changing Brightness

1. Turn **ON** the mask using the remote.
2. Turn the dial to view **Settings** and click the dial button to select.
3. Select **Screen Brightness**, then use the dial to increase or decrease the brightness.
4. Click the dial button to make your selection.
5. After selecting a brightness level, you will be sent back to the main menu.

### Changing Language

1. Turn **ON** the mask using the remote.
2. Turn the dial to view **Settings** and click the dial button to select.
3. Select **Language** and use the dial to select the language of your choice.
4. After selecting a language, you will be sent back to the main menu.

### Reset Progress

1. Turn **ON** the mask using the remote.
2. Turn the dial to view **Settings** and click the dial button to select.
3. Select **Reset Progress**, then click the dial button to select **Confirm**.
4. After confirming, you will be sent back to the main menu.

### Factory Reset

1. Turn **ON** the mask using the remote.
2. Turn the dial to view **Settings** and click the dial button to select.
3. Select **Factory Reset**, then click the dial button to select **Confirm**.
4. After confirming, you will be sent back to the main menu.

**Pro Tip:** You can replay the demo video on the remote and access our support page through settings.

## LED SKIN ROUTINES

The mask offers **3 LED routines** to choose from, whether your goals are treating your skin or maintaining it.

### TREATMENT ROUTINES

The Shark CryoGlow LED mask emits energy in the red and infrared light spectrum and in the red, blue and infrared light spectrum for the treatment of mild-to-moderate inflammatory acne. Select either the **Red + Infrared Light Therapy or Blended Blue Light Therapy** routine below and use your selection **exclusively once daily for a minimum of 8 weeks**. For both the treatment of acne and to reduce redness and improve roughness, we recommend using the Blended Blue Light Therapy mode exclusively for 8 weeks. Do not stack modes.

#### Red + Infrared Light Therapy (6 MINUTES)

Reduce redness and improve roughness.\*

\*Results from 12-week clinical study with 59 subjects and 1x/day use in Red + Infrared Light Therapy treatment mode. Individual results may vary.

\*\*Results from 12-week clinical study with 60 subjects and 1x/day use in Blended Blue Light Therapy treatment mode. Individual results may vary.

#### Blended Blue Light Therapy (8 MINUTES)

Reduce acne for a clearer, smoother complexion.\*\*

### MAINTENANCE ROUTINE

Once you have achieved your target results with a treatment routine, use Skin Sustain **once daily** to help maintain your results.

#### Skin Sustain (4 MINUTES)

Designed to maintain glowing skin after an 8-week cycle.

**Did you loose your streak? Don't worry!**  
**Readjust your goals and continue with your skincare journey.**

No UV light is emitted from this mask.

## UNDER-EYE COOLING

Use cooling in the morning to help soothe your under-eyes and in the evening to relax.

### LED + COOLING

Under-eye cooling will automatically start with any LED routine and has the option to be turned off. Use the remote's chill button to cycle through 3 temperatures to quickly soothe your under-eye area.

### COOLING ONLY

#### Under-Eye Revive (5, 10 or 15 minutes)

An adjustable under-eye cooling **ONLY** routine (no LEDS) for a soothing and refreshing experience designed for optimal relaxation.

### FIND YOUR PERFECT ROUTINE

Scan the QR code to take our quiz to find the best routine for your skin concerns.

## ERROR MESSAGES

Error message	Solution
Remote Over Temperature. Let remote cool down before using.	Let remote cool down before using. The mask should not be used in temperatures higher than (30°C).
Remote Under Temperature. Let remote warm up before using.	Let remote warm up before using. The mask should not be used in temperatures lower than (0°C).
Remote error A3. Call Customer Service.	Call Customer Service.
Charger not recognized.	Use the provided Shark USB-C charging cable and wall plug to charge the mask.
Remove remote from charger.	The mask's battery is charged; remove it from the charger.
Low Battery. Charge before use.	The mask's battery is low; please charge it.
Battery error B1. Call Customer Service.	Call Customer Service.
Remote error B2. Charge the remote before next use.	Charge the remote before next use using the provided Shark USB-C charging cable and wall plug.
Charging error C1. Disconnect and reconnect to charge.	Disconnect and reconnect the remote to charge using the provided Shark USB-C charging cable and wall plug provided.
Fan error F1. Call Customer Service.	Call Customer Service.
Error A1. Call Customer Service.	Call Customer Service.
Error A2. Call Customer Service.	Call Customer Service.

## SYMBOLS

Symbol	Explanation
	Type BF Applied Part
	Class II Equipment
	Direct Current
	Refer to User Guide
<b>IP22</b>	Protection against ingress: IP2X: Protected against solid foreign objects of 12.5 mm Ø and greater IPX2: Protection against vertically falling water drops when ENCLOSURE tilted up to 15°
AS TO ELECTRICAL SHOCK, FIRE AND MECHANICAL HAZARDS ONLY IN ACCORDANCE WITH AAMI ES 60601-1:2005/A2:2021; CAN/CSA-C22.2 No. 60601-1:14(Reaffirmed 2022); IEC 60601-1-11:2015/AMD1:2020; IEC 60601-2-83:2019/AMD1:2022	
	Manufacturer
	This Device complies with the Medical Device Regulation (EU) 2017/745
	Authorized Representative In the European Union
	Medical Device

## TWO (2) YEAR LIMITED GUARANTEE

When a consumer buys a product in Europe, they get the benefit of legal rights relating to the quality of the product (your statutory rights). You can enforce these rights against your retailer. We give you an additional manufacturer's guarantee of two years. These terms and conditions relate to our manufacturer's guarantee only – your statutory rights are unaffected.

The conditions below describe the prerequisites and scope of our guarantee. They do not affect your statutory rights or the obligations of your retailer and your contract with them.

### Shark Guarantees

Every Shark product comes with a free parts and labour guarantee. You'll also find online support at [sharkclean.eu](http://sharkclean.eu)

### How do I register my Shark guarantee?

You can register your guarantee online within 28 days of purchase. To save time, you'll need the following information about your product.

- Date you purchased the unit (receipt or delivery note).
- To register online, please visit [sharkclean.eu](http://sharkclean.eu)

### IMPORTANT:

- The guarantee will only cover your product from the date of purchase.
- The free Shark guarantee is only valid in the country where the product was purchased.
- Please keep your receipt at all times. Should you need to use your guarantee we will need your receipt to verify the information you have supplied to us is correct. The inability to produce a valid receipt may invalidate your guarantee.

### What are the benefits of registering my free Shark guarantee?

When you register your guarantee we'll have your details to hand if we ever need to get in touch. You can also receive tips and advice on how to get the best out of your Shark product and hear the latest news about new Shark technology and launches. If

you register your guarantee online, you'll get instant confirmation that we've received your details.

### What is covered by the free Shark guarantee?

Repair or replacement of your Shark product (at Shark's discretion), including all parts and labour. A Shark guarantee is in addition to your legal rights as a consumer.

### What is not covered by the free Shark guarantee?

1. Normal wear and tear of wearable parts (such as accessories). Replacement parts are available for purchase at [sharkclean.eu](http://sharkclean.eu)
2. Damage caused by misuse, abuse, negligent handling, failure to perform required maintenance or damage due to mishandling in transit.
3. Damage caused by maintenance not authorised by Shark.

### Where can I buy genuine Shark spares and accessories?

You'll find a full range of Shark spares and replacement parts/accessories for all Shark machines at [sharkclean.eu](http://sharkclean.eu)

Please remember that damage caused by the use of non-Shark spares may not be covered under your guarantee.

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SHARK and SHARK CRYOGLOW are registered trademarks in the European Union of SharkNinja Operating LLC.
© 2025 SharkNinja Operating LLC. All rights reserved.
SharkNinja Operating LLC 89 A St #100, Needham, MA 02494, USA
Emergo Europe B.V. Westervoortsedijk 60, 6827 AT Arnhem, The Netherlands