CURRENTBODY SKIN™

LED LIGHT THERAPY MASK - SERIES 2

USER GUIDE AND MANUAL

CONTENTS

CurrentBody Skin LED Light Therapy Mask	04
What's in the box?	05
Important safety information	06
Warnings	07
Cautions	09
Contraindications	09
Device overview	11
How to use	12
Warranty	16
Specification	17
Symbols Explained	18
Safety and EMC	19
Manufacturer information	22
Français	23
Deutsch	37
Español	51
Italiano	65
繁體中文	79
简体中文	93
日本語	107
한국어	121

Welcome to CurrentBody Skin LED Light Therapy Mask

Developed through our expertise in beauty technology, CurrentBody Skin is a range of high-performance beauty tech for home use.

Every CurrentBody Skin product pairs the gold standard of clinical technology with unique design & advanced software to maximise your results.

Veritace® - Quality, Safety & Results in Beauty Tech

Choosing Veritace® guarantees that the technology in your hands is one of the best, safest and most reliable on the market.

We achieve this standard through our unique Veritace® product development process. This allows us to ensure that every product sold worldwide is produced to the highest standard.

The Veritace® Mark on the box is our guarantee of quality. Every product can be traced back through individual product testing to ensure it meets the Veritace® standard.

WHAT IS IN THE BOX?

System Includes:

CurrentBody Skin LED Light
Therapy Mask

Multi-strap x 3

USB-C Charging Cord

Storage Bag

User Manual

Eye Inserts

NFC card

IMPORTANT SAFETY INFORMATION

- DO NOT operate the device where it could come in contact with water.
- DO NOT use or charge the device with a power adapter that is damaged or faulty.
- DO NOT harshly bend or abuse the USB-C Charging Cord.
- DO NOT attempt to take the device apart. The device is not servicable by the user. This will void the warranty.
- STOP using the device immediately if you think it is defective.

CUSTOMER CARE

Please contact CurrentBody via the website if you have any questions or need any advice.

WARNINGS

- This device should not be used by children 12 years or younger. People between the ages of 13 and 18 must always use the device under adult supervision.
- Please consult your doctor prior to use if you are, or suspect you are pregnant or you are breastfeeding.
- Please consult your doctor prior to use if you have or suspect you have cancer.
- DO NOT use the device if it is wet. If the device has fallen into water, you must unplug it from any power supply before you touch it.
- DO NOT look directly into the light as this may cause temporary discomfort. Ensure that you place the device on the face so that you can see through the eye holes and the device is not covering the eyes.
- DO NOT use if your skin is red or swollen from a previous treatment.
- DO NOT use the device to treat skin that is painful to touch, raised (higher than the surrounding skin), crusting or bleeding. This may be a sign of a more serious condition and should be checked by your doctor.
- No modification of this equipment is allowed.
- This device contains a Lithium Ion Battery. Before disposal of this device, please review your local laws and requirements surrounding Lithium Ion Battery disposal. The preferred method of disposal is through recycling.

- DO NOT let children play with the device. Cleaning and charging the device must not be done by children.
- DO NOT use the device to treat moles, warts, open sores, cancerous lesions or any skin condition.
- DO NOT use the device on spots without pigmentation or areas on the skin that have vitiligo.
- DO NOT use the device if you have any open wounds.
- DO NOT use the device outdoors to prevent damage to the product.
- DO NOT drop or insert any object into any openings of the device.
- DO NOT use or store the device near hot surfaces or in direct sunlight.
- DO NOT use the device with non-manufacturer-approved accessories.
- DO NOT use solvents or abrasive cleaners on the device.
- DO NOT use the device in any other way than specified in this manual.
- DO NOT use the device if you have had a recent surgery or other surgical procedures.
- DO NOT use the device over a bacterial infection as this may result in the worsening of symptoms.
- Close supervision is necessary when this device is used, cleaned or stored near children or those with reduced physical or mental abilities.
- DO NOT use the device with any other USB-C cable than provided by CurrentBody.
- DO NOT scratch the light emitting surface.

CAUTIONS

- Discontinue use if the device overheats or becomes uncomfortably hot.
- Always unplug the device from the outlet after charging is complete.
- Do not store the device in direct sunlight or in a location with temperatures over 100 Fahrenheit as this may cause damage to the plastic coating.
- Use only the supplied Universal USB-C charging cord. Use of an incorrect power cord could damage the device.
- Do not attempt to open or repair the device or power supply.
- Do not drop the device on hard surfaces, as this could damage it. Using a damaged device could result in injury.
- If you have any concerns about the skin you are treating, you should cease use and check with your doctor.
- Contact your doctor if you have any medical questions regarding use of this device.

CONTRAINDICATIONS

Certain medical conditions, medications and chemicals are known to induce photosensitivity. Photosensitivity is a common side effect of various medications. These can include many commonly used medications, skincare, supplements, antibiotics, diuretics and steroids.

For an extensive list of common causes of photosensitivity, please see here:

<https://www.currentbody.com/pages/led-contraindications>

St John's Wort, certain perfumes, antiperspirant, antibacterial soaps, artificial sweeteners, naphthalene (mothballs), petroleum products, brightening agents found

in laundry detergent and cadmium sulphide (a chemical injected into the skin during tattooing).

Anyone with the following medical conditions should not use this device:

- Epilepsy or seizures triggered by light
- Metabolic disorders such as Porphyria
- Autoimmune disorders such as Systemic Lupus Erythematosus (SLE)
- Photosensitive disorders
- Active cancers
- Use of photosensitive medication or skincare
- Eye disease
- Glaucoma
- Cataracts
- Recent laser eye surgery
- Light induced headaches or migraines.

Test yourself for light sensitivity prior to using:

1. Ensure your skin is clean and dry.
2. Place your LED product in contact with the skin. Ideally inner thigh, arm or back of hand.
3. Turn on your device and complete a full treatment.
4. Wait 24 hours. If there are no reactions, carry out the treatment protocol indicated for your device.
5. If you experience any reactions, please contact our Customer Support team.

DEVICE OVERVIEW

The CurrentBody Skin LED Light Therapy Mask is powered by 3 highly precise LED wavelengths, proven in scientific testing to be most effective for anti-ageing and skin rejuvenation.


The flexible silicone mask molds to your face to give maximum skin coverage for the most consistent, effective results. The LED light penetrates the skin, energising the cells in the deepest layers where collagen is formed, reducing the appearance of fine lines and wrinkles and resulting in a brighter complexion.

HOW IT WORKS

The CurrentBody Skin LED Mask is powered by 236 Veritace® red (633nm), near-infrared (830nm) and deep near-infrared (1072nm) LEDs.

The device is programmed to deliver 10-minute treatments and with consistent use, 3-5 times per week can result in firmer, more radiant skin in just 4 weeks.

Light Wavelength Diagram

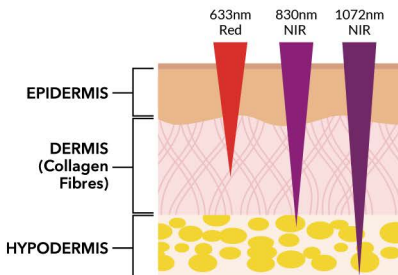HOW TO USE

Set up:

1. Unpack your device. The device may come with some charge, but we recommend always fully charging before first use.
2. Insert the USB-C power cord into the controller and then into a suitable USB charger (not supplied).
3. Insert the USB charger into a mains outlet and switch on.
4. The white power indicator lights will flash to indicate charging. The device will take 4.5 hours to charge fully. When the battery is fully charged the white light will stop flashing and turn solid white.
5. Ensure you unplug the device from the USB-C charger before using. For your safety, the device will not work when charging.
6. Place the device in the protective storage bag when not in use.
7. When the device needs recharging, the power indicator light will flash orange. To recharge the device, please follow the instructions above.

INSTRUCTIONS FOR USE:

You must conduct a sensitivity test before first use.

1. Wash and fully dry your face prior to use so your skin is clean and free from skincare products. Do not use any other skincare with the LED as it may cause irritation or block the light.
2. Attach the straps, place the mask to your face and adjust the straps around your head.
3. Press and hold the power button to turn on your device. Your treatment will automatically run for 10-minutes.
4. When your treatment has finished you can apply your skincare as normal.

FOR BEST RESULTS:

- Use 3-5 times weekly.
- DO NOT carry out more than one treatment on the same area within a 24 hour period.
- Use the device on clean, dry skin.
- Ensure skincare is applied after using the device as it can cause irritation and block light penetration, not before.
- Wipe device clean after each use.

MAINTENANCE:

- This device has no serviceable parts. The battery is not accessible and not serviceable. Users should not attempt to access or service the battery as this will damage the device and void the warranty.
- Device includes only an ON/OFF button. All settings are set to the factory setting.
- Charge the device only with the supplied USB-C power charging cord and a universal USB-C power adapter (not included).
- The power indicator light will begin to flash orange when battery levels are low. This indicates the device should be recharged.
- Compatible Power Supply: 100-240V
- The lifetime of the Lithium battery can be maximized by allowing it to fully charge and then fully discharge during regular use.

CLEANING:

1. Turn device off and unplug from any power source.
2. Wipe the tops of the LEDs with a damp or microfibre cloth to remove any residue after each use.
3. Clean the device body by using damp or microfibre cloth. Do not use chemical or alcohol based cleaners on the plastic coating.
4. Allow the device to dry completely prior to use.

TROUBLESHOOTING

System is not turning on:

- a. Is the device fully charged? (power indicator light is on and solid).
- b. Is the power indicator light illuminated?
- c. Is the device disconnected from the mains? The device will not work while charging.

System is not charging:

- a. Is USB-C power cord firmly inserted into the device? (does power indicator light turn on?).
- b. Try a different wall outlet to determine if there is power to the outlet.

System shuts off abruptly, prior to completion of 10-minute treatment:

- a. System is designed to shut off automatically after 10 minutes.
- b. System may shut off after multiple repeated uses if the device exceeds desired temperature. Allow the device to cool for several minutes and try again.

Some lights appear turned off:

The near- infrared (NIR) light spectrum is not visible to the human eye. These lights will look as if they are not on, but they are operating properly. If the visible lights are on, then the device is operating as designed.

The device does not switch off after a 10 minutes treatment:

The device is faulty. Please contact CurrentBody.

The LEDs flicker during use:

The device is faulty. Please contact CurrentBody.

WARRANTY

Under normal use this product is guaranteed for 2 years after its original purchase date against manufacturing defects in materials or workmanship. This warranty extends only to the original purchaser. If the product should become faulty within the warranty period, contact CurrentBody.

- CurrentBody reserves the right to repair a defective product.
- CurrentBody reserves the right to replace a defective product with the most compatible product currently available.

This warranty does not cover damage caused by the following:

- Accident, misuse, abuse or alteration
- Servicing by unauthorised persons
- Connection to incorrect power sources
- Use with unauthorised accessories

CurrentBody shall not be responsible for any incidental, special or consequential damages resulting from the use of this product. All implied warranties, including but not limited to implied warranties of fitness and merchantability are limited in duration to two years from date of original purchase. This warranty gives you specific legal rights, and you may also have other rights that exclusion or limitation of incidental or consequential damages, or limitations on how long an implied warranty lasts, so some of the above limitations and exclusions may not apply to you. This warranty is given in lieu of all other warranties, written or oral, whether expressed by affirmation, promise, description, drawing, model or sample. Any and all warranties other than this one, whether expressed or implied, including implied warranties of merchantability and fitness for a particular purpose are hereby disclaimed. If the product is defective within the warranty period please contact CurrentBody.

SPECIFICATION

Specification	
Produce name	LED Light Therapy Mask
Model	MK-90H
Power source	Adapter: Input: 100-240V/50-60Hz Output: 5.0V=1.0A
Battery capacity	2600mAh
Weight	1220g±10g
Dimensions	190 × 98 × 207 mm (L × W × H)
Wavelengths	Red 633nm + NIR 830nm + Deep NIR 1072nm
Energy Source	LEDs
Power Density	30mw/cm2
Number of LEDs	236 LED diodes: 110 × 633nm, 110 × 830nm, 16 × 1072nm
Treatment time	3-5 times weekly for 4 weeks
Location for Use	Face
Environment of Use	OTC
Safety and EMC	IEC 60601-1, IEC 60601-1-2, IEC 60601-1-6, IE C60601-1-11, IEC 60601-2-57
Biocompatibility	ISO 10993-1, ISO 10993-5, ISO 10993-10
Operating Environment	Temperature: -10~35 C Humidity: 30-80%RH
Storage Environment	Temperature: -20~60 C Humidity: 30-90%RH

THE SYMBOLS AND DESCRIPTION

No.	Symbols	Description
1		Manufacturer.
2		"WEEE (Waste Electrical and Electronic Equipment)". This means that the product including its battery must not be disposed of in household waste.
3		Consult instructions for use.

THE SYMBOLS AND DESCRIPTION

No.	Symbols	Description
4		The date of manufacture.
5	IP22	Protection against ingress of solids, it means the volume of the solids should less than 12mm. Protection against ingress of water, it means the de vice could protected against vertical falling water drop 15°.
6		Type BF Applied Part.
7	CE	This symbol means: Conforms to certain EC Directives on health, safety, and environmental protection standards. CE stands for 'Conformité Européenne'.
8	UK CA	This symbols implies conformity to UK Directives on health, safety and environmental protection standards. UKCA stands for UK Conformity Assessed.
9		This symbol means: Class B, Civil or household products using light sources
10		Batch code.
11		Symbol for "SERIALNUMBER". This symbol shall be accompanied by the manufacturer's serial number.
12		Temperature limit.
13		Humidity limit.
14		This symbol means keep dry.
15		This symbol means keep the de vice away from direct sunlight.

SAFETY AND EMC

- 1) It complies with Medical Electrical Safety Standards (IEC 60601-1, IEC 60601-1-11, IEC 60601-2-57).
- 2) This device also complies with Medical EMC Standard (IEC 60601-1-2).

Guidance and manufacturer's declaration - electromagnetic emissions		
LED Light Therapy Mask is intended for use in the electromagnetic environment specified below. The customer or the user of LED Light Therapy Mask should assure that it is used in such an environment.		
Emissions	Compliance	Electromagnetic environment-- guidance
RF emissions CISPR 11	Group 1	LED Light Therapy Mask uses LED energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	LED Light Therapy Mask is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	$P_{st} \leq 1.0$; $P_{lt} \leq 0.65$; $d(t) \text{ sec} \leq 0.5$; $d_{max} \leq 4\%$ $d_c \leq 3.3\%$	

Guidance and manufacturer's declaration - electromagnetic immunity			
LED Light Therapy Mask is intended for use in the electromagnetic environment specified below. The customer or the user of LED Light Therapy Mask should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment --guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 8 kV contact ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV air	± 8 kV contact ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrical fast transient/burst IEC 61000-4-4	± 2 kV 100 kHz repetition frequency	± 2 kV 100 kHz repetition frequency	Mains power quality should be that of a HOME HEALTHCARE ENVIRONMENT
Surge IEC 61000-4-5	$\pm 0,5$ kV, ± 1 kV Line-to-line	$\pm 0,5$ kV, ± 1 kV Line-to-line	Mains power quality should be that of a HOME HEALTHCARE ENVIRONMENT

Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % UT; 0,5 cycle (At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°) 0 % UT; 1 cycle and 70 % UT; 25/30 cycles (Single phase: at 0°) 0 % UT; 250/300 cycle	0 % UT; 0,5 cycle (At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°) 0 % UT; 1 cycle and 70 % UT; 25/30 cycles (Single phase: at 0°) 0 % UT; 250/300 cycle	Mains power quality should be a HOME HEALTHCARE ENVIRONMENT. If the user Of LED Light Therapy Mask requires continued operation during power mains interruptions, it is recommended that LED Light Therapy Mask be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a HOME HEALTHCARE ENVIRONMENT
NOTE UT is the a.c. mains voltage prior to application of the test level			

Guidance and manufacturer's declaration – electromagnetic immunity			
LED Light Therapy Mask is intended for use in the electromagnetic environment specified below. The customer or the user of LED Light Therapy Mask should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment --guidance
Conducted RF IEC 61000-4-6	3 Vrms 0,15 MHz and 80 MHz 6 Vrms in ISM and amateur radio bands between 0,15 MHz and 80 MHz	3 Vrms 0,15 MHz and 80 MHz 6 Vrms in ISM and amateur radio bands between 0,15 MHz and 80 MHz	Portable and mobile RF communications equipment should be used no closer to any part of LED Light Therapy Mask, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $d = 1,2 \sqrt{P}$ $d = 1,2 \sqrt{P} \quad 80 \text{ MHz to } 800 \text{ MHz}$ $d = 1,2 \sqrt{P} \quad 800 \text{ MHz to } 2,7 \text{ GHz}$ Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation Distance in meters (m).

Radiated RF IEC 61000-4-3	10 V/m 80 MHz to 2.7GHz 80 % AM at 1 kHz	10 V/m 80 MHz to 2.7GHz 80 % AM at 1 kHz	Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, should be less than the compliance level in each frequency range. b Interference may occur in the vicinity of Equipment marked with the following symbol: 
NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies. NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			
a. Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which LED Light Therapy Mask is used exceeds the applicable RF compliance level above, LED Light Therapy Mask should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating LED Light Therapy Mask. b. Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 10 V/m.			

Recommended separation distances between portable and mobile RF communications equipment and the LED Light Therapy Mask			
LED Light Therapy Mask is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of LED Light Therapy Mask can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and LED Light Therapy Mask as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter (W)	Separation distance according to frequency of transmitter (m)		
	150 kHz to 80 MHz $d = 1,2/\sqrt{P}$	80 MHz to 800 MHz $d = 1,2/\sqrt{P}$	800 MHz to 2.7 GHz $d = 2,3/\sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23
For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer. NOTE1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies. NOTE2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

MANUFACTURER INFORMATION

The Beauty Tech Group Ltd

Q17, Crossley Road, Stockport, Greater Manchester, SK4 5BB

EU Authorised Representative: Certification Company,
Veluwezoom 42, 1327 AH, Almere, The Netherlands

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CURRENTBODY SKIN™ LED Light Therapy Mask is manufactured for The Beauty Tech Group Ltd.