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











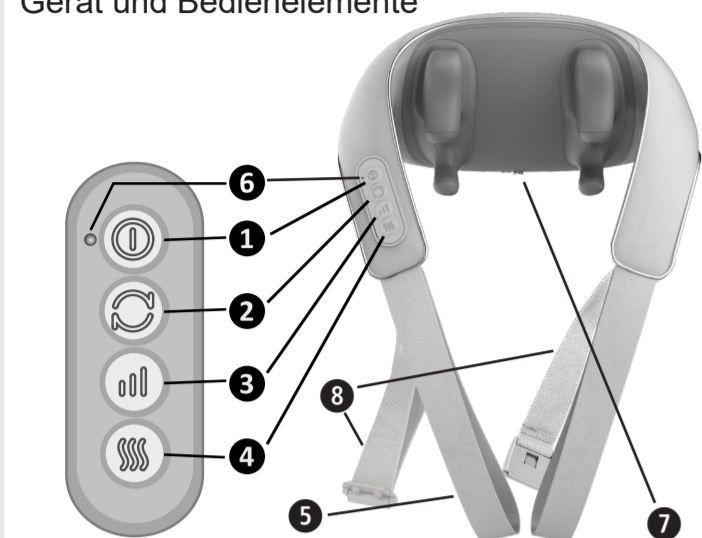
## GB Instruction Manual Cordless Neck Massager NM 990

### Device and controls

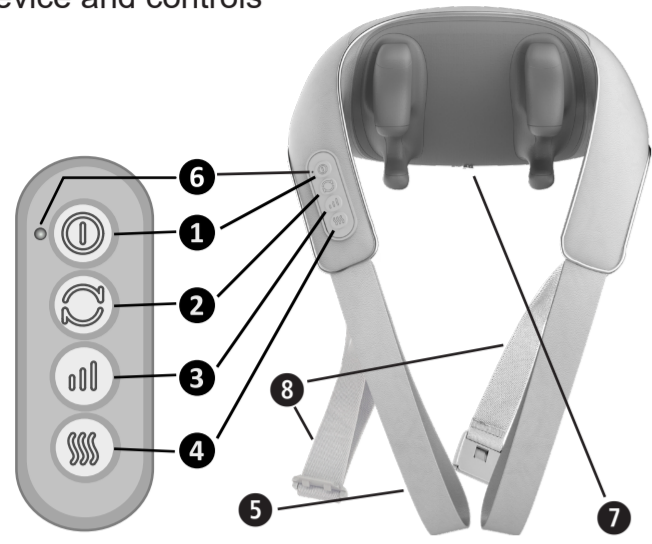- 1 On/Off button
- 2 Button to select between 2 massage directions
- 3 Button to select between 3 massage intensity levels
- 4 HEAT button, switches heat function on/off
- 5 Hand strap (each 1 on right and 1 on left side)
- 6 Light indicator (LED)
- 7 Removable neck cover
- 8 Back strap

### Explanation of symbols

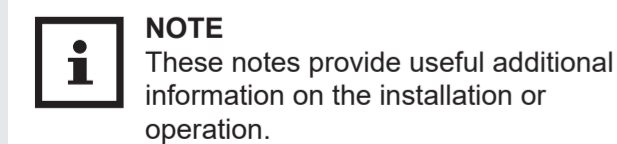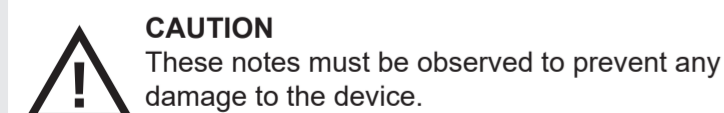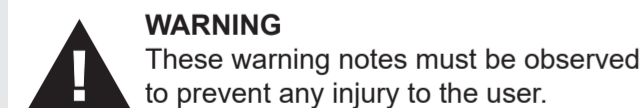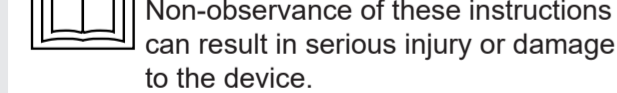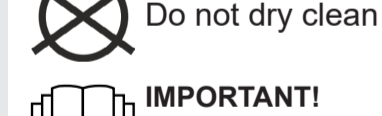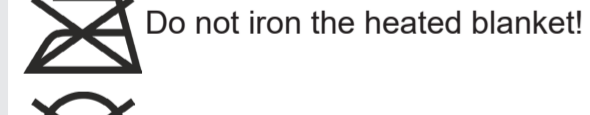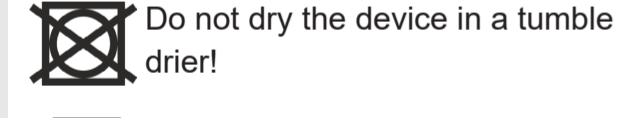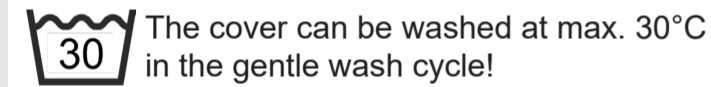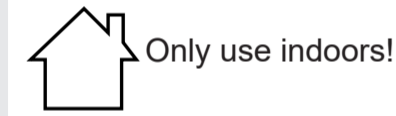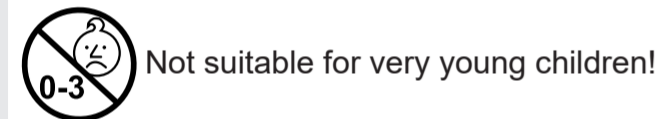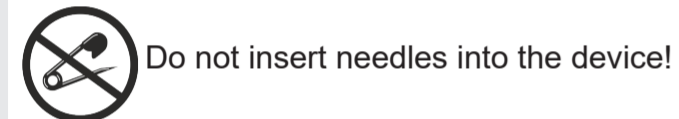**LOT** LOT number

**Manufacturer**

Recycling symbols/codes:  
These are used to provide information about the material and its proper use and recycling.

## GB Safety instructions

**Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.**

### About the power supply

- Only charge the unit's battery using the USB charging cable supplied. Never leave the unit unattended when it is connected to the mains (for charging) or in operation.
- Never charge the device in a flammable environment (e.g. flammable objects, gasoline, flammable gases, paint, etc.)
- Only connect the unit for charging so that the cable and power connector are freely accessible.
- The built-in batteries are not intended for replacement or removal by the end user. If problems occur with the power supply, contact the customer service.
- This device must only be supplied at safety extra low voltage (SELV) corresponding to the marking on the appliance.
- To disconnect the unit from the mains, always pull the plug out of the power connection. Do not pull on the cable!

### Specific groups of people

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children must be supervised to ensure that they do not play with the device.
- Do not use this unit to supplement or replace medical treatment. Chronic complaints and symptoms may be made worse.
- Do not use the **NM 990** neck massager, or only after consulting your doctor, if:
  - you are pregnant,
  - you have a cardiac pacemaker, artificial joints or electronic implants,
  - you suffer from one or more of the following illnesses and complaints: Circulation disorders, varicose veins, open wounds, bruises, broken skin, inflammation of the veins.
- Do not use the appliance near your eyes or near any other sensitive areas.
- If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
- The surface of the appliance is hot. People who are insensitive to heat must be careful when using the appliance.
- Before using the massager, consult your doctor regarding any unexplained pain, if you are undergoing medical treatment or using medical equipment.
- Children under the age of three must not use this device due to their inability to react to overheating.
- This appliance contains non-replaceable batteries.

### Before using the unit

- Check the massager carefully for damage before each use. Do not operate a defective device.
- Do not use the device if any damage is visible, if it is not working properly, if it has been dropped or if it has become wet. To avoid hazards, send the unit to the service centre for repair.

### Operating the unit

- The device is not to be used for commercial use or medical purposes.**
- If you have health concerns, consult your doctor before using the device.**
- Only use the massager indoors.**
- Do not use the massager in damp rooms (e.g. when bathing or showering).**

- The device must only be used for its intended purpose as described in the instruction manual.
- Using the device for any other purpose invalidates the warranty.
- Only use the device in the proper place, that is, around the neck and on the shoulders.
- Do not use the device while you are asleep or when you are lying in bed.
- Do not use the device before going to sleep. The massage has a stimulating effect.
- Never use the device while driving a vehicle or operating a machine.
- The maximum operating time for one session is 15 minutes. Prolonged use reduces the service life of the device.
- Very long periods of use can lead to overheating. Always allow the device to cool down before using it again.
- Never put the unit down next to an electric stove or other sources of heat.
- Do not allow the massage cushion to come into contact with pointed or sharp objects.
- Do not cover up the appliance when it is switched on. Never use it under blankets or cushions. There is a risk of fire, electric shock or injury.

### Intended Use

The **NM 990** massager is intended for targeted treatment, mainly of the neck and shoulder area, but can also be used for the upper or lower legs, abdomen or back. The pleasant massage has 3 levels of intensity and can be enhanced with heat treatment. The device is not intended for medical-therapeutic use. It is not suitable for clinics, massage practices, wellness studios, etc. The massage device is not a medical device and is not intended for the treatment of pain. In case of doubt or uncertainty, consult your doctor before use.

### Maintenance and cleaning

- The unit is maintenance-free. You may only carry out cleaning work on the unit yourself. In the event of malfunctions or damage, do not repair the unit yourself, as this will invalidate any warranty claims. Consult your specialist dealer and have repairs carried out only by authorised service centres in order to avoid hazards.
- Only use accessories intended for the unit and supplied by the manufacturer.
- Cleaning and user maintenance must not be carried out by children without supervision.

### Items supplied and packaging

Please check first of all that the device is complete and is not damaged in any way.

In case of doubt, do not use it and contact your supplier or your service centre.

The following parts are included:

- 1 **medisana** Neck massager **NM 990**
- 1 Charging cable, USB-C
- 1 Short instructions

The packaging can be reused or recycled. Please dispose properly of any packaging material no longer required. If you notice any transport damage during unpacking, please contact your supplier without delay.


**WARNING**  
Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!

### Charging the battery

- Use a 5V mains adapter with an output current of 2A. Connect the supplied USB charging cable to the adapter and to the USB-C port of the device.
- The LED **6** indicates the charging process by flashing slowly. When the battery is fully charged, the LED lights up constantly. The device cannot be used during the charging process. You can recharge the battery at any time.
- If the LED **6** flashes quickly the remaining charge is low – the battery should be recharged as soon as possible.

### Using the device

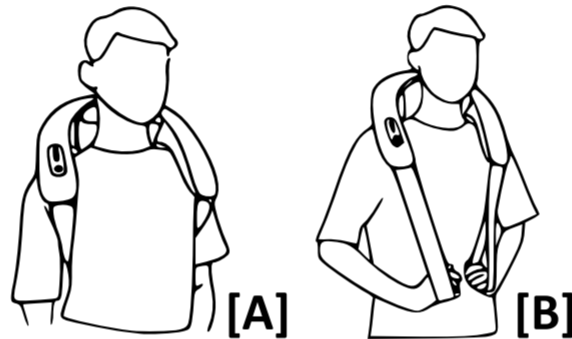Place the massager around the neck **[A]** so that it has sufficient hold on the shoulders. If you (optionally) fasten the back strap **8**, like shown in Fig. **[A]**, then you can use the **NM 990** hands-free.

- Press the On/Off button **1** for approx. 2 seconds to switch on the massage function. The device starts working straight away at low intensity. Press button **3** to select one out of 3 massage intensity levels, and button **2** to change the massage direction. Press the On/Off button **1** for approx. 2 seconds once again to switch all functions off.
- Press the HEAT button **4** to switch on the heat function. The LED indicator **6** changes its colour to red. To switch the heat function off, press the HEAT **4** button again.
- Put your hands through the hand straps **5** on the right and left sides **[B]**. By pulling on the straps you may stabilize the device in its position and adjust the massage intensity (pulling stronger increases the massage intensity).
- After 15 minutes of constant use the device will switch off automatically.

**CAUTION**  
Make sure that the maximum operating time of 15 minutes is not exceeded.

### Troubleshooting

#### • Device cannot be switched on or cannot be charged:

- Only use the appliance in rooms with the appropriate operating temperature. Failure to do so may result in functional impairments. If the device has been outside the permitted environment, place it in the correct environment for approx. 2 hours without further use and try charging or switching it on again.
- Also check whether the battery is sufficiently charged before switching the device on again.
- Device does not work or does not work normally, battery charge level too low:**
  - Always charge the battery in good time. Otherwise, the electronics may switch off the protection.

### Cleaning

- Before cleaning the appliance, make sure that the appliance is switched off and the power supply unit is disconnected. Allow the unit to cool down.
- Clean the massager with a slightly damp sponge only. Never use brushes, aggressive cleaning agents, petrol, thinners or alcohol. Rub the seat dry with a clean dry cloth.
- The neck cover **7** can be removed (zipper) and washed in the washing machine with mild detergent in the gentle wash cycle at 30°C.
- May not be bleached, ironed or dry cleaned!
- Never immerse the appliance in water to clean it and ensure water cannot get into the appliance.
- Wait until the device is completely dry before using it again.
- You should ideally store the device in its original packaging and keep it in a clean, dry place.
- Store the device in a safe, clean, cool and dry, well-ventilated place. Secure it against falling and shocks.

**Disposal:** This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your local authority or your supplier for information about disposal.

### Technical specifications

Name and model:	<b>medisana</b> Neck massager <b>NM 990</b>
Power supply:	Li-Ion battery 7.4V, 2200mAh Input: USB-C, 5V= 2.0A
Max. operating time:	15 minutes, with automatic switch-off
Operating conditions:	Only use in dry rooms at –5°C to +40°C
Storage conditions:	In a cool, dry place
Dimensions (L×W×H):	Approx. 23.5 × 15.5 × 19.5 cm
Weight:	Approx. 1.2 kg ± 0.1 kg
Article No.:	88963
EAN No.:	4015588 88963 9

### Warranty and repair terms

Your statutory warranty rights are not restricted by our guarantee below. Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is.

The following warranty terms apply:

- The warranty period for **medisana** products is 3 years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
- Defects in material or workmanship will be removed free of charge within the warranty period.
- Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
- The following is excluded under the warranty:
  - All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
  - All damage which is due to repairs or tampering by the customer or unauthorised third parties.
  - Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
  - Accessories which are subject to normal wear and tear.
- Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

**medisana GmbH,**  
Carl-Schurz-Str. 2, 41460 NEUSS, GERMANY

**In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.**

The current version of this instruction manual can be found under  
<https://docs.medisana.com/88963>.

Information about service can be found here:  
[www.medisana.com/servicepartners](http://www.medisana.com/servicepartners)

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