
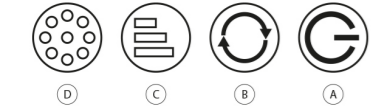
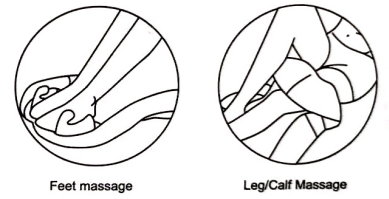
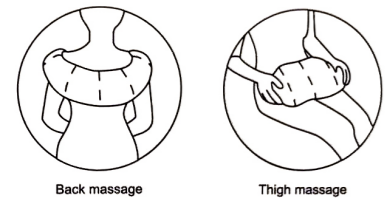
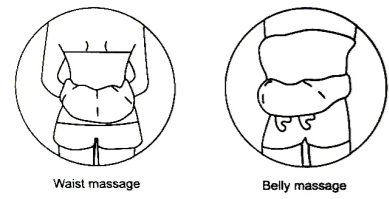
Rechargeable 6D Neck and Back Massager InnovaGoods

InnovaGoods®

User's Manual

Use the device for its intended purpose only. Fully charge the battery before its first use. Make sure that the voltage of the power source corresponds to that of the appliance. Do not pull on the cable and ensure it does not tangle. Keep the cable away from hot and sharp surfaces. Clean with a damp cloth. Do not use abrasive products for cleaning. Do not immerse the appliance or cable in liquid. Do not handle it with wet hands. Do not use near water or damp areas. Suitable for domestic use only. Use the device for its intended purpose only. Store the device in a dry, dust-free place. Do not leave the appliance unsupervised when it is turned on. If the appliance, the cable or the pieces looked damaged or do not work correctly, do not use them or attempt to repair them; they should only be repaired by a qualified technician. This device is not intended for use by children or persons with reduced physical, sensory or mental capabilities, unless they are supervised by an adult responsible for their safety. This device is not household waste. Please follow the environmental regulations in place for its recycling. Discard the product in accordance with the European Directive 2002/96/EC. Enquire about your local recycling legislation for electrical and electronic products and follow these steps. Do not dispose of this product with household waste. Following these instructions will help prevent harmful consequences to personal and environmental safety.

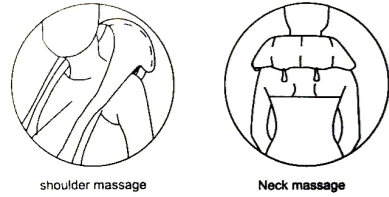CONTENT / TECHNICAL SPECIFICATIONS

- InnovaGoods Rechargeable 6D Neck and Back Massager
- DC In: 5V / 1A / 5W
- Battery: 2000mAh / 7.4V / 18.4W Li-ion
- Charging time: up to 180 min
- Usage time: 90 min approx.
- Charging and power LED
- 3 massage modes, 6 bi-directional nodes
- 3 speeds
- Heat function: 2 intensities
- Wavelength (IR): 650-700 nm
- Automatic shutdown: 10 min
- Hands-free straps
- Removable zip cover
- USB-C cable
- Materials: PU, PVC, PA, PP, ABS, PA, PP, ABS, silicone, iron, nylon, elastane, polyamide, polyester

INSTRUCTIONS

- To charge the InnovaGoods Rechargeable 6D Neck and Back Massager, plug the USB-C cable into the charging connector of the massager and connect to a USB power supply; the LED indicator will show red while charging and change to green when charging is complete.
- Place InnovaGoods Rechargeable 6D Neck and Back Massager on the area of the body you want to massage: neck, shoulders, back, waist, abdomen, legs or feet. Use the straps for optimal placement. You should remain in a comfortable and relaxed position during the massage.
- Press and hold button A on the control panel for a few seconds to turn on the device. After 10 minutes of use, the massager will turn off automatically. You can turn it off manually with the A button.
- Press button B to choose between the 3 massage modes: continuous clockwise rotation, continuous counter-clockwise rotation, alternating rotation with short or quick pause.
- Press the C button to select the massage speed: low, medium, high. It is advisable to choose a pleasant level of intensity.

ENGLISH

WARNINGS

IMPORTANT: KEEP THE INSTRUCTIONS FOR FUTURE REFERENCE. If you feel any discomfort when using the device, discontinue use and consult a doctor. To avoid risks, do not use it if any parts are damaged and do not attempt to disassemble or repair it. Do not use the device on areas where the skin is infected or sore. Suitable for domestic use only. Do not use abrasive products or materials to clean it.

- Press the D button to select the heat level (high, low) or to deactivate the heat function.

MAINTENANCE

The massage area cover has a zipper for easy removal and washing. You can clean the plastic outer area with a slightly dampened cloth.

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