# Rechargeable Portable Oral Irrigator InnovaGoods®


waste. Following these instructions will help prevent harmful consequences to personal and environmental safety.

## InnovaGoods®

### TECHNICAL SPECIFICATIONS

- DC In: 5 V / 1 A
- Battery: Li-ion 2000 mAh / 7.4 Wh
- Power: 4 W
- Pressure: 30-110 PSI
- Capacity: 220 ml
- IPX6 protection

## User's Manual
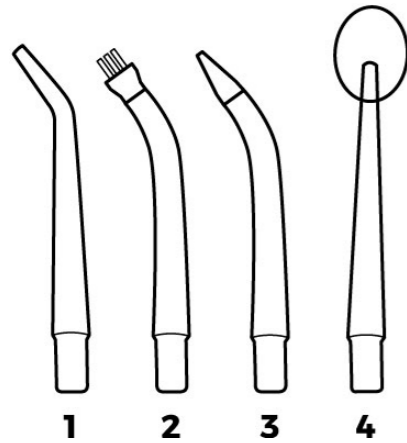
## ENGLISH

### WARNINGS

**IMPORTANT: KEEP THE INSTRUCTIONS FOR FUTURE REFERENCE.** Do not use chargers of more than 1A or fast charging chargers. Mineral water is recommended to prevent problems with limescale. Fully recharge the battery the first time. Make sure that the voltage of the power supply is compatible with the device. Do not pull on the cord and make sure that it does not become tangled. Keep the cord away from hot or sharp surfaces. Clean the device with a slightly damp cloth. Do not use abrasive products or materials to clean it. Suitable for domestic use only. Use the device for its intended purpose only. Store the product in a dry place with no dust. Do not leave the device unattended when it is switched on. If the device, the pieces or the cord show defects or they do not work properly, do not use them or attempt to repair them. Only a qualified technician should do it. This device is not designed for use by children or people with physical and mental disabilities, unless they are fully supervised by an adult who is responsible for their safety. This product is not household waste, so follow the environmental standards for recycling. Read the instructions carefully before first use. Keep the manual for further queries. Discard the product in accordance with the European Directive 2002/96/EC. Enquire about your local recycling legislation for electrical and electronic products and follow these steps. Do not dispose of this product with household

### INSTRUCTIONS

- Remove the rear silicone protector and charge the Rechargeable Portable Oral Irrigator InnovaGoods® using a USB source. The LEDs will flash while it is charging and will remain constant when charging is complete (approx. 3 hours). The device will not work while it is being charged.
- Open the container and insert water; you can add one part of your usual mouthwash (alcohol-free).
- Position the accessory you wish to use: (1) standard cleaning; (2) orthodontic cleaning; (3) periodontal cleaning; (4) tongue cleaning.
- Insert the nozzle into the mouth, above the teeth and close your lips. Press the ON button to switch the Rechargeable Portable Oral Irrigator InnovaGoods® on. Press again to change between settings: LOW (low pressure, for sensitive teeth or gums); HI (daily cleaning); PULSE (intermittent jet).
- Normally, a container of water is needed for the top teeth and another for the bottom teeth. With practice, you can avoid splashes and use the water more efficiently.

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