
InnovaGoods®

User's Manual

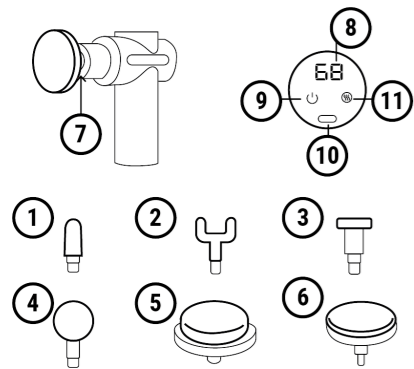ENGLISH

WARNINGS

IMPORTANT: KEEP THE INSTRUCTIONS FOR FUTURE REFERENCE. Clean and dry the skin before using the device. A maximum daily use of 15 minutes per muscle group is recommended. Fully recharge the battery the first time. Make sure that the voltage of the power supply is compatible with the device. Do not pull on the cord and make sure that it does not become tangled. Keep the cord away from hot or sharp surfaces. Clean the device with a slightly damp cloth. Do not use abrasive products or materials to clean it. Do not put the device or cord in liquid. Do not handle it with wet hands. Do not use the device close to a source of water or in damp places. Suitable for domestic use only. Use the device for its intended purpose only. Store the product in a dry place with no dust. If the device, the pieces or the cord show defects or they do not work properly, do not use them or attempt to repair them. Only a qualified technician should do it. This device is not designed for use by children or people with physical and mental disabilities, unless they are fully supervised by an adult who is responsible for their safety. This product is not household waste, so follow the environmental standards for recycling. Discard the product in accordance with the European Directive 2002/96/EC. Enquire about your local recycling legislation for electrical and electronic products and follow these steps. Do not dispose of this product with household waste. Following these instructions will help prevent harmful consequences to personal and environmental safety.

TECHNICAL SPECIFICATIONS

- DC In: 5 V / 1 A
- Battery: Li-ion 7.4 V / 1500 mAh
- Power: 25 W
- 6 Speeds (2600 to 3200 bpm.)
- 6 interchangeable heads
- LCD Screen
- Automatic switch off
- USB-C/DC cable

INSTRUCTIONS

- Charge the Hot & Cold Massage Gun InnovaGoods® in a USB source (10) for at least 5 hours the first time you use it. The LCD indicator (8) will show the charging status; when it reaches 99, it will be charged.
- Fully insert the head you wish to use into the hole at the front of the Hot & Cold Massage Gun InnovaGoods®. Do not insert or remove the massage heads with the device switched on.
- Various heads are included for different areas of use and functions:
- 1 - Bullet tip: deep massage, area of the joints and soles of feet.
- 2 - U-shaped: massage for the neck area, trapezius muscle and spine area.
- 3 - Flat tip: relaxing massage for any part of the body.
- 4 - Round: general massage for large muscle groups.
- 5 - Cold-effect tip: muscle relief.
- 6 - Heat-effect tip: muscle decongestion.
- Press and hold the «ON» button (9) for 2 seconds to switch the device on/off.
- Select the speed you require by pressing the «ON» button (9).
- Heat effect: Insert head (6) and press the «Heat» button (11). The screen (8) will display «HH» when the heat is active and «HO» when the heat is deactivated. The head will heat up until it reaches a temperature of approximately 40 °C.
- Cold effect: place head (5) in the freezer for 2 hours. Insert it into the hole (7). If it feels too cold, wait for a few minutes before use.
- Apply the massager vertically to the area wish to massage.
- We recommend you don't massage for more than 60 seconds on the same body area and that sessions do not exceed 15 minutes each day.

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