
Body Sculpting Belt InnovaGoods® -

InnovaGoods®

User's Manual

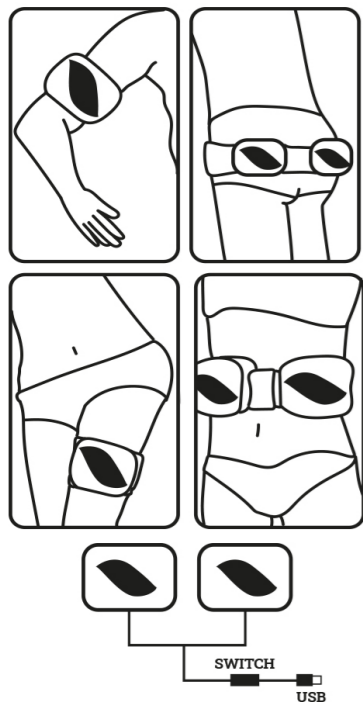www.innovagoods.com

ENGLISH

TECHNICAL SPECIFICATIONS

- DC In: 5 V / 600 mA
- Power: 2 x 18 W
- 2 «Body Sculpting Belt InnovaGoods®» modules
- Vibration: 4500 rpm
- 2 elastic bands: 44 cm
- 1 elastic band: 77 cm
- 1 USB-2DC cable

INSTRUCTIONS

- Pass the elastic band that you want to use through the metal side rings and behind each module of the Body Sculpting Belt InnovaGoods®.
- Plug a DC connector into each of the modules and the other end of the cable into a USB source.
- Activate the switch on the cable to switch the device on or off. You can use one module separately or both modules

together.

- Fit the Body Sculpting Belt InnovaGoods® using the elastic bands onto the area you want to massage (arms, legs, waist, hips, buttocks, etc.). Recommended daily use is between 10 and 15 minutes each day on each area of the body.

WARNINGS

Do not use the product if you are suffering from any inflammation or injuries in the area of application. If you feel any discomfort when using the device, discontinue use and consult a doctor. Do not use the device if you have heart problems, epilepsy or if you are pregnant. Do not use the device on areas where the skin is infected or sore. Make sure that the voltage of the power source corresponds to that of the appliance. Unplug the appliance whenever you are not using it. Do not pull on the cable and ensure it does not tangle. Keep the cable away from hot and sharp surfaces. Before cleaning, make sure that appliance is disconnected from the power supply and has cooled completely. Clean with a damp cloth. Do not use abrasive products for cleaning. Do not immerse the appliance or cable in liquid. Do not handle it with wet hands Do not use near water or damp areas. Suitable for domestic use only. Use the device for its intended purpose only. Store the device in a dry, dust-free place. Do not leave the appliance unsupervised when it is turned on. This appliance is not a toy. Keep out of reach of children. Do not modify the appliance or any of its components. If the appliance, the cable or the pieces looked damaged or do not work correctly, do not use them or attempt to repair them; they should only be repaired by a qualified technician. This device is not intended for use by children or persons with reduced physical, sensory or mental capabilities, unless they are supervised by an adult responsible for their safety. This device is not household waste. Please follow the environmental regulations in place for its recycling Discard the product in accordance with the European Directive 2002/96/EC. Enquire about your local recycling legislation for electrical and electronic products and follow these steps. Do not dispose of this product with household waste. Following these instructions will help prevent harmful consequences to personal and environmental safety.

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ALL RIGHTS RESERVED. MADE IN P.R.C.