
Body Shaping EMS Massager InnovaGoods®

InnovaGoods®

User's Manual

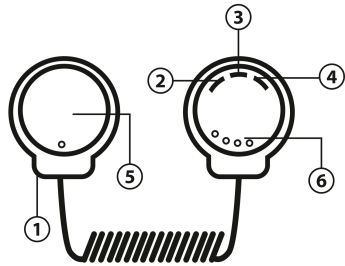ENGLISH

TECHNICAL SPECIFICATIONS

- DC In: 5 V / 500 mA
- Battery: Li-ion 150 mAh / 0.55 Wh
- 10 intensities.
- 6 massage modes
- 4 magnetic patches
- USB/DC cable

INSTRUCTIONS

- Plug (1) the Body Shaping EMS Massager InnovaGoods® into a USB source; the red LED (5) will switch on while charging and will switch off when charging is complete.
- Clean the skin on the area to be treated. Peel the transparent protective film from each electrode and keep it. Place both patches on the area of the body to be treated. Place each Body Shaping EMS Massager InnovaGoods® module on a patch; they will connect magnetically.
- Press and hold the «ON » button (2) to switch the device on. By default, it will switch on to the last programme selected. Change the programme from the following list (6) using the «MODE» button (3): **LED I:** Standard EMS (legs and arms) / **LED II:** Deep EMS (legs and arms) / **LED I + II:** Combined EMS (legs and arms)/ **LED III:** Standard EMS (buttocks) / **LED IV:** Definition EMS (neck and shoulders)/ **LED III + IV:** Combined EMS (abdomen).
- Press the «ON » button (2) to increase the intensity. Press the «OFF » button (4) to reduce the intensity of the massage or to switch the device off.
- The maximum recommended daily use is 10 minutes. When finished, stick the transparent film back onto the electrodes. If the protective film loses its adherence

over time, clean it gently with water and neutral soap and dry it with a clean cloth.

- Always start with the lowest intensity and increase it gradually as your body gets used to it (this may take several weeks).
- If the intensity is too high, reduce the intensity or remove one of the modules from the patch to stop the massage.

WARNINGS

Do not use the product if you are suffering from any inflammation or lesions in the area, high blood pressure, heart problems, if you have a pacemaker, if you are menstruating, if you are pregnant or suspect you may be pregnant. Fully recharge the battery the first time. Make sure that the voltage of the power supply is compatible with the device. Do not pull on the cord and make sure that it does not become tangled. Keep the cord away from hot or sharp surfaces. Clean the device with a slightly damp cloth. Do not use abrasive products or materials to clean it. Do not put the device or cord in liquid. Do not handle it with wet hands. Do not use the device close to a source of water or in damp places. Suitable for domestic use only. Use the device for its intended purpose only. Store the product in a dry place with no dust. If the device, the pieces or the cord show defects or they do not work properly, do not use them or attempt to repair them. Only a qualified technician should do it. This device is not designed for use by children or people with physical and mental disabilities, unless they are fully supervised by an adult who is responsible for their safety. This product is not household waste, so follow the environmental standards for recycling Discard the product in accordance with the European Directive 2002/96/EC. Enquire about your local recycling legislation for electrical and electronic products and follow these steps. Do not dispose of this product with household waste. Following these instructions will help prevent harmful consequences to personal and environmental safety.

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ALL RIGHTS RESERVED. MADE IN P.R.C.