
InnovaGoods®

User's Manual

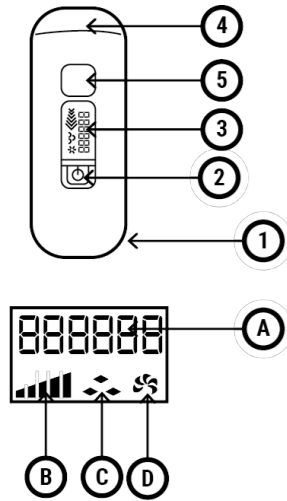ENGLISH

WARNINGS

IMPORTANT: KEEP THE INSTRUCTIONS FOR FUTURE REFERENCE.

Avoid looking directly at the pulsed light. Always wear the IPL safety glasses. Do not use the device if your skin is naturally dark or is very tanned. Do not use on energy levels from 2 to 5 if your skin has recently been exposed to the sun or artificial UVA rays. Do not expose the skin to sunlight without any protection after using it on energy levels from 2 to 5. For two weeks after each session, you should protect the skin with SPF 30 sunscreen or higher. Do not use around or near to the eyes. Exercise caution when dealing with facial hair. Do not use on tattoos, permanent makeup, skin spots, moles, warts, nipples, genitals or lips. It may cause burns, blisters or changes in skin colour. Do not use on areas with eczema, psoriasis, herpes, wounds or any skin condition. If in doubt, consult a doctor. Do not use on areas of the skin with a history of herpes outbreaks. Do not use if you are pregnant or while breastfeeding. There is a risk of optical radiation, of possible damage

to the sight or even loss of vision if you do not follow the instructions. If you apply the pulsed light for too long on the same area, it may damage the skin. If you suffer from photosensitivity, are taking photosensitising medication or take herbal remedies, consult your doctor before using the device. Do not remove grey or white hair. Do not use on areas that have been injected with fillers or toxins. Remove makeup, jewellery, earrings, piercings, etc. before each use. Avoid contact with liquids, humidity, switch off and unplug the device when you have finished using it if used in the bathroom. Do not immerse in the bathtub. This device has not been designed for medical or health purposes. Make sure the light emitter is in good condition and does not have any leaks, otherwise, do not use it. If you feel any discomfort when using the device, discontinue use and consult a doctor. If you notice excessive or prolonged redness, stop the treatment until the skin regains its normal appearance. If in doubt, consult your doctor. Clean the skin contact area regularly with a cloth and hot water to avoid the accumulation of oils. Suitable for domestic use only. Make sure that the voltage of the power source corresponds to that of the appliance. Unplug the appliance whenever you are not using it. Do not pull on the cable and ensure it does not tangle. Keep the cable away from hot and sharp surfaces. Before cleaning, make sure that appliance is disconnected from the power supply and has cooled completely. Clean with a damp cloth. Do not use abrasive products for cleaning. Do not immerse the appliance or cable in liquid. Do not handle it with wet hands. Do not use near water or damp areas. Suitable for domestic use only. Use the device for its intended purpose only. Store the device in a dry, dust-free place. Do not leave the appliance unsupervised when it is turned on. This appliance is not a toy. Keep out of reach of children. Do not modify the appliance or any of its components. If the appliance, the cable or the pieces looked damaged or do not work correctly, do not use them or attempt to repair them; they should only be repaired by a qualified technician. This device is not intended for use by children or persons with reduced physical, sensory or mental capabilities, unless they are supervised by an adult responsible for their safety. This device is not household waste. Please follow the environmental regulations in place for its recycling. Discard the product in accordance with the European Directive 2002/96/EC. Enquire about your local recycling legislation for electrical and electronic products and follow these steps. Do not dispose of this product with household waste. Following these instructions will help prevent harmful consequences to personal and environmental safety.

TECHNICAL SPECIFICATIONS

- AC: 100-240 V / 50-60 Hz - 0.8 A
- DC: 12 V / 2 A
- Wavelength: ≥ 530 nm
- Treatment surface: 3 cm²
- Pulses: 999999 (+/-1%)
- Working temperature: 5-45
- Accessories included: IPL protective glasses and shaver

PRIOR CONSIDERATIONS

- The duration of a cycle of hair growth is from 15 to 25 months. To achieve long-lasting hair removal, multiple sessions of IPL Hair Removal Device InnovaGoods® will be required. This depends on factors such as the area of the body the hair is

being removed from, the strength and colour of the hair and the way the IPL Hair Removal Device InnovaGoods® is being used.

- In a hair removal plan, it is advisable to leave two weeks between the first 4 treatment sessions, and 4 weeks between the following sessions until you achieve the desired results.
- The IPL Hair Removal Device InnovaGoods® has 5 energy levels: level 1 (4 j/cm²), level 2 (5 j/cm²), level 3 (6 j/cm²), level 4 (7 j/cm²), level 5 (8 j/cm²).
- As the energy level increases, so do the hair removal results, as well as the risks of possible side effects. Select the energy level that is best suited to your skin type.
- The IPL Hair Removal Device InnovaGoods® incorporates an eye protection system that prevents the emission of the light pulses if the treatment surface is not in full contact with skin. Always wear the protective glasses.
- Light-based hair removal is not suitable for people with dark or very tanned skin. On very dark skin, it may cause burns, blisters and changes in skin colour. Always use at the minimum power.
- IPL Hair Removal Device InnovaGoods® has 2 settings for use: manual and automatic.

INSTRUCTIONS

- Before starting each session, the hair should be shaved and skin should be clean and dry. Hair should not be removed beforehand using wax, epilator or tweezers. Use the shaver.
- Skin must not have been exposed to the sun for at least 2 weeks before treatment.
- Plug the adaptor into a power socket and to the connector (1) of the device; the screen will flash for a few seconds.
- To switch the device on or off, press and hold button (2) for 3 seconds; this will activate the fan and screen. The screen will display the number of flashes remaining (A), the icon of the power selected (B), the icon of the programme type (C) (manual or automatic) and the fan icon (D). Press and release button (2) to select the next energy level (B) on the screen (3).
- Put on the IPL protective glasses.
- Automatic setting: Press and hold button (5) for 3 seconds and a screen will be displayed with an icon showing 3 dots (C). Place the device head (4) in a fully vertical position on the skin and glide it slowly without losing contact with the skin. Avoid treating the same area more than once or overlapping areas of skin. The pulses will vary between 1 and 3 seconds depending on the power selected.
- Manual setting: Press and hold button (5) for 3 seconds and the screen with the icon showing 3 dots (C) will disappear. Place the device head (4) in a fully vertical position on the skin and press button (5) once to emit a pulse; if you don't feel any discomfort, increase the energy level. Try different energy levels until you reach the highest that doesn't cause discomfort. Move the head to the next area of skin and repeat the process.
- If, 48 hours after starting the treatment, you don't experience any adverse effects, you can continue with the full treatment at the same energy level or increase it, repeating this test on each part of the body to be treated.
- It is normal to feel a slightly warm feeling and tingling during the treatment, as well

- as having a slight redness of the skin. If you observe a severe redness of the skin, blisters or burns, stop the treatment immediately.
- Avoid exposure to the sun on any areas of the skin treated with energy levels from 2 to 5. Ensure you carefully protect the treated skin with sunscreen throughout the hair removal period and for at least 2 weeks after the last session.
- Switch the device off by pressing button (1); a countdown from 10 seconds will appear on the screen; do not unplug the cable until the screen displays «OFF».
- Clean the head with a slightly damp cloth and then dry it.

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