NEXXT

Model: CO250924-1

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NL

EN

DE

FR

ES

NEXXT

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ES


NEXXT

Thank you for purchasing the foot massager!


- Before using, please read the instructions carefully, and use it correctly.
- Please keep this "User Guide" properly, in order to keep readying anytime for reference.
- This product is designed for home scenes, and is also suitable for portable use in other indoor scenes (such as business trips, etc.)

Components & How to use


1. Plug the AC adapter into the socket.
2. Connect the machine to the power cord and make sure it is powered on.
3. Select your favorite massage program and start it.
4. Please disconnect the power supply after use.

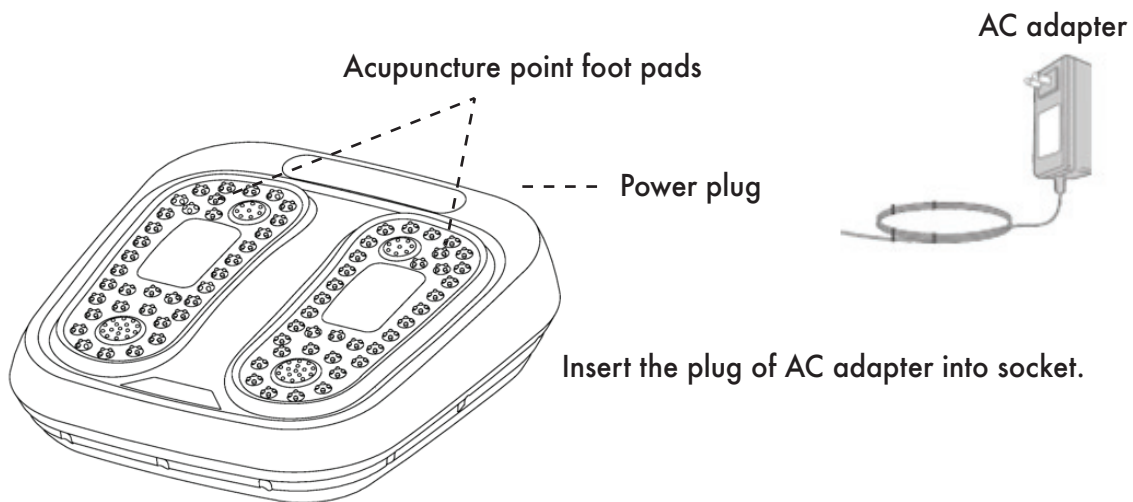Demonstration:

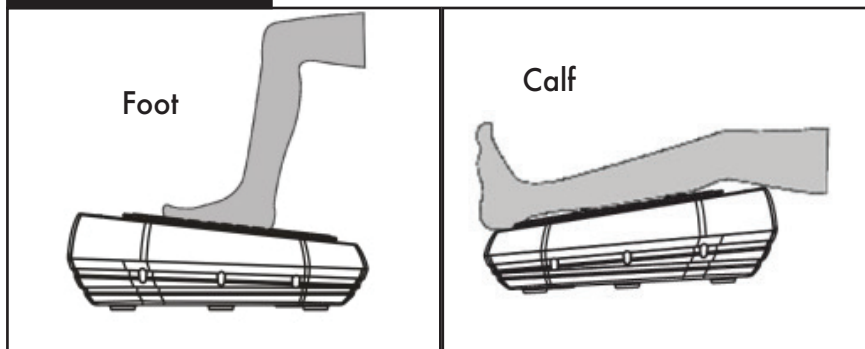Please place your calves or feet on the machine.
The acupuncture point foot pads can be removed for cleaning.
Please use it after it is completely dried after cleaning.

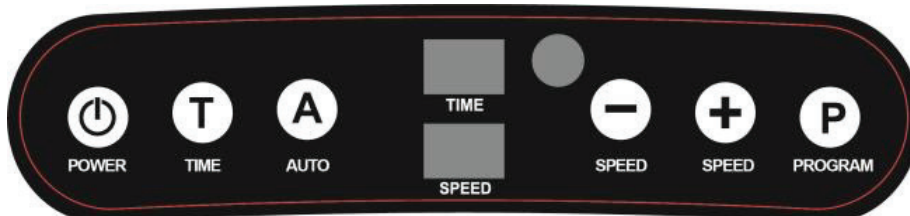SAFETY WARNINGS


- Read the instruction manual before using the NEXXT FOOT MASSAGER WITH SHIATUS for the first time.
- Always place the NEXXT FOOT MASSAGER WITH SHIATUS on a flat surface.
- Ensure that children do not use the NEXXT FOOT MASSAGER WITH SHIATUS without adult supervision.
- Do not allow the NEXXT FOOT MASSAGER WITH SHIATUS to make contact with water.
- If the power cable is damaged, it should be replaced with a special cable by the manufacturer or authorized technical service.
- Unplug the device after use.
- As with any other physical exercise device, consult your doctor before you start using the NEXXT FOOT MASSAGER WITH SHIATUS.
- This device may be used by children aged eight and over, by people with reduced physical sensory or mental abilities, and by people lacking in experience and knowledge if they are appropriately supervised or given suitable instructions on the safe use of the device and understand the dangers its use entails.
- Children should not be allowed to play with the device.
- Children should not clean or do maintenance on the device without supervision.
- Consult your physician before using this product if:
 - You are pregnant
 - You have a pacemaker, an implant or other aids
 - You suffer from diabetes
 - You have a circulatory disorder
 - You have concerns about your health
- The elderly and infirm, should consult their doctor before using the product, even if there are no specific pre-existing conditions. Because skin friction can be painful, it is recommended that the device not be directly applied to the skin.
- Stop using immediately if you feel discomfort or pain.
- Do not use this product at bedtime.
- The massager has a stimulating effect and can delay sleep.
- To avoid a short circuit, connect the device correctly and keep it away from metal (needles), waste and water.
- Remove the adapter from the power outlet when not in use.
- If you suspect a problem, turn off the device as soon as possible. When not in use for extended periods store the device in its packing in a safe, dry cool place.
- Do not place any other objects on the massager.
- If the power cord is damaged, it must be replaced, but only by an authorized and qualified person.
- To reduce the risk of electric shock, do not use the device on wet skin.
- Keep the product away from moisture.
- Do not use the product if it or its parts are wet or damp.
- Children should not be allowed to play with the device.
- Children should not be responsible for cleaning or maintaining the device, unless they are over the age of eight and supervised by an adult.
- Keep the device and its cable out of the reach of children under the age of eight and of pets and animals.
- To avoid injury and irreparable damage always sit while using the massager.
- To avoid dangerous injury, never stand on the unit or use it forcefully.
- Do not attempt to disassemble or repair the product.
- Repairs may only be carried out by suitably qualified technicians.

Forbidden

- Do not connect the device to a voltage transformer or other accessories.
- Do not use it if the power supply is damaged.
- Do not run the device idle. Please disconnect the power supply when not in use.
- Do not pull or twist the power cord.
- Do not use the device in a wet environment.
- Do not use the device with more than one person. Please be careful not to use the device while holding a baby.
- Do not stand on the device.
- Do not use the device during thunderstorms.
- The device is only suitable for feet and legs. Do not use it on other parts of body.

Specification

Model No.	: CO250924-1
Rated voltage	: AC100V-240V
Rated power	: 30W
Rated frequency	: 50-60Hz
Rated time	: 10Min
Product setup size	: L363mmXW396mmXH125mm
Product weight	: 4.2kg
AC adapter	: 1PC
Carton package	: 1pc/carton

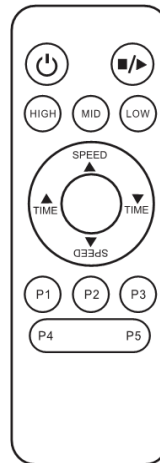	: On / Off
	: Start / Stop
	: Speed increase
	: Speed decrease
	: Time increase
	: Time decrease
	: select 16 speeds
	: select 10 speeds
	: select 6 speeds

P Mode (P1,P2,P3,P4,P5) : Time and speed can not be adjustable under P mode.

Manual mode (H0) : Time and speed are adjustable under manual mode.

Power-on display : Default speed :00, Default time :10 minutes.

Press the Start/Stop button : The device will stop operating while it is running.

- In shut down state, the display will light up by pressing any button.
- In power-on state, press the start/stop button to start/stop the device.
- After starting, press the +/- button to increase/decrease the speed, the speed range is from 1 to 20.
- In pause state, under manual mode, press +/- to increase/decrease the setting time. The time can be set from 1 to 10 minutes and the default is 10 minutes.

Programs

- P1: Vibration + massage, time and speed can be adjusted manually.
- P2: Vibration only, speed cycle according to the P2/P3 table.
- P3: Vibration + massage speed cycle according to the P2/P3 table.
- P4: 1 minute of vibration + 1 minute of massage, speed cycle according to the P4/P5 table.
- P5: 2 minutes of vibration + 3 minutes of massage, speed cycle according to the P4/P5 table.

- Automatic program P4/P5 table, the speed changes once per second, and the cycle per minute is as follows:
1,1,2,3,4,5,6,7,8,9,10,11,12,13,14,15,16,17,18,19,20,19,18,17,16,17,18,19,20,20,19,18,17,16,17,18,19,20,20,19,18,17,16,15,14,13,12,11,10,9,8,7,6,5,4,3,2,2,1,0

- Automatic program P2/P3 table, the speed changes once per second, and the cycle per minute is as follows:
1,1,2,3,4,5,6,7,8,9,10,11,12,13,14,15,16,17,18,19,20,19,18,17,16,17,18,19,20,20,19,18,17,16,17,18,19,20,20,19,18,17,16,15,14,13,12,11,10,9,8,7,6,5,4,3,2,2,2,1

Warranty

- During the warranty period, if you encounter the following faults, you can repair them free of charge:
 - If you find any damage or defect when purchasing, please contact the nearest after-sales service.
 - Damage occurs when operating according to the instructions.
 - Non-consumable items have been damaged.

- If the following situations occur, the user needs to cover the cost of repairing:
 - Damage caused by improper operation.
 - Damage caused by fire, natural disasters, abnormal voltage.
 - Damage caused by modification or repair by unofficial personnel.
 - The product name and purchase date are not indicated in the warranty card held by the user.

- The warranty card is only valid in the country of purchase.
- The warranty card can not be reissued if missing, please keep it properly.

Correct Disposal of this product

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

NEXXT

is the innovative brand with high-quality gadgets for in and around the home. Go for sustainable comfort with our inspiring items and make your life easier, simpler or more relaxed. Always at a competitive price. Because we believe in affordable innovation for modern life. NEXXT, the new ease.

Thank you for choosing one of our NEXXT products!

Before you get started, we strongly recommend that you read this entire manual carefully to avoid possible bodily harm and damage to the product. We cannot be held responsible for the consequences of incorrect or faulty installation if the instructions in this manual are not followed carefully. If you have any questions, problems with the installation or have any doubts about safety, please feel free to contact us at info@conceptonlinebv.nl.

Importer:
Concept Online BV
De Kreek 28
NL -3871 MM Hoevelaken
info@conceptonlinebv.nl

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