BR3ND

BR-2230 Massagematras

Massage mat / Massagematte / Matelas de massage
Materasso da massaggio / Colchón de masaje / Mata do masažu

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Polish manual 45

Icons used in this manual

The following warnings are used in this manual and on the packaging

Warning symbol: This symbol indicates important safety information. When you see this symbol, pay close attention to the accompanying instructions to prevent injury or damage. Follow the precautions carefully to ensure the safe operation of the device.

Information symbol: This symbol provides additional information and helpful tips that can assist you in using your device. When you see this symbol, read the accompanying text for further explanation and guidelines to enhance functionality and ease of use.

Introduction

Congratulations on the purchase of your new massage mat! This comfortable and versatile product provides the ideal solution for relaxation and relief from muscle tension, wherever and whenever you desire. This manual contains instructions for the correct installation, use, and maintenance of your massage mat. Please read this guide carefully to make the most of your purchase and ensure the long-lasting and safe operation of your massage mat. Thank you for choosing our product!

Component inspection


Before you begin, it is important to check that all components of your massage mat are present. Carefully inspect the contents of the packaging to ensure you have received everything necessary for the optimal operation of your device.

- 1x Massage mat
- 1x Massage pillow
- 1x Handheld control panel

Description of parts


Pillow controls

1. On/off Button
2. Direction
3. Intensity
4. Heating

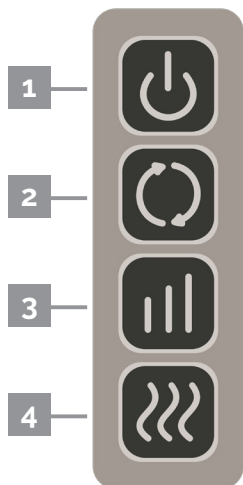Mat controls

- | | |
|--------------|--------------|
| 1. Shoulders | 6. Timer |
| 2. Back | 7. Mode |
| 3. Buttocks | 8. Intensity |
| 4. Legs | 9. Heating |
| 5. Standby | |

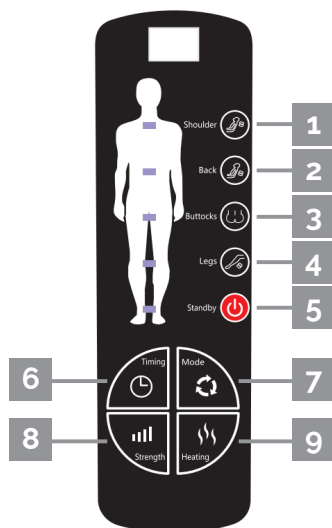Safety warnings

ATTENTION! Please read all instructions and warnings carefully before using the massage mat. Use the device only for its intended purpose.

This product is **not** suitable for...

- Elderly, injured, ill, weakened, or bedridden individuals, seriously ill patients, or those undergoing cancer treatment, or suffering from conditions such as acute spinal cord injuries, heart problems, stomach ulcers, stomach bleeding, lung diseases, severe heart conditions, liver diseases, malignant tumours, or acute abdominal issues: consult a doctor before using the product.
- Pregnant women, infants, and young children.
- Individuals with blood disorders, an increased risk of bleeding, in the early stages of closed soft tissue injuries, bruises (purple spots), or fresh bone fractures. Also, do not use on an empty stomach, immediately after a heavy meal, after excessive alcohol consumption, or after intense physical exercise. Avoid using if recovering from nerve or vascular surgeries.

Safety warnings:

- Unplug the device after use or before cleaning.
- This device should not be used in damp environments or bathrooms. Store it in a dry place, out of direct sunlight, after use.
- Do not use the device while sleeping to avoid prolonged stimulation of a specific body part.
- Do not use the massage pillow within half an hour after a meal, as various bodily reactions may cause undesirable effects.
- This device is not intended for use by individuals (including children) with physical, sensory, or mental impairments, or by those lacking experience or knowledge, unless supervised by a person responsible for their safety.
- Children and individuals with limited mobility should use the product under appropriate supervision.
- During pregnancy, or if suffering from back pain, acute illnesses, heart problems, or other physical abnormalities, only use under the guidance of a physician.
- Drink 500 ml of water or warm tea after the massage to help eliminate toxins from the body.
- If the power cord of the device is damaged, contact the manufacturer; do not use unsuitable replacements.
- If the surface of the device becomes hot during use, persons insensitive to heat should exercise extra caution.
- This product is not recommended for individuals with disabilities, sensory or neurological disorders, or those lacking experience and common sense (including children, individuals with intellectual disabilities, and people with physical paralysis). When used, it must always be under the supervision of a responsible adult. Ensure that children do not play with the product.

Instructions

NOTE:

Keep this manual in a safe place, as it contains important information.

Instructions for the neck control panel

Standard operating time for the pillow: 15 minutes


1. **On/off:** Briefly press this button to start the pillow. The neck rollers will begin to operate at high intensity, and the heating function will be activated. Press this button again briefly to turn off the pillow.
2. **Direction:** Press this button to change the massage direction of the rollers. Choose between forward kneading or reverse "Tuina" massage.
3. **Intensity:** The device has 3 intensity levels: low, medium, and high. Press this button to select the desired intensity level.
4. **Heating:** Press this button to toggle the heating function on or off.

Instructions for the handheld control panel

Standard operating time for the pillow: 10 minutes

- Standby:** Briefly press this button to start the pillow. The default operating time is 10 minutes in automatic mode (AU). Press this button again briefly to turn off the pillow.
- Time:** Press this button to adjust the duration between 10, 20, or 30 minutes.
- Mode:** Press this button to switch between the vibration modes of the pillow, including 5 different modes from "F1" to "F5", and the automatic mode "AU".
- Intensity:** Press this button to switch between 3 intensity levels: "S1" (low), "S2" (medium), and "S3" (high).
- Heating:** Press this button to turn the heating function of the pillow on or off. There are 3 heating levels: "H1" (low), "H2" (medium), and "H3" (high). When the display shows "H0", the heating function is turned off.
- Shoulders:** Press this button to direct the vibration to the shoulders. The display will show "P1".
- Back:** Press this button to move the vibration to the back. The display will show "P2".
- Buttocks:** Press this button to direct the vibration to the buttocks. The display will show "P3".
- Legs:** Press this button to move the vibration to the legs. The display will show "P4".

Operating modifications

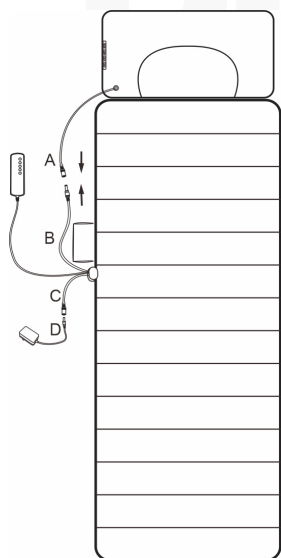Option 1 Connect A and B, and connect C and D to operate both the pillow and seat cushion simultaneously.

Option 2 Connect only A and D to operate the pillow independently.

Option 3 Connect only C and D to operate the seat cushion independently.

*A is the cable from the pillow
B and C are the cables from the mat
D is the adapter*

Troubleshooting

First, check the following issues yourself. If problems persist, please contact us or your local distributor for assistance.

Problem	Cause	Solution
Entire device not working	Power outage	Wait for the power to be restored
	Adapter is not plugged into a socket	Properly connect the adapter to the socket
	Power switch is off	Turn on the power switch
	DC plug from the adapter is not properly connected to the mat	Firmly insert the DC plug into the DC socket of the mat
Noise	Device is not positioned correctly	Position the device correctly and securely
Low intensity	Mains voltage is too low	Use a voltage regulator
Mat not working or defective	Possible defect in the mat	Unplug the adapter, do not use the device, and contact the manufacturer

Note: Users should not disassemble or repair the massage mat themselves. Self-repair may void the warranty, and the company will not be liable for any consequences.

Disposal of used product

Do not dispose of this product/device with normal household waste at the end of its lifespan. Take the product to a public collection point for the recycling of used electrical and electronic devices. The symbol indicating this can be found on the product, the user manual, or the packaging.

Technical data

Product name	Brend BR-2230 Massage Mat
Voltage	12V
Current	3A
Power	36 Watts
Input	100-240V ~ 50-60Hz
Output	12V DC / 3A
Materials	Fleece fabric; foam filling; ABS plastic
Separate use	Yes
Heat compression	Yes
Number of vibration positions	4 (shoulders, back, buttocks, and legs)
Number of heating levels	Mat: 3 heating levels Pillow: 1 heating level
Operating times	10, 20, and 30 minutes
Intensity	5 adjustable modes
Massage methods (pillow)	Kneading and heating
Massage methods (mat)	Vibration and heating
Cable length	2 meters
Pillow dimensions (LxWxH)	45 x 28 x 10 cm
Dimensions (LxW) (mat + pillow)	170 x 55 cm
Country of origin	China
Rovus Trademark Owner, Importer to the EU, Distributor	BREND Grote Ohé-weg 3, 6081PT Haelen Netherlands

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