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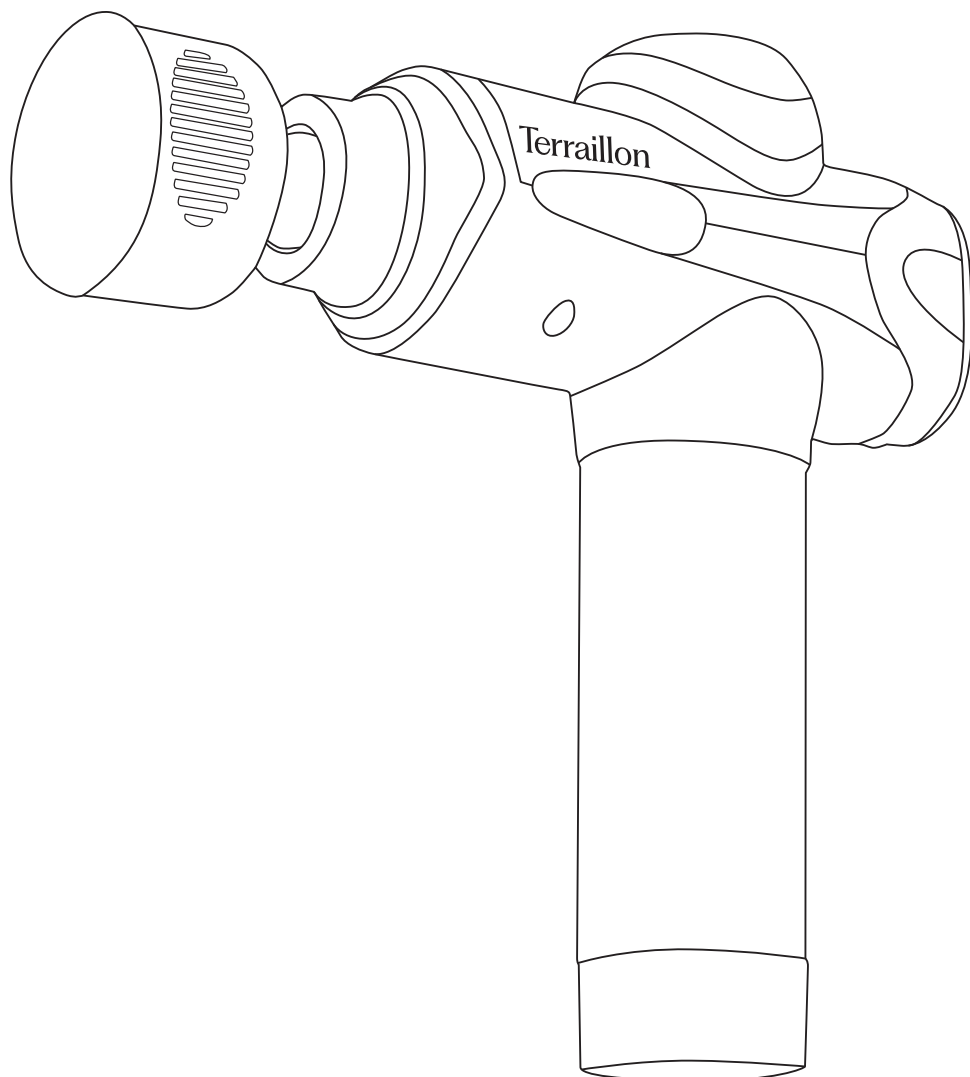
Dear Customer,

Thank you for purchasing the Gun Massager Expert. Please read the instructions before using the massager as this manual provides additional information including safety precautions and warnings.

If you have any difficulties or need more information, you can find our help centre by scanning the QR code opposite with your smartphone. If you cannot find the answer to your question, you can also contact our After Sales Service: support.uk@terraillon.com

DISCLAIMER

CAUTION

If you have the following conditions or any health issues, please consult with your doctor before use.

1. Pregnancy, diabetes including complications such as neuropathy and retinal damage.
2. Use of cardiac pacemaker, recent surgery, epilepsy or migraines, protrusion of lumbar intervertebral discs and anterior displacement of bone vertebrae.
3. Recent joint replacement surgery, intrauterine device placement, metal implants and other concerns of the body.

WARNING: consult doctor before using this product to reduce the risk of injury.

PRECAUTIONS

1. Minors are not allowed to use this product. This product is not a toy, please store it out of reach of children.
2. Do not immerse the product in water. If the product falls into the water, please submit it for inspection in time.
3. Do not place the product near the fire source or into the fire.
4. Do not use the gun massager on the head, near the heart or on the genitals.
5. Do not use the gun massager on the bones and joints. The massager is designed for use on soft tissues.
6. Do not drop the product intentionally.
7. Do not remove screws or attempt to disassemble the product.
8. It is forbidden to replace or supplement medical treatment with the device, if you feel pain or are receiving treatment, please consult before using the device your doctor.
9. Do not put fingers, hair or other parts of the body near the axis and back of the massage head to avoid crushing injuries.
10. If the equipment is damaged, stop operation immediately.
11. Do not place fingers, hair or other body parts near the axis or back of the massager to prevent injury.

PACKAGE CONTENTS

Gun Messenger x 1

Massage Head x 4


Heat & Cool Massage Head x 1

Charging cable x1

Carrying case x 1

Instruction manual x 1

PRODUCT OVERVIEW

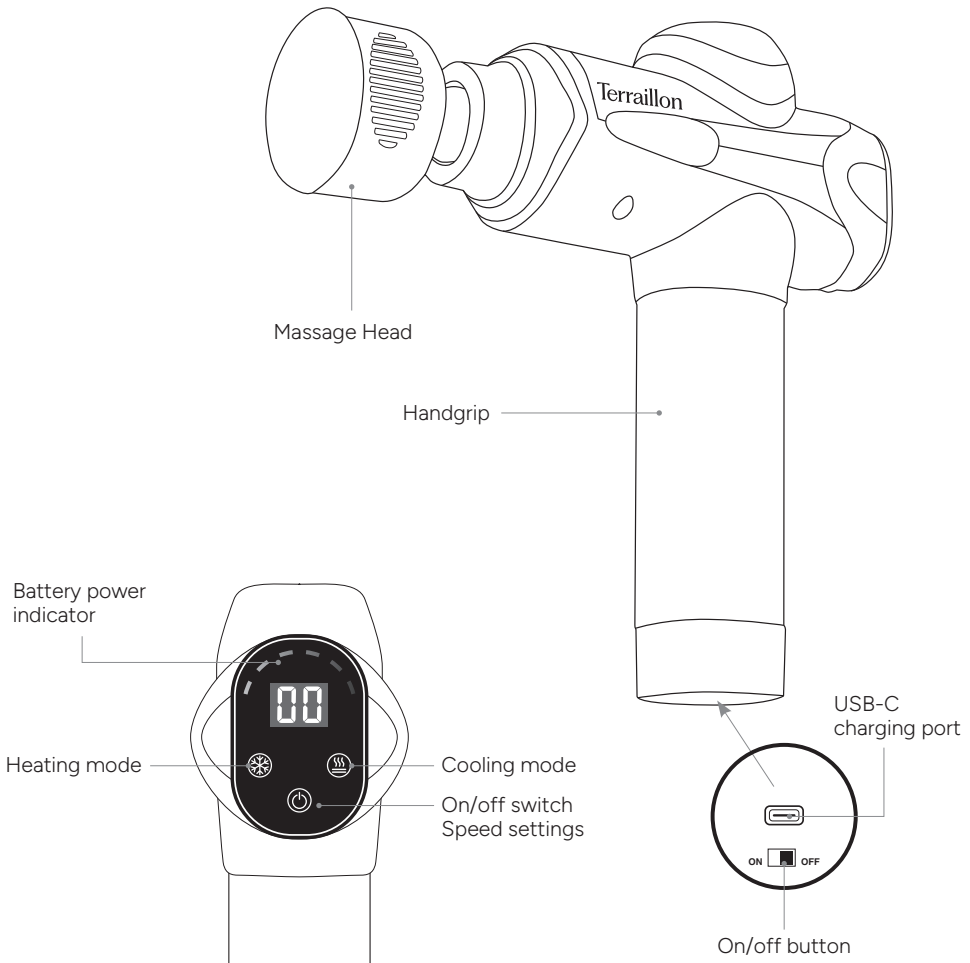CHARGING THE PRODUCT


⚠ For your safety, always use an appropriate 5V input charger/adaptor to charge the massager. DO NOT USE a quick-charging USB adapter as this may damage the product.

Important: The Gun Massager must not be turned on while charging.

1. Use the USB-C cable included with this massager.
2. Plug the USB-C plug into the USB-C charging port located under the massage gun handle.
3. Connect the USB plug to a USB port on a standard wall adapter to begin charging.
4. It is recommended to fully charge the battery for five hours before its first use.
5. When the Gun Massager is charging, the battery level indicator light flashes. When fully charged, all battery indicator lights are on.
6. Put the massage gun on charging, if the battery charging is at level 2.

NOTE:

- It is not recommended to drain the battery completely.
- You must select a qualified charger from a reputable manufacturer.
- You must charge the massager if it has not been used for more than 6 months.

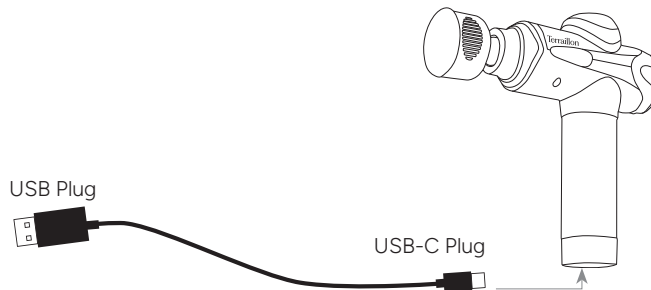⚠ Safety instructions for charging the product's lithium battery:

- The input voltage is DC 5V 2A.
- In order to charge the product safely and reliably, please use a wall adapter that is certified and in good condition.
- Only charge the product in ambient temperatures between 10°C /50°F and 40°C /104°F.
- Do not heat above 70°C /158°F. Do not dispose of in a fire.
- Keep out of direct sunlight.
- Turn off the product after use.
- Please note: we will not take any responsibility for the misuse of the product, such as using voltage higher than indicated. Misuse may result in serious injury or loss of property.
- In extreme cases, abuse or misuse of a lithium battery can lead to an explosion, heat generation, fire or smoke development.

MESSAGE HEADS

- **Spherical Head:** suitable for Relaxing Large Muscle Groups, including your arms, waist, back, buttocks, and thighs.

- **Flat Head:** suitable for Relaxation and Shaping of Various Muscle Groups and especially the abdominal belt, the pectoral muscles of the chest or even the biceps and triceps femoris in the leg.

- **Cone Head:** suitable for deep tissue impact recovery sessions. Ideal for targeting trigger points and small specific areas: such as the soles of the feet, palms, Achilles heel and rhomboid muscles in the back.


- **U-Shaped Head:** suitable for trapezius muscles, areas around the achilles tendon muscle, ankle muscle and calf muscles.

- **Heat & Cool massage head:** to help relax muscles and promote blood circulation (warmth) and support muscle regeneration and soothe painful tissues (freshness).

CONNECTING THE ATTACHMENT

- Place the massage head you want on the end of the gun and snap it on. Then turn on the gun using the on/off button under the massage gun handle to start the massage.
- For the Heat & Cool tip, fit the vertical bar of the tip into the notch on the gun. Make sure it is inserted all the way to the bottom.
- To change the Massage Head, turn off the gun by holding down the ⏻ button: you can remove the used head and insert the new one

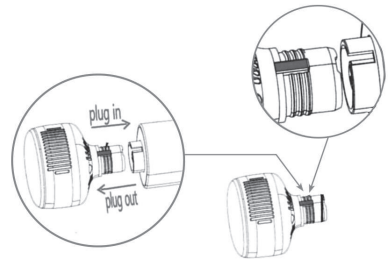USING THE PRODUCT

1. Long press the ⏻ button of the gun massager to start the massager. The massage will start at the slowest speed and the light will turn blue.
2. Briefly press the ⏻ button to change the intensity.
3. The gun massager features twenty speeds.
4. The Gun Massager switches off automatically after each 15-minute session, corresponding to the recommended massage time. To restart, simply long press the ⏻ button.
5. If you wish to stop the session before the 15-minute timer, simply long press the ⏻ button until the massager powers off.

USING THE HEAT & COOL TIP

- Hot and cold therapy can only be used with the Heat & Cool tip.
- For hot therapy on/off mode, press the  hot therapy mode button. The heating function will get start. Hot therapy has 3 temperature levels:

- the initial temperature will be about 38°C and the indicator turns blue.
- when the button is pressed again for the second time, then the temperature will turn about 40°C and the indicator will turn to green.
- press the button again for the third time. Then the temperature will increase to 42°C and the indicator will turn red.

	blue	38°C
	green	40°C
	red	42°C

According to this operation cycle there will be a need to long press the hot therapy mode key for two seconds to turn it off.

- For cold therapy on/off mode, long press the  cold therapy mode button for two seconds, the cold function will get start. Cold therapy has 3 temperature levels :

- the initial temperature will be about 20°C and the indicator turns red.
- When the button will be pressed again for the second time then the temperature will turn about 18°C and the indicator will turn to green.
- press the button again for the third time. Then the temperature will increase to 16°C and the indicator will turn blue.

	red	20°C
	green	18°C
	blue	16°C

- Vibration, hot therapy and cold therapy functions can be operated independently.

FAQ

• Why is my Gun Messenger Expert not turning on?

If the Gun Messenger Expert will not power on, please ensure the battery is fully charged and that the on/off button on the underside of the appliance is in the on position.

• Why is my Gun Messenger Expert not charging?

If your Gun Messenger Expert will not charge, please ensure each end of the charging cable is securely attached to the Gun Messenger and the power source. Check to ensure the wall outlet has power. If using a computer USB port, please ensure that the computer has power.

• Why does my Gun Messenger Expert stop suddenly?

This function activates when the head is locked. Remove the obstruction or cause of the lock and check the head can moves freely. Or the battery is discharged, charge the device.

• Why does my Gun Messenger Expert display an E1 message?

This message appears when the battery is at its minimum charge level. The product must be charged before it can be used again.

TECHNICAL SPECIFICATION

Product name	Gun Massager Expert
Vibration mode	20 levels
Speeds	2800 rpm - 3200 rpm
Charging time	about 5 hours
Charging interface	USB-C
Usage time	2-4 hours
Battery specifications	Li-Ion / 11,1 V / 2 000 mAh / 22,2Wh
Product size	240 x 180 x 70 mm
Hot mode temperature	38°C, 40°C, 42°C
Cool mode temperature	20°C, 18°C, 16°C

The temperature is indicative and may vary depending on skin and environmental temperature.

CARE AND MAINTENANCE

To reduce the risk of failure, electric shock, injury, fire and equipment damage, please read all the instructions and warnings in this manual.

1. Do not dispose of, store or use the massager near or in fire or open flames.
2. Do not remove screws or attempt to disassemble the massager.
3. Please consult with your doctor before use. This massager is not intended to replace or supplement prescribed medical treatments.
4. Clean the device with a slightly wet towel, wipe and dry it with dry soft cloth.
5. Please charge the device at least once a month, if it is not in use for a long time.

WARRANTY

This product is guaranteed against defects in materials and workmanship. This warranty does not cover damage resulting from accidents, improper use or negligence. If you have a complaint, first contact the store where you made your purchase.

ENVIRONMENTAL PROTECTION

When the product is no longer usable, please take it to a collection point for the recycling of waste electrical and electronic equipment.

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Terrailon SAS France & Headquarters
Terrailon Asia Pacific Ltd
Terrailon Corp USA