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TERRAILLON CONNECTED BODY COMPOSITION ANALYZER

Dear Customer;

Thank you for purchasing this product. Please read these instructions carefully to take full advantage of your product

- Don't use this product if you wear a pacemaker or implantable electronic device.
- Any information provided by this device is in no way meant to treat, cure or prevent, any disease or illness from happening.

Safety instructions for charging the product's lithium Ion Polymer battery:

For your safety, always use an appropriate 5V charger/adaptor to charge the scale. DO NOT USE a quick-charging USB adapter as this may damage the product.

1. Use the Type-C cable included.
2. Connect the USB-C plug to the charging port on the side of the scale.
3. Connect the USB end to a USB socket on a computer or standard wall adapter to start charging.
4. When the scale is charging, the battery picto on display will blink. When fully charged, the battery indicator will remain steady.

NOTE:

- You must select a qualified USB charger from a reputable manufacturer.
- You must charge the scale every 3 month.
- Only charge the product in ambient temperatures between 10°C /40°F and 40°C /104°F.
- Do not heat above 70°C /158°F. Do not dispose of in a fire.
- Keep out of direct sunlight.
- Terrailon will not take any responsibility for the misuse of the product, such as using voltage higher than indicated. Misuse may result in serious injury or loss of property.
- In extreme cases, abuse or misuse of a lithium battery can lead to an explosion, heat generation, fire, or smoke development.

NOTES ON SAFETY

Please read this section carefully to familiarize yourself with features and operations before using the unit.

- The warning signs and the sample icons shown here are listed in order for you to use this product safely and correctly as well as to prevent product damage, risk and injury to you or others.

 PRECAUTION WARNING	Indicate the right condition to use the product and prevent damage risk and injury.
 IMPORTANT FOR USERS	Indicate the important elements users should read before using the product.
 CARE AND MAINTENANCE	Indicate matters in which the possibility of damage may happen as a result of incorrect handling and improper maintenance.

TERRAILLON CONNECTED BODY COMPOSITION ANALYZER

⚠ PRECAUTION WARNING

- Do not use the scale on subjects who have a pacemaker or implantable electronic device. It may cause the device to malfunction. When in doubt, please consult your physician. For pregnant women and people who have body implants such as artificial limbs, contraceptive devices, metal plates or screws, the results given by the Body Fat Analyser may be inaccurate. However there is no health risk associated with using the Body Fat Analyzer.
- Do not disassemble the scale as incorrect handling may cause injury.


⚠ IMPORTANT FOR USERS

- Body fat percentage estimates will vary with the amount of water in the body, and can be affected by dehydration or over-hydration due to such factors as alcohol consumption, menstruation, illness, intense exercise, etc.
- Do not use on pregnant women. The result is inaccurate.

⚠ CARE AND MAINTENENCE

1. Do not disassemble the scale; it contains no user serviceable parts. Damage to the scale may occur as a result of improper handling.
2. Clean the scale after use with a dampened cloth. Do not use solvents or immerse the unit in water.
3. Avoid excessive impact or vibration to the scale, such as dropping it onto the floor.
4. For body fat and body water estimates, subject must always be barefoot on the electrodes.

PRODUCT OVER VIEW

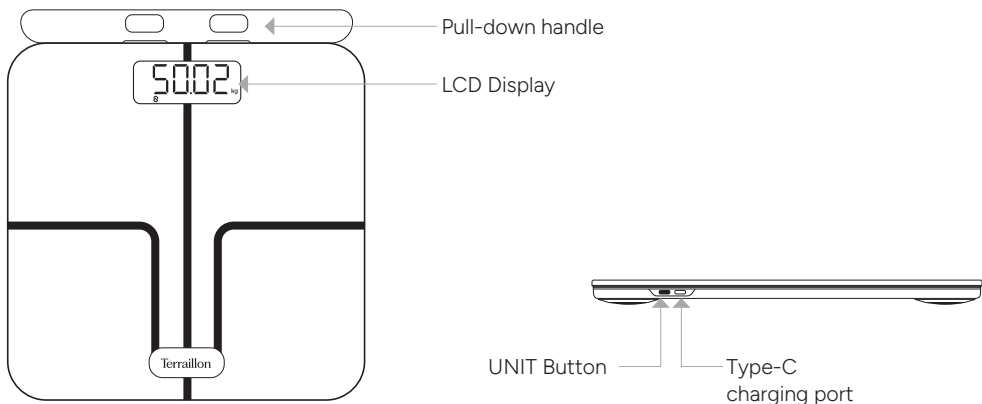TERRAILLON CONNECTED BODY COMPOSITION ANALYZER

⚠️ PREPARATION BEFORE USE

Charging the battery for first-time use.

To charge the scales, plug the supplied cable (to be used with a USB charger, not supplied) into the connector on the side of your scales.

⚠️ Charge your scale when **Lo** is displayed.

Once charging is complete, unplug the charger. Leaving the product charging for too long is not recommended as this may damage the battery.

If you have any difficulties or need more information, you can find our help centre by scanning the QR code opposite with your smartphone. If you cannot find the answer to your question, you can also contact our After Sales Service: support.uk@terraillon.com

FIRST USE

- You will need your smartphone to install your bathroom scales.
- Download the MyBodyCheck Terraillon, then open it.


- Follow the instructions that appear on your smartphone. Ensure you enter all the information requested.
- Your smartphone must remain near to the bathroom scales throughout the installation process.
- If you encounter any problems, please consult the help centre:

TERRAILLON CONNECTED BODY COMPOSITION ANALYZER

USE OF THE APPLIANCE


To use the device, stand on the scale barefoot and pull the handle. Make sure your fingertips touch the electrodes as on the pictures.

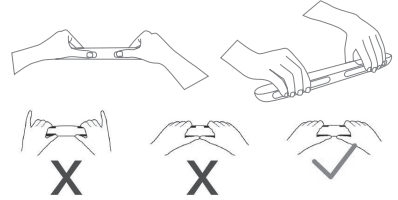Straighten your arms and keep them 45° from your body.

Standing up straight with your feet apart, make sure that the front of your feet cover upside electrodes and that your back of the feet cover downside electrodes.

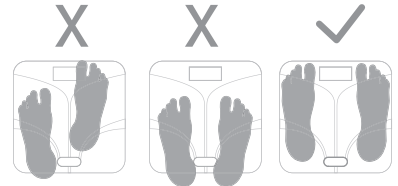When using the scale for the first time after pairing, it is essential to keep the application open in order to synchronise the first data. To obtain accurate readings, keep your posture stable on the scale for about 30 seconds while you take the measurement.


At the end of the measurement, you can check your weight (1) on the scale screen, as well as your BMI (2) and body fat (3). You can also find more details on your MyBodyCheck app. From the 2nd use, the device supports offline mode. No need to open the app for each measurement.

TERRAILLON CONNECTED BODY COMPOSITION ANALYZER

MULTI-USER MANAGEMENT

At the end of the first pairing process, your MyBodyCheck application will ask you to become the owner of the device. There can only be one owner of the device. The owner has the rights to modify certain parameters, such as the device configuration. To invite a new user while remaining the owner of your Master Coach Expert, you can follow the steps below to display the QR code for pairing the invited user. The guest user can add a device by scanning the QR code (Add device / Scale / Connection via Scan code). The guest user's profile is then associated with the bathroom scale.

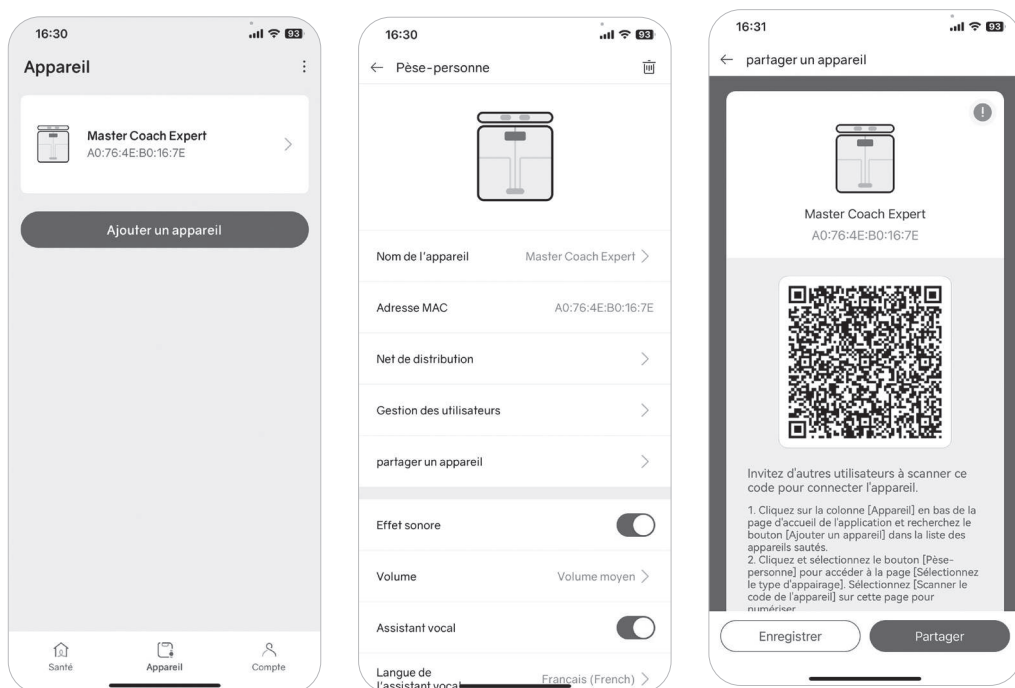PROTECT THE ENVIRONMENT

After life, dispose of this product at a designated waste recycling point.

E.I.R.P. : 20 dBm

GUARANTEE

This product is guaranteed against defects in materials or manufacturing. During the guarantee period, any such defects will be repaired free of charge (proof of purchase must be shown in the case of a complaint under guarantee).

The guarantee does not cover damage resulting from improper use or negligence. In the event of a claim, first contact the shop where you purchased the product.

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iPhone is registered trademark of Apple Inc.
Android is a trademark of Google Inc.
The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.
Other trademarks and trade names are those of their respective owners.
Terrailon SAS France & Headquarters
Terrailon Asia Pacific Ltd