
QUICK START GUIDE

Terraillon®

FR GUIDE DE DÉMARRAGE RAPIDE
IT GUIDA RAPIDA INTRODUTTIVA

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
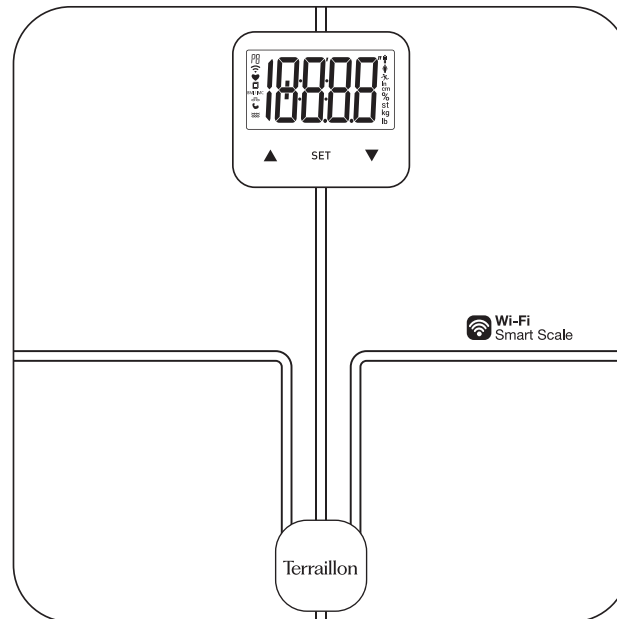
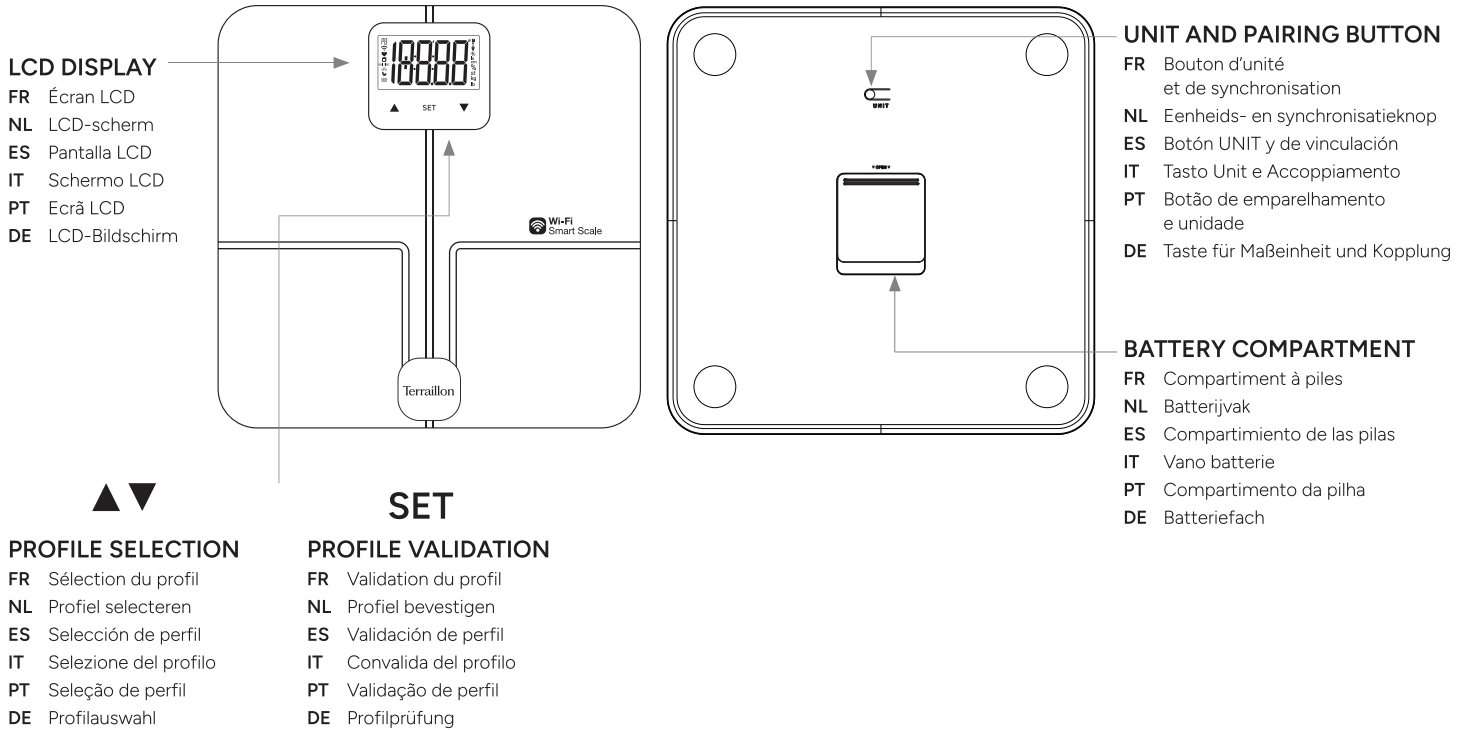
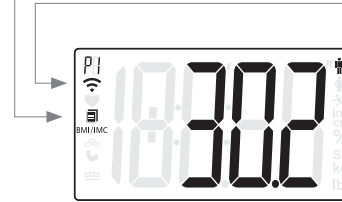
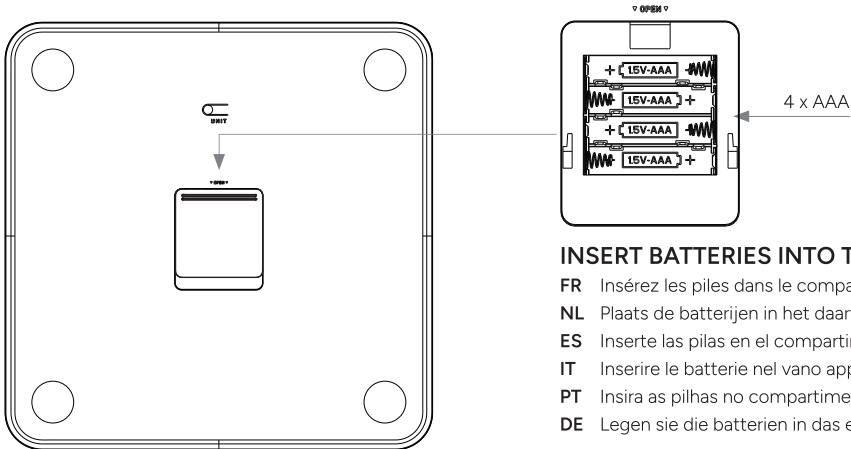
INSERT BATTERIES INTO THE DEDICATED COMPARTMENT.

- FR Insérez les piles dans le compartiment dédié.
- NL Plaats de batterijen in het daarvoor bestemde vak.
- ES Inserte las pilas en el compartimento específico.
- IT Inserire le batterie nel vano apposito.
- PT Insira as pilhas no compartimento dedicado.
- DE Legen sie die batterien in das entsprechende fach ein.

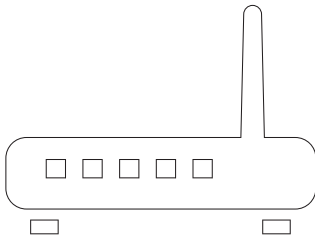GO OVER TO YOUR INTERNET ROUTER AND MAKE SURE YOU KNOW YOUR Wi-Fi PASSWORD. THEN CONNECT YOUR SMARTPHONE TO YOUR Wi-Fi NETWORK.

- FR Rapprochez-vous de votre box internet et assurez-vous de connaître votre mot de passe Wi-Fi. Connectez ensuite votre smartphone à votre réseau Wi-Fi.
- NL Zorg dat u zich in de buurt van uw internetkastje bevindt en uw wifiwachtwoord bij de hand hebt. Verbind uw smartphone vervolgens met uw wifinetwerk.
- ES Acérquese a su router y asegúrese de tener a mano la contraseña del Wi-Fi. A continuación, conecte su smartphone a la red inalámbrica.
- IT Avvicinarsi al modem Internet e assicurarsi di conoscere la password della rete Wi-Fi. Quindi, connettere il proprio smartphone alla rete Wi-Fi.
- PT Aproxime-se da sua box de internet e assegure-se de que sabe a sua palavra-passe do Wi-Fi. Em seguida, ligue o seu smartphone à sua rede de Wi-Fi.
- DE Sie müssen sich in der Nähe Ihres Routers befinden und Ihr WLAN-Passwort kennen. Verbinden Sie Ihr Smartphone mit Ihrem WLAN-Netz.

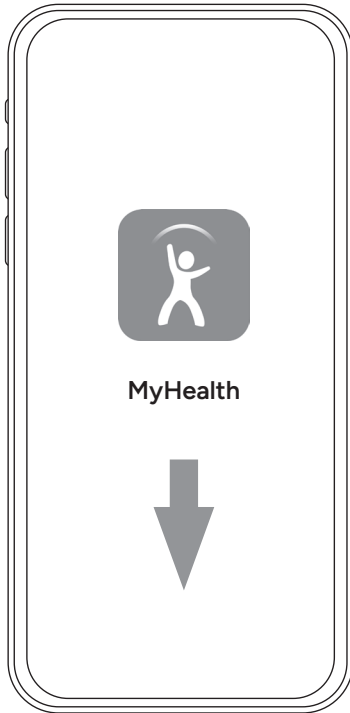DOWNLOAD MyHealth MOBILE APPLICATION.

- FR Téléchargez l'application mobile MyHealth.
- NL Download de mobiele applicatie MyHealth.
- ES Descarga la aplicación móvil MyHealth.
- IT Scaricare l'app mobile MyHealth.
- PT Descarregue a aplicação móvel MyHealth.
- DE Laden Sie die mobile App MyHealth herunter.

CREATE YOUR ACCOUNT OR LOG IN WITH YOUR CREDENTIALS.

- FR Créez votre compte ou connectez-vous avec vos identifiants.
- NL Maak een account aan of meld u aan met uw inloggegevens.
- ES Cree una cuenta o conéctese con sus credenciales.
- IT Creare un account o accedere con i propri dati di accesso.
- PT Crie a sua cunta ou ligue-se com os seus identificadores.
- DE Erstellen Sie Ihr Konto oder melden Sie sich mit Ihrem Login an.

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DURING THE CREATION OF YOUR ACCOUNT, AND WHEN THE APPLICATION INVITES YOU TO DO SO, SELECT THE WI-FI BATHROOM SCALES TO START THE PAIRING PROCEDURE.

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NOTES: IF YOU ARE CONNECTED WITH A PRE-EXISTING ACCOUNT, PRESS THE BUTTON AS INDICATED, THEN "ASSOCIATE A NEW DEVICE". WHEN THE APPLICATION INVITES YOU TO DO SO, SELECT THE WI-FI BATHROOM SCALES TO START THE PAIRING PROCEDURE.

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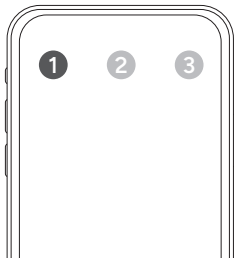
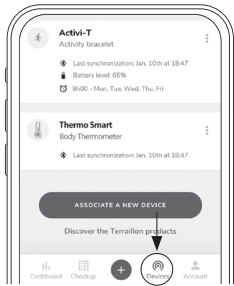
FOLLOW THE INSTRUCTIONS ON THE MOBILE APPLICATION.

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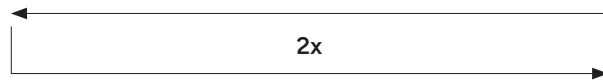
HEART RATE - ACTIVATION AND DESACTIVATION

FR Fréquence cardiaque – Activation et désactivation de la fonction - NL Hartslag - Activeren en uitschakelen van de functie

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EN To disable the heart rate: press the plate to turn on the display (0.0kg), hold SET until OFF appears, SET to confirm the choice. If the function is disabled, the heart rate will not be displayed during the body composition measurement. Note that if the function is deactivated, it is deactivated for all profiles using the product. (Note that when the heart rate is activated, the measurement time is longer.) The heart rate function is only available on the product.

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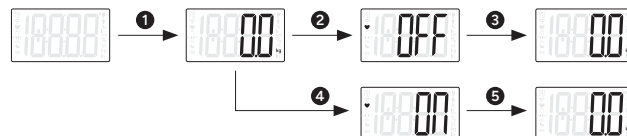
EN ① Return to zero after starting up ② Hold set button ③ Heart rate is enabled by default ④ When the heart rate is off ⑤ Display 2 seconds after return to zero.

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IT Note sulla sicurezza - PT Notas sobre segurança - DE Sicherheitshinweise

EN

Please read this section carefully to familiarize yourself with features and operations before using the unit. The warning signs and the sample icons shown here are listed in order for you to use this product safely and correctly as well as to prevent product damage, risk and injury to you or others. This appliance can be used by children aged from 8 years (body fat analysis results are valid only from 18 years) and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

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EN

PRODUCT SPECIFICATION

Wi-Fi 802,11 b/g/n (2,4 Ghz compatible) - 16.51dBm - 17.94dBm
 Weight & BMI (5 kg-180 kg)
 Body fat analysis (body fat, muscle mass, body water and bone mass*)
 Heart rate measurement function
 8 user memories (with automatic recognition)
 Battery powered (4 x AAA alkaline batteries (not included))
 Dimensions: 32,7 x 32,7 x 2,19 cm

MAINTENANCE

Do not disassemble the scale: it contains no user serviceable parts. Do use a wet soft cloth, dipped into water and wrung out, to wipe the dirt. Then use a dry soft cloth to dry up the device. DON'T use propellant, abrasive or other chemicals.

PROTECT THE ENVIRONMENT

Dispose of used batteries in a designed container so that they can be collected and recycled. After life, dispose of this product at a designated waste recycling point.

WARNING

Do not expose batteries to high temperatures on account of the risk of liquefaction, vaporisation and explosion. Please ensure the batteries are replaced with batteries of the same type (otherwise your device may be damaged).

DECLARATION OF CONFORMITY

Hereby, Terrailon SAS declares that MASTER COACH ULTRA complies with the essential health and safety conditions. This product meets the essential requirements and other relevant provisions of the directive RED 2014/53/UE. The complete declaration of conformity is available on <https://terrailonhelp.zendes.com/hc/en-us> (Declaration of Conformity section).

WARRANTY

This product is guaranteed against defects in materials or manufacturing for 2 years. During this period, any such defects will be repaired free of charge (proof of purchase must be shown in the case of a complaint under warranty). The warranty does not cover damage resulting from improper use or negligence. In the event of a claim, first contact the shop where you purchased your product.

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EN *available via the MyHealth application only

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EN If you have any difficulties or need more information, you can find our help centre by scanning the QR code opposite with your smartphone. If you cannot find the answer to your question, you can also contact our After Sales Service: support.uk@terraillon.com.

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Terrillon SAS France & Headquarters

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Terrillon Asia Pacific Ltd

4/F, Eastern Centre
1065 King's Road
Quarry Bay – Hong Kong
Tel: +852 (0)2960 7200
customerservice@terrillon-asia.com

Terrillon Corp USA

contact@terrillon.fr

Terrillon UK

support.uk@terrillon.com

Terrillon®

MASTER COACH ULTRA