Vibrosit Gimnasia pasiva

Manual de instrucciones

ES

Contenidos

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EN

Contents

Intended Use	3
Safety Notices	3
Items Supplied	5
Device Overview	5
Use	6
Changing the Battery of the Remote Control	6
Cleaning and Storage	7
Troubleshooting	7
Technical Data	8
Warranty	8
Disposal	8

**IMPORTANT NOTICES!
PLEASE KEEP FOR REFERENCE!**

Dear Customer,

We are delighted that you have chosen this **massage device**.

Before using the device for the first time, please read the operating instructions carefully and keep them for future reference.

The manufacturer and importer do not accept any liability if the information in these operating instructions is not complied with.

Safety Notices

- The device is intended for massaging and stimulating your feet and leg muscles and is used for relaxation. It is not a replacement for medically required treatment.
- The device is only suitable for use when sitting or lying down.
- The device is not designed to be operated with an external timer or a separate remote control system.
- The device is for personal use only and is not intended for commercial applications.
- The device should only be used as described in the operating instructions. Any other use is deemed to be improper.
- The warranty does not cover faults caused by incorrect handling, damage or attempted repairs. The same applies to normal wear and tear.
- This device should not be used by people (including children) with impaired sensory or mental faculties or who do not have adequate experience and / or knowledge of such devices, unless they are supervised by a person responsible for their safety, or have received instructions from such a person on how to use the device.
- Children** must not play with the device.
- Cleaning and **user maintenance** may not be performed by **children**.
- Do not use the mains adapter with another device and never use the device with another mains adapter.
- Check the device for damage each time before you use it. The device may only be used if the device and connecting cable have no visible signs of damage! If the device or the connecting cable are damaged.

Dangers for the Following Groups of People

The following people may use the device **only after consultation with their doctor**:

- people who wear a pacemaker or similar medical implants.
- people with heart disease or who suffer from epilepsy.
- people who have back complaints or spinal disorders.
- women who are pregnant.
- people who suffer from circulation problems or thrombosis.
- people who have undergone an operation for which the treatment is not yet complete.
- people with physical abnormalities.
- people who have a malignant tumour.
- Do not use the device near to or on injuries, swelling, skin irritations, infected skin, warts, fungal infections, abnormal changes in the skin or body parts with no feeling.

DANGER – Danger of Electric Shock

- Use and store the device only in closed rooms. Do not operate the device in rooms with high humidity.
- Never immerse the device in water or other liquids! Ensure that the device and the connecting cable cannot fall into water or become wet.
- Should the device fall into water, switch off the power supply immediately. Do not attempt to pull the device out of the water while it is still connected to the mains power!
- Never touch the device and the connecting cable with wet hands when these components are connected to the mains power.
- Do not use the device if it has malfunctioned or been dropped in water. Have it checked in a specialist workshop before using it again.
- The device must always be disconnected from the mains power when left unattended and before cleaning it.

WARNING – Danger of Injury

- Danger of suffocation!** Keep children and animals away from the packaging material.
- Children younger than 8 years of age and animals must be kept away from the device and the connecting cable.
- Danger of burns from battery acid!**
 - Batteries can be deadly if swallowed. The battery and remote control should therefore be kept out of the reach of children and animals. If a battery is swallowed, immediate medical attention must be sought.
 - If the battery has leaked, do not allow your skin, eyes and mucous membranes to come into contact with the battery acid. Wear protective gloves if necessary. If you touch battery acid, rinse the affected areas with abundant clear water straight away and seek medical assistance immediately.
- Danger of tripping!** Make sure that the connecting cable does not present a trip hazard.

WARNING – Danger of Fire

- Do not use the device in rooms containing easily ignitable or explosive substances.
- Do not cover the device during operation in order to prevent it from catching fire.
- Unwind the connecting cable fully before connecting to the mains power.

NOTICE – Risk of Damage to Material and Property

- Only connect the mains adapter to a properly installed plug socket. The plug socket must also be readily accessible after connection so that the connection to the mains can quickly be isolated. The mains voltage must match the rating plate of the device and the mains adapter.
- The device must be switched off when the mains adapter is pulled out or inserted into the plug socket.
- Make sure that the connecting cable is not squashed, bent or laid over sharp edges and does not come into contact with hot surfaces.
- Switch off the device and pull the mains adapter out of the plug socket if an error occurs during operation or before a thunderstorm.
- When the mains adapter is removed from the plug socket, always pull on the mains adapter and never the mains cable.

- Do not pull or carry the device by the connecting cable. Never wind the mains cable around the device or mains adapter.
- Use only the battery type (CR2025).
- Remove the battery when it is dead or if you will not be using the device for a prolonged period. This will prevent any damage caused by leaks.
- Remove a leaking battery immediately.
- Batteries must not be charged or reactivated in any way, disassembled, thrown into a fire, immersed in liquids or short-circuited.
- Use only original accessories from the manufacturer in order to guarantee that there is no interference that may prevent the device from working and to avoid possible damage.
- Never place heavy objects on the device or the connecting cable.
- Protect the device from heat sources, naked flames, sub-zero temperatures, persistent moisture, wet conditions and impacts.


Items Supplied

- 1 x massage device
- 1 x mains adapter
- 1 x remote control (not included cell battery CR2025)
- 1 x operating instructions

Check the items supplied for completeness and the components for transport damage. If you find any damage, do not use the device but contact our customer service department.


Remove any possible films, stickers or transport protection from the device. **Never remove the rating plate and any possible warnings!**

Device Overview

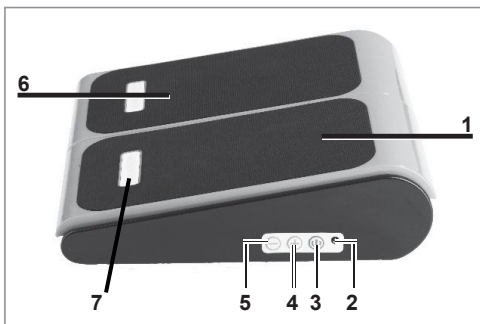Main Machine:

- 1 Vibrating surface on right
- 2 Connecting socket for the mains adapter
- 3 Button 
- 4 Speed Button +
- 5 Speed Button -
- 6 Vibrating surface on left
- 7 Display (Time and Temperature)

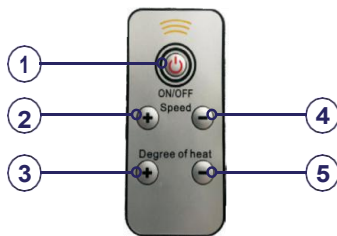Remote control:

1. Power
2. Speed up
3. Increase temperature
4. Speed down
5. Reduce temperature

→ Adapter

Use

PLEASE NOTE!

-No pressure should be exerted on the device as it may be damaged. Never stand on the device.

-CAUTION – Danger of Crushing! Never place your hands and feet between the frame and vibrating surfaces while the device is in operation.

-Do not use the device for longer than 30 minutes without stopping.

-Do not fall asleep while you are using the device!

-A tingling sensation / itching during or after the massage is normal due to the circulation-boosting effect.

1-Pull the insulation strip out of the battery compartment of the remote control to activate the battery.

2-Place the massage device on the floor. The surface should be firm, dry and flat.

Please Note: Protect any delicate surfaces (laminated, parquet floor, etc.) from scratches etc. with a protective mat.

3-Connect the barrel connector of the mains adapter to the connecting socket (2) on the device and insert the mains adapter into an easily accessible plug socket.

4-The device can be used either while sitting or lying down:

Sitting: Sit down on a chair, armchair or similar place in front of the device. Place your feet down on the vibrating surfaces (1, 6).

Lying down: Lie down in front of the device in such a way that your calves can be placed on the vibrating surfaces.

5-Switch on the device by pressing the button  (3) on the device or on the remote control.

6-Set the vibration level (1 – 10) you want by pressing the buttons + (4) and – (5) on the device or on the remote control.


The device switches off automatically after approx. 30 minutes. **Please Note:** It can be switched off manually at any time using the button  on the device or on the remote control.

7-Set the temperature level (30 – 50°C) you want by pressing the degree buttons + (3) and – (5) on the remote control. (Buttons on the side of the machine can't control the temperature.)

8-Pull the mains adapter out of the plug socket and disconnect the barrel connector from the device.

Install and Changing the Battery of the Remote Control

Before inserting the battery, check whether the contacts on the remote control and on the battery are clean, and clean them if necessary.

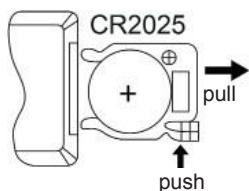1. The battery compartment is located on the narrow side of the remote control. Pull it out as shown.
2. Take out the used battery and insert a new battery (3 V, CR2025). The positive pole must be facing upwards.
3. Push the battery compartment back into the remote control to the limit stop. You should feel and hear it click into place.

Replace the battery when the power of the remote control starts to wane.

Cleaning and Storage

PLEASE NOTE!

- Do not use any caustic or abrasive cleaning agents or cleaning pads to clean the device. These may damage the surfaces.
- Make sure that the device is not connected to the mains power.
- Wipe down the device with a slightly damp cloth if necessary. Then wipe it down with a dry cloth.
- Keep the device in a dry place which is protected from sunlight and is not accessible to children or animals.

Troubleshooting

If the device does not work properly, first check whether you are able to rectify the problem yourself. If the problem cannot be solved with the steps listed, contact the customer service department.

Do not attempt to repair an electrical device yourself!

Problem	Possible cause	Solution
The device does not work.	The mains adapter is not inserted correctly in the plug socket.	Push the mains adapter into the plug socket as far as it will go.
	The device is not connected to the mains adapter correctly.	Correct the connection.
	The plug socket is defective.	Try another plug socket.
	There is no mains voltage.	Check the fuse of the electrical distributor (fuse box).
	The button  (3) has not been (correctly) pressed.	Press the button  (3) again.
The device does not react to the remote control.	The battery of the remote control is too weak.	Replace the battery of the remote control.
	There are obstructions between the device and remote control (no visual contact).	Point the remote control directly at the device.

Technical Data

XJ-F-26D	
Power	48W
Voltage	24V 2A
Max weight	40KG
Amplitude	4mm
Motor Speed	1400 rpm - 3200 rpm
Level	0-10
Heating	30 °C – 50°C
Product Size	35.5 x 33 x 10.5 cm
Packing Size	40.5 x 38.5 x 15.5 cm
N.W. / G.W.	5 kg / 5.8 kg

Warranty

In the following cases, the guarantee according to current regulations will not be effective:

- Any damage resulting from repairs or manipulation unauthorized by the company.
- Any damage resulting from misuse, negligence or accident.
- Any damage caused by natural disasters or abnormal voltage.
- Any damage resulting from cases of use of the machine in contexts other than those described in this manual.

Recycling

IMPORTANT INFORMATION ON THE CORRECT DISPOSAL OF THE PRODUCT REGARDING THE STANDARD OF THE EUROPEAN UNION 2002/96/EC

The waste materials of the electrical products must not be dumped with the rest of rubbish waste of the household.

These materials shall be returned to the shop of taken to a recycling point.

Dear Customers, for any doubt related to this product, you may contact our Customer

Service. We are more than happy to help you.

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