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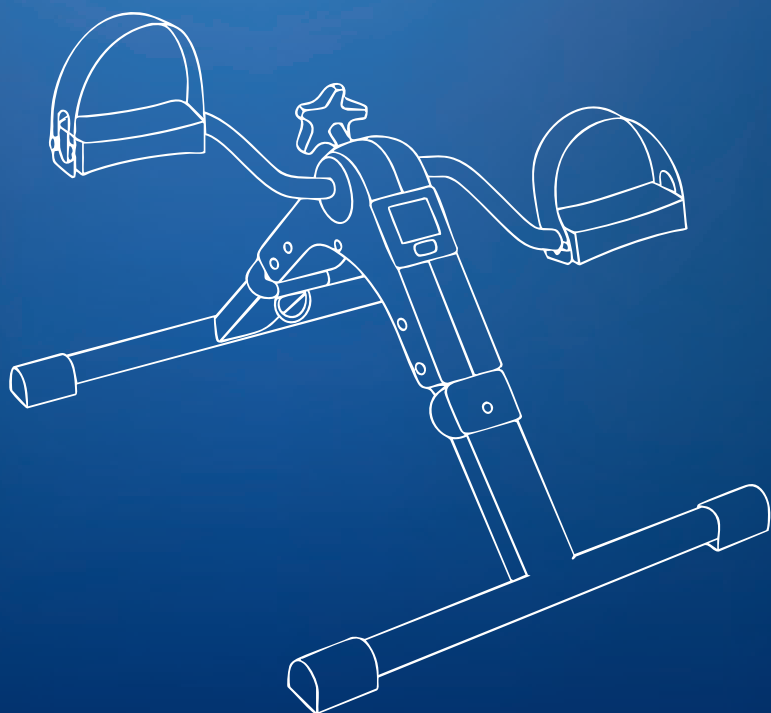
## IMPORTANT

Keep the instructions for future reference. Read these instructions before assembling and using this product. Adult assembly is required.

Thank you for choosing our RONDA folding electronic pedal exerciser! In this manual, you will find all the information you need to assemble, use and maintain this product correctly and safely.

This device has been designed to provide a comfortable, efficient pedalling experience that can be adapted to a variety of needs.

Before you begin, we recommend that you read this manual carefully to ensure optimal and safe use of the product.

As a guarantee of the quality of our products, the manufacturer has a management system certified under ISO 13485. This product complies with the requirements of Regulation (EU) 2017/745 on medical devices.

Notice to the user and/or patient: any serious incident that has occurred in relation to the product must be reported to the manufacturer and to the competent authority of the Member State in which the user and/or patient is established.

## 1. WARNINGS AND PRECAUTIONS

- 1. WARNING!** Do not use without medical supervision in specific cases: Persons with heart problems, recent injuries or reduced mobility should consult a doctor before use.
- 2. WARNING!** Do not use on unstable or sloping surfaces.
- 3. WARNING!** Do not overload the device: The bottom bracket is designed for home use, not for heavy-duty training.
- 4. WARNING!** Avoid contact with water and moisture.
- 5. WARNING!** To avoid injury, wear closed footwear.
- 6. WARNING!** Keep out of reach of unsupervised children: This is not a toy and its misuse can cause accidents.

## 2. PURPOSE OF THE PRODUCT

- 1. Physical exercise and health maintenance:** It allows physical activity to be carried out anywhere, as well as contributing to the improvement of cardiovascular and muscular endurance.
- 2. Personalised training:** It offers different levels of resistance to adapt the exercise to each user.
- 3. Ergonomic and accessible use:** It can be used in different positions: sitting or lying down, and its compact design makes it easy to transport and store.
- 4. Improved general well-being:** Reduces stress and improves mental health through regular physical activity and helps to combat sedentary lifestyles in people who work long hours sitting down.

### 3. BENEFITS AND CONSIDERATIONS

There are multiple benefits to using the bike trainer, such as improved mobility, blood circulation and muscle strengthening, and it's a great way to enjoy a moderate cardio workout.

Its compact and versatile design allows it to be used at home or anywhere, adapting to different fitness levels. In addition, with the electronic panel you can enjoy monitoring functions to follow your training progress, contributing to a more active and healthy lifestyle.

For safe and efficient use, it is important to place it on a stable surface and adjust the resistance progressively.

### 4. PRODUCT FEATURES

<b>PRODUCT NAME</b>	Foldable electronic pedal set
<b>MODEL</b>	RONDA
<b>MODEL REFERENCE</b>	A-00053/04
<b>MAXIMUM SUPPORTED WEIGHT</b>	264,55 lb
<b>NET PRODUCT WEIGHT</b>	4,40 lb
<b>PRODUCT DIMENSIONS</b>	19,48"*13,97"* 9,44"
<b>MATERIAL</b>	Steel

### 5. PARTS AND COMPONENTS

## 6. INDICATIONS AND CONTRAINDICATIONS

### INDICATIONS

This product is intended for:

- People in rehabilitation who need to regain mobility in legs or arms after an injury, surgery or accident.
- Older adults to help maintain physical activity, improve circulation and strengthen muscles.
- Users with a sedentary lifestyle: perfect for those who work long hours sitting down and want to activate circulation and prevent muscle problems.
- People looking for risk-free exercise.

### CONTRAINDICATIONS

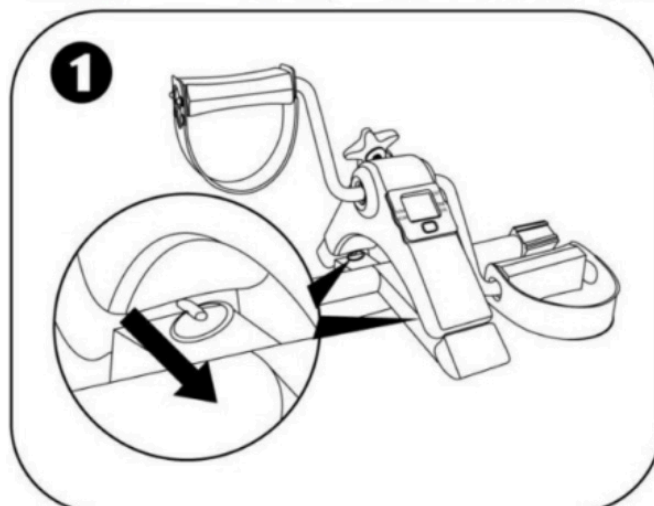
This product is not intended for:

- Persons with recent fractures or serious injuries.
- Patients with severe cardiac or respiratory conditions in a relaxed manner.
- Users with severe neurological disorders without medical control.
- Unsupervised children.


## 7. INSTRUCTIONS FOR USE

You will find your bottom bracket folded. To unfold it:

**Step 1:** Pull firmly on the ring shown in the picture to allow the product to unfold.

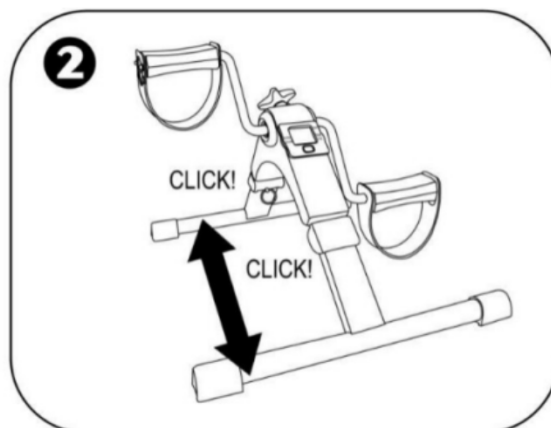
**Step 2:** Then pull the bottom bracket outwards until it clicks into place.


**COMMISSIONING OF THE ELECTRONIC PANEL:**

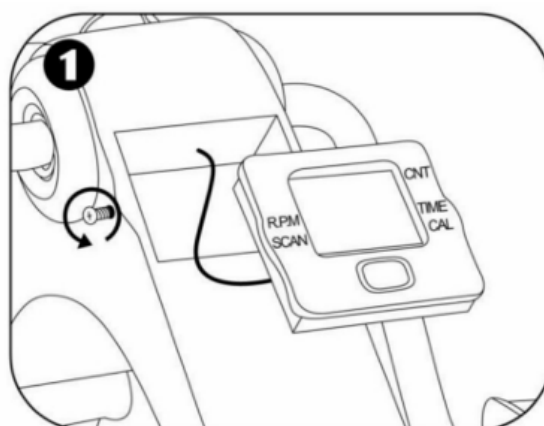
**Step 1:** Loosen and remove the screws on each side of the display and then carefully remove the display from the hole in which it is located.

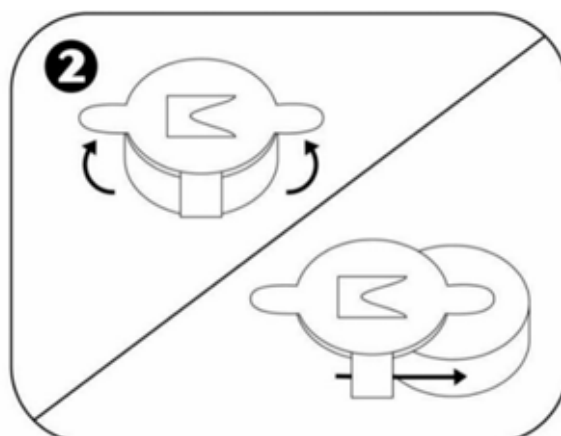
**Step 2:** On the back of the display unit, lift up the battery tabs.

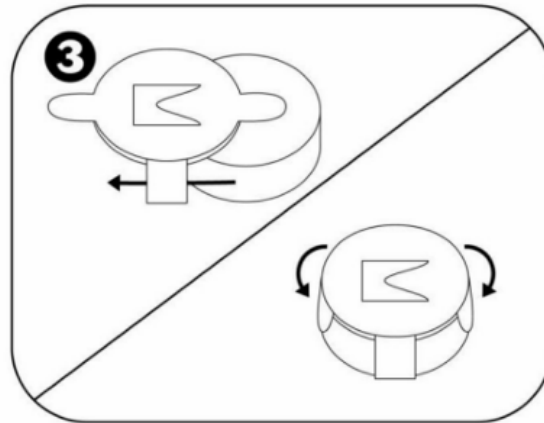**Step 3:** Insert a new 1.5 V AG13 battery, positive side down, and close the tabs.



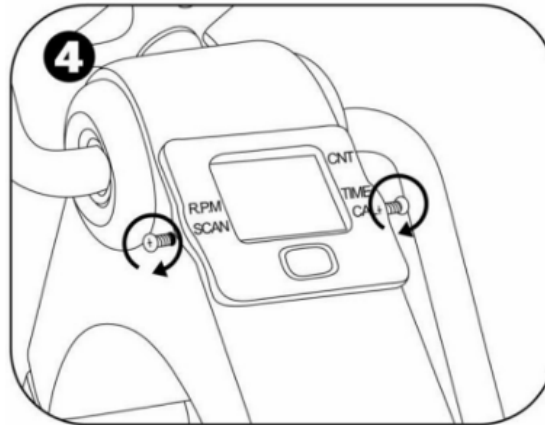
**Step 4:** Put the display back in place and secure it again with the screws.

**FUNCTIONS OF THE ELECTRONIC PANEL:**

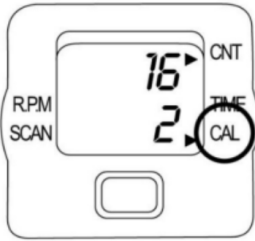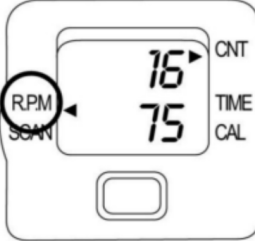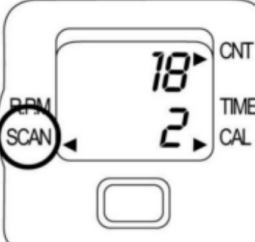	<p>CNT: Number of completed rotations.</p>
	<p>TIME: Duration of training.</p>

	<p>CAL: Estimated calorie expenditure.</p>
	<p>R.P.M.: Rotations per minute.</p>
	<p>SCAN: Toggle the TIME - CNT - RPM - CAL functions on the display.</p>

## 8. INSTRUCTIONS FOR USE

To turn on the display, press the centre button once or start pedalling.

The display switches off automatically after 4 minutes of inactivity.

All display functions are activated as soon as you start pedalling and stop as soon as you stop pedalling.

## 9. MAINTENANCE AND CLEANING

- Clean the structure with a soft cloth to remove dust and dirt.
- Do not use abrasives or solvents that could damage the materials.
- Constant exposure to moisture should be avoided.

## 10. PROBLEM SOLVING

**Problem: The screen does not turn on.**

- Possible cause: the batteries have stopped working
- Solution: Check their condition and replace them if necessary.

**Problem: Unusual noises when pedalling.**

- Possible cause: The bottom bracket is not fully extended.
- Solution: Check if the product has been deployed correctly and is fixed in its position of use

**Problem: Pain or discomfort when pedalling.**

- Possible cause: Incorrect posture or incorrectly adjusted resistance.
- Solution: Make sure you are in a suitable position (sitting or lying down) and adjust the resistance to a comfortable level.

## 11. CONDITIONS OF GUARANTEE

This product is covered by a 36 month warranty against manufacturing defects.

Thank you very much for choosing Mobiclinic for this RONDA folding electronic bottom bracket!

We value your trust in our products and are committed to always offering you the best quality. If you have any questions or need additional assistance, please do not hesitate to contact us through our customer service at [clientes@mobiclinic.com](mailto:clientes@mobiclinic.com). We appreciate your comments and suggestions so that we can continue to improve.

Thank you for being part of Mobiclinic!

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**XIAMEN NEAQ IMP.& EXP. CO., LTD**
**SUNGO EUROPE B.V.**
**GRUPO R. QUERALTÓ S.A.**