
MANUAL DE INSTRUCCIONES

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MOBIDAL

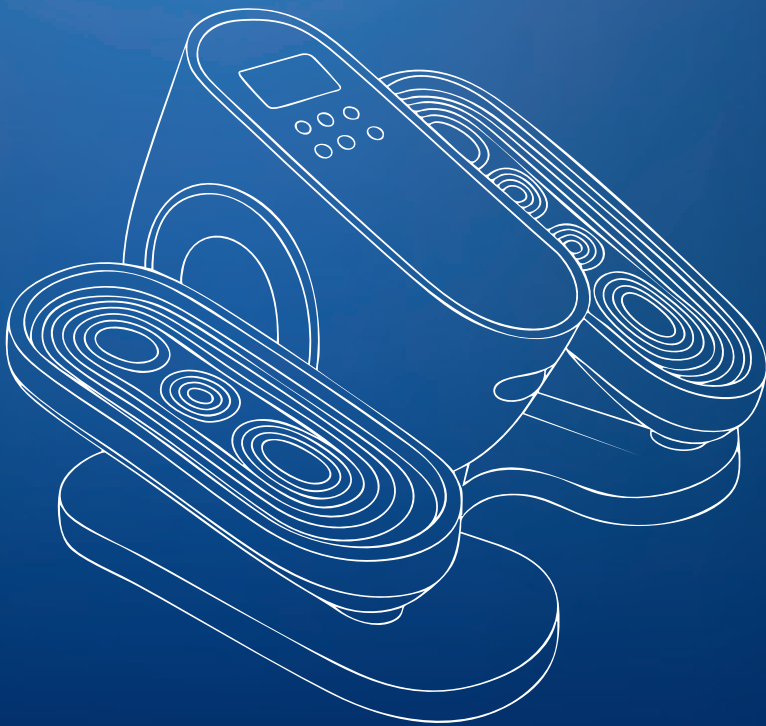Manual de instrucciones	ES
Instruction manual	EN
Manuel D'Utilisation	FR
Manuale D'Istruzioni	IT
Anweisungen	DE
Instruções	PT
Gebruiksaanwijzing	NL
Bruksanvisning	SW
Instrukcja obsługi	PL
Brugsanvisning	DA

ES - Lea las instrucciones antes de usar el producto.
EN - Please read instructions before use.
FR - Lisez les instructions avant d'utiliser le produit.
IT - Leggere le istruzioni prima di usare il prodotto.
DE - Lesen Sie die Gebrauchsanweisung, bevor Sie das Produkt verwenden.
PT - Leia as instruções antes de utilizar o produto.
NL - Lees de instructies voordat u het product gebruikt.
SW - Läs instruktionerna innan du använder produkten.
PL - Läs instruktionerna innan du använder produkten.
DA - Læs vejledningen, før du bruger produktet.

IMPORTANTE

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**IMPORTANT**

Keep these instructions for future reference. Read these instructions before assembling and using this product. Adult assembly is required.

Thank you for choosing our MOBIDAL mini elliptical trainer! In this manual, you will find all the information you need to assemble, use and maintain this product correctly and safely.

The mini elliptical trainer is designed to improve overall physical condition through low-impact exercise that activates the lower body and improves blood circulation. This product allows you to train comfortably at home, strengthening muscles, improving joint mobility and helping to burn calories, all without overloading the joints.

Before you begin, we recommend that you read this manual carefully to ensure optimal and safe use of the product.

1. PRODUCT PURPOSE

1. Improve overall physical condition through low-impact exercise.
2. Stimulate blood circulation in the legs and feet.
3. Strengthen and activate the lower body muscles.
4. Promote joint mobility in the hips, knees, and ankles.
5. Allow training even while seated (office or home).
6. Preventing a sedentary lifestyle and facilitating daily physical activity.
7. Serve as a complementary tool in mild rehabilitation processes.

2. PRODUCT FEATURES

PRODUCT NAME	Mini elliptical trainer
MODEL	MOBIDAL
MODEL REFERENCE	A-00407/11
MAXIMUM WEIGHT SUPPORTED	55,11 lb
NET WEIGHT OF THE PRODUCT	13,66 lb

PRODUCT DIMENSIONS	17,71"*14,96"*10,62"
MATERIAL	ABS plastic and stainless steel
OTHER FEATURES	<ul style="list-style-type: none"> • Power: 60 W • Input voltage: 110-240 V • Nominal frequency: 50/60 Hz • Noise level: < 60 dB • Manual mode (HR) and 3 automatic programmes (P1, P2, P3) • LED display: speed, time, distance, temperatures • Remote control included

3. PARTS AND COMPONENTS

Remote control (function control):

- **ON/OFF:** Switching the device on and off.
- **HR:** Activate manual mode.
- **P1 / P2 / P3:** Selection of automatic programmes.
- **+ / - SPEED:** Speed adjustment
- **+ / - TIME:** Adjust the usage time (5 to 30 minutes).
- **MODE:** Switch between available modes.
- **STOP:** Stop the movement at any time.


- Power cable
- User manual

4. INSTRUCTIONS FOR USE

Place the mini elliptical on a flat, stable surface.

1. Connect the mini elliptical to the power supply using the power cable.
2. Sit opposite her on a chair or sofa with a straight back.
3. Place your feet on the pedals.
4. Press the power button (POWER) on the panel or remote control.
5. Use the speed arrows (+/-) to adjust the intensity.
6. Choose between manual mode (HR) or automatic programmes (P1, P2, P3).

Tip: Maintain an upright posture and perform the movement at a steady pace. You can combine its use with other tasks (reading, watching television, working).

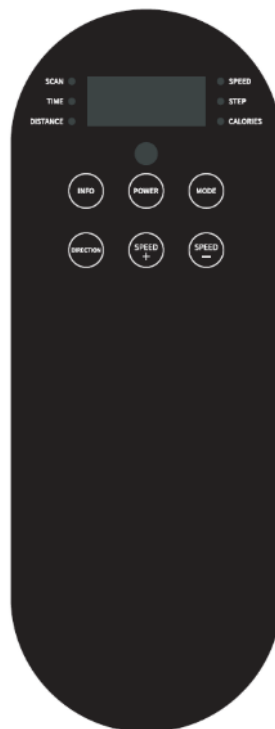MANUAL EXERCISE MODE (HR):

By default, the product starts in manual mode (HR).

1. To activate manual mode, press the MODE button until "HR" appears on the display.
2. Use the +/- SPEED buttons to adjust the speed and +/- TIME to set the exercise time.

3. Press POWER to start the movement.
4. Once in operation, you can press DIRECTION to change the direction of the pedals (forward/backward).

AUTOMATIC PROGRAMMES P1, P2, P3:

1. Press MODE once to select programme P1, twice for P2 and three times for P3.
 2. In these programmes, the speed and direction of the pedals change automatically according to the predefined pattern. Only the time can be adjusted.
- P1 (Daily exercise):
 - Recommended for elderly people or those with reduced mobility.
 - Direction: forwards and backwards (changes every 60 seconds).
 - Speed: cycle 1-2-3-4-5-4-3-2-1-3-5.
 - P2 (Light exercise):
 - Ideal for users who work sitting down or with little activity.
 - Direction: forwards and backwards (changes every 90 seconds).
 - Speed: cycle 2-3-4-5-6-7-8-4-6-8-4.
 - P3 (Intensive training):
 - Aimed at users who want to burn calories and strengthen muscles.
 - Direction: forwards and backwards (changes every 120 seconds).
 - Speed: cycle 5-6-7-8-9-10-8-6-8-10-8.

Note: With the appliance switched off or unplugged, the pedal can be moved forward easily. Moving it backwards will require a little more effort.

5. MAINTENANCE AND CLEANING

- Disconnect the device before cleaning it.
- Clean with a dry or slightly damp cloth.
- Do not use abrasive products or immerse it in water.
- Store in a dry, dust-free place.

6. WARNINGS AND PRECAUTIONS

GENERAL

1. To avoid short circuits, connect the device correctly and keep it away from metal objects, water, and rubbish.
2. Always use the device on a level and stable surface. If it is not, it could tip over or cause injury.
3. Unplug the product when not in use.
4. Store the equipment in its original packaging or in a cool, dry, secure place if you are not going to use it for a long period of time.
- 5.
6. If the power cable is damaged, it must be replaced by qualified personnel.
7. **Do not** use the product or any of its parts if they are damp or wet.
8. Keep the appliance and cable out of reach of children under 8 years of age.
9. To avoid injury, do not use the product while standing or apply excessive force during use.
10. **Do not** attempt to disassemble or repair the device yourself. Repairs must be carried out by a qualified technician.

BEFORE USE

1. **WARNING!** Use the product indoors, on a stable surface.
2. **WARNING!** Do not use outdoors or near sources of heat or moisture.
3. **WARNING!** Do not use the product while standing; it is only to be used while seated, otherwise the product may be damaged.

WHILE IN USE

1. **WARNING!** Do not stand on the pedals.
2. **WARNING!** Do not use the product if you are dizzy, injured, or ill.
3. **WARNING!** Keep away from children and pets.

IMMEDIATE SUSPENSION OF USE

1. **WARNING!** If the appliance emits abnormal sounds or overheats.
2. **WARNING!** If it gets wet or comes into contact with liquids.
3. **WARNING!** If you experience pain, discomfort, or dizziness during use.

ELECTRICAL WARNINGS

1. **WARNING!** Do not touch the plug with wet hands.
2. **WARNING!** Do not pull on the cable to unplug it.
3. **WARNING!** Disconnect the appliance after use or if it is not going to be used for a long period of time.

7. PROBLEM SOLVING

Problem: It does not turn on.

- Possible cause: Cable not properly connected.
- Solution: Check the connection to the mains.

Problem: The remote control is not responding.

- Possible cause: Battery is flat.
- Solution: Replace the battery in the remote control.

Problem: Abnormal noise.

- Possible cause: Use on an unstable surface.
- Solution: Place the equipment on flat ground.

Problem: Stops on its own.

- Possible cause: Scheduled time has ended.
- Solution: Restart the cycle from the panel.

8. WARRANTY CONDITIONS

This product is covered by a 36-month warranty against manufacturing defects.

9. EC DECLARATION OF CONFORMITY

The manufacturer declares that this device complies with the essential requirements and other relevant provisions of the following European directives:

- Directive 2014/30/EU on Electromagnetic Compatibility (EMC).
- Low Voltage Directive 2014/35/EU (LVD).
- Directive 2011/65/EU and its amendment 2015/863/EU (RoHS 2.0).

Thank you very much for choosing Mobiclinic for this MOBIDAL mini elliptical trainer! We value your trust in our products and are committed to always offering you the best quality.

If you have any questions or require further assistance, please do not hesitate to contact us via our customer service department at **clientes@mobiclinic.com**. We appreciate your comments and suggestions so that we can continue to improve.

Thank you for being part of Mobiclinic!

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KIMS INDUSTRY & TRADE CO. LTD
GRUPO R. QUERALTÓ S.A.