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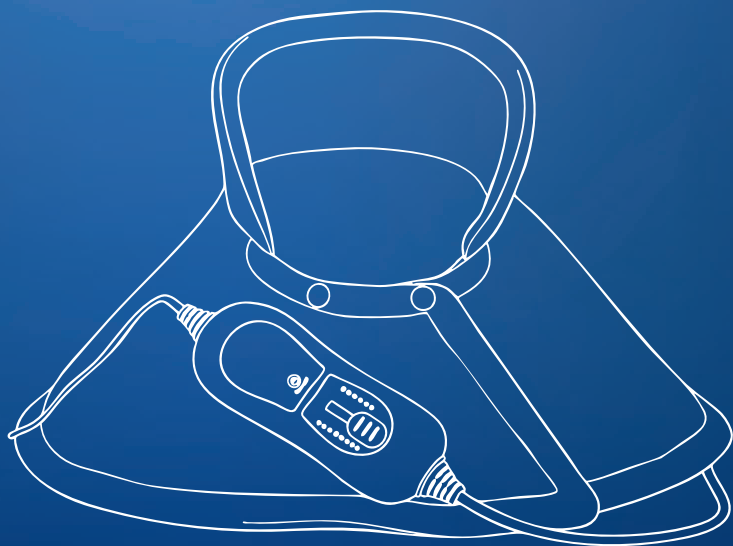
**IMPORTANT**

Keep these instructions for future reference. Read these instructions before assembling and using this product.

Thank you very much for choosing our 24,40"\*23,62" electric neck pad! In this instruction manual, you will find all the information you need to use and maintain this product safely. Please read all instructions carefully before using the product.

Below, you will find a detailed description of their features, purposes, and recommendations for use:

**PRODUCT PURPOSE**

The electric neck pad is designed to provide a comforting source of heat to the neck, shoulders and upper back. Its purpose is to relieve muscle tension, combat the cold and offer a feeling of general well-being by emitting controlled and safe heat.

**1. PRODUCT CHARACTERISTICS**

<b>PRODUCT NAME</b>	Electric neck pillow
<b>MODEL</b>	24,40"*23,62"
<b>REFERENCE</b>	VD-00055/21
<b>PRODUCT WEIGHT</b>	1,65 lb
<b>MATERIAL</b>	Soft and comfortable fleece
<b>OTHER FEATURES</b>	<ul style="list-style-type: none"> <li>• Power: 100 W</li> <li>• Voltage: 220–240 V</li> <li>• Minimum temperature: 53 °C</li> <li>• Maximum temperature: 65 °C</li> <li>• Temperature levels: 3</li> <li>• Automatic switch-off time: 90 minutes</li> <li>• Protection: Overheating safety system</li> <li>• Washable: Hand or machine wash, maximum 30°C (without electric control)</li> <li>• Includes: Electric pad with cable and control unit</li> </ul>

## 2. PARTS AND COMPONENTS

- Electric heating pad.
- Control panel:

- 1) Power off button (OFF).
- 2) Position 1 – 1st temperature level.
- 3) Position 2 – 2nd temperature level.
- 4) Position 3 – 3rd temperature level.

## 3. INSTRUCTIONS FOR USE

- Connect the plug to the socket on the bottom of the pad and the other end to a power supply (AC220-240V).
- To preheat the electric pad, select the third setting for 15 minutes. The temperature will rise gradually.
- After preheating, turn the control knob to select the desired temperature level.

Note: The heating pad is set to operate continuously for 90 minutes. If you wish to extend this time, it is advisable to allow it to cool for 5 minutes before using it again.

## 4. MAINTENANCE AND CLEANING

Follow these instructions when you are not using the electric pad and before cleaning it:

- Ensure that the control button is in the "OFF" position.
- Disconnect the power supply.
- Disconnect the plug that connects the control to the device.
- Allow the heating pad to cool before cleaning it.
- Small stains can be cleaned with a damp cloth or sponge and, if necessary, with a little liquid detergent for delicate fabrics.

**Warning:** The electric pad must not be dry cleaned, folded, tumble dried or machine pressed/ironed.

- If the electric pad is very dirty, it can be washed in a washing machine. You should use the delicate cycle at a temperature of 30°C.
- Use a detergent for delicate fabrics.
- It should not be washed in the washing machine more than 5 times.
- Never plug in the electric pad to dry it.

## 5. WARNINGS AND PRECAUTIONS

1. **WARNING!** Children under three years of age (0 to 3 years) should not use this device.
2. **WARNING!** Do not use the heating pad if it is folded or creased.
3. **WARNING!** Do not use if wet.
4. **WARNING!** The electric and magnetic fields from this heating pad may interfere with the operation of pacemakers.
5. **WARNING!** Do not cover the switch or place it on the pad when in use.
6. **WARNING!** Do not pull, twist or bend the cables.
7. **WARNING!** This pad should be checked regularly for signs of wear or damage.

## 6. WARRANTY CONDITIONS

This product is covered by a 3-year warranty against manufacturing defects.

Thank you very much for choosing Mobiclinic for this 24,40"\*23,62" electric neck pad.

We value your trust in our products and are committed to always offering you the best quality and safety. If you have any questions or need additional assistance, please do not hesitate to contact us at [clientes@mobiclinic.com](mailto:clientes@mobiclinic.com). We appreciate your comments and suggestions so that we can continue to improve.

Thank you for being part of Mobiclinic!

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CO.,LTD**
**GRUPO R. QUERALTÓ S.A.**